DEMONSTRATED PROTOCOL

Nuclei Isolation for Single Cell Multiome ATAC + Gene Expression Sequencing

Overview

This protocol outlines how to isolate, wash, and count nuclei suspensions for use with the Chromium Next GEM Single Cell Multiome ATAC + Gene Expression (GEX) protocol (CG000338). Cryopreserved primary cells (PBMCs) and cell lines (GM12878 cells; 3T3 cells) were used to develop this protocol. PBMCs were cryopreserved in IMDM + 40% FBS + 15% DMSO. Cell lines were cryopreserved in RPMI + 15% FBS + 5% DMSO.



For optimal assay performance, nuclei isolation should be performed using this protocol and not the standalone protocols for nuclei isolation for ATAC or RNA sequencing only. The recommended buffer compositions, final nuclei suspension concentration, and the wash step guidelines presented in this protocol for nuclei sample preparation are critical for optimal Chromium Single Cell Multiome ATAC + GEX assay performance. Failure to adhere to these guidelines may result in suboptimal assay performance.

DO NOT use this protocol for isolating nuclei from tissue-derived cells. Refer to Demonstrated Protocols-Nuclei Isolation from Embryonic Mouse Brain for Single Cell Multiome ATAC + Gene Expression Sequencing (CG000666) & Nuclei Isolation from Complex Tissues for Single Cell Multiome ATAC + Gene Expression Sequencing (CG000375) for isolating nuclei from tissue-derived cells.

Additional Guidance

Consult Demonstrated Protocol Cell Preparation Guide (Document CG00053) for Tips & Best Practices.

Cells carry potentially hazardous pathogens. Follow material supplier recommendations and local laboratory procedures and regulations for the safe handling, storage and disposal of biological materials.

Cell Sourcing

Cell Type	Species	Supplier
GM12878	Human	Coriell Institute
3T3	Mouse	ATCC
Normal Peripheral	Human	AllCells
Blood MNC (PBMC)		

Optimization Recommendations

The following demonstrated protocol was performed using the indicated sample types. Optimization of some protocol steps may be needed for other cell types.

· Lysis time:

Perform a lysis timeline to determine appropriate lysis incubation time for a specific cell type. For optimization experiments, RNAse inhibitor may be omitted from the buffer and instead of the 10x Genomics' Nuclei Buffer, PBS may be used for nuclei resuspension. However for the actual experiment, ensure that RNAse inhibitor and the 1x Nuclei Buffer are used as recommended.

· Lysis buffer strength:

If nuclei quality is poor at short lysis times, buffer strength can be decreased for a gentler lysis

· Sample cleanup steps:

Additional cleanup steps such as washes, filtering, density gradient centrifugation, and FACS may be necessary to clean up excess debris present in the sample

See Appendix for additional optimization and troubleshooting guidance



Preparation – Buffers

Diluted Nuclei Buffer Prepare fresh, maintain at 4°C	Stock	Final	1 ml
Nuclei Buffer* (20X) DTT RNase inhibitor (check vendor-specific stock concentration)	20Χ 1000 mM 40 U/μl	1Χ 1mM 1 U/μl	50 μl 1 μl 25 μl
Nuclease-free Water	-	-	924 µl

Wash Buffer Prepare fresh, maintain at 4°C	Stock	Final	4 ml
Tris-HCl (pH 7.4)	1 M	10 mM	40 µl
NaCl	5 M	10 mM	8 µl
MgCl ₂	1 M	3 mM	12 µl
BSA	10%	1%	400 µl
Tween-20	10%	0.1%	40 µl
DTT	1000 mM	1 mM	4 μl
RNase inhibitor	40 U/µl	1 U/μl	100 μl
Nuclease-free Water	-	-	3.40 ml

Lysis Buffer Prepare fresh, maintain at 4°C	Stock	Final	2 ml
Tris-HCl (pH 7.4)	1 M	10 mM	20 μl
NaCl	5 M	10 mM	4 μl
MgCl ₂	1 M	3 mM	6 µl
Tween-20	10%	0.1%	20 μl
Nonidet P40 Substitute	10%	0.1%	20 μl
(alternatively, use IGEPAL CA-630)			
If using Sigma 74385 or i8896,			
prepare a 10% stock			
Digitonin	5%	0.01%	4 μl
(incubate at 65°C to dissolve			
precipitate before use)			
BSA	10%	1%	200 μl
DTT	1000 mM	1 mM	2 μl
RNase inhibitor 40 U/µl	40 U/µl	1 U/μl	50 µl
Nuclease-free Water	-	-	1.67 ml

RPMI + 10% FBS (maintain at 4°C, pre-warm at 37°C before use)

PBS + 0.04% BSA (maintain at 4°C)

Additional Buffers

Specific Reagents & Consumables

Vendor	Item	Part Number
10x Genomics	Nuclei Buffer* (20X) 2000	153/2000207
Thermo Fisher	Digitonin Tubes, 0.2 ml, flat cap tube**	BN2006 AB0620
	Sorvall Microtube Adapters**	76003750
Sigma-Aldrich	Trizma Hydrochloride Solution, pH 7.4 Sodium Chloride Solution, 5 M Magnesium Chloride Solution, 1M Nonidet P40 Substitute (alternatively, use IGEPAL CA-630) Sigma Protector RNase inhibitor (substitution is not recommended) DTT	T2194 59222C M1028 74385 i8896 3335402001 646563
Miltenyi Biotec	MACS BSA Stock Solution	130-091-376
Bel-Art	Flowmi Cell Strainer, 40 µm	H13680-0040
Bio-Rad	Tween 20	1662404
Corning	1X Phosphate-Buffered Saline, pH 7.4	21-040-CV

^{*}Included in the 10x Genomics Single Cell Multiome ATAC Kit A

See Appendix for DNase Treatment specific reagents & buffers

^{**}ONLY for Low Cell Input Nuclei Isolation protocol


Protocol Overview


Frozen Cells:

Thaw cells using the protocol for thawing cell lines or primary/fragile cells (step 1)

Fresh Cells:

Wash cells with PBS + 0.04% BSA, determine cell count, and proceed to Nuclei Isolation (step 2)

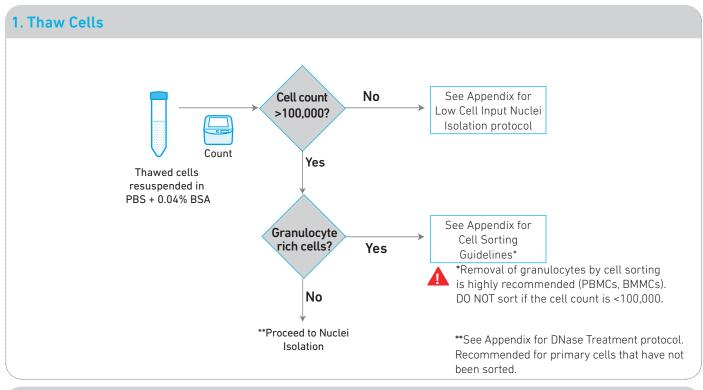

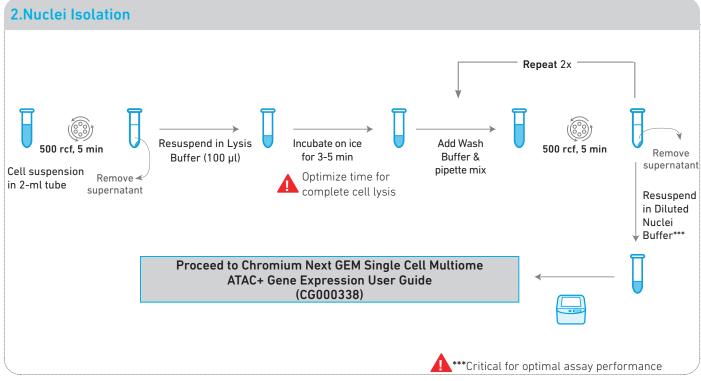

Protocol

If using fresh cells, perform 1-2 washes with PBS + 0.04% BSA, determine cell count, and proceed to Nuclei Isolation (step 2).

1. Thaw Cells (if using frozen cells)

For cell lines (used for GM12878 and 3T3 cells):

- a. Remove cryovials from storage, thaw in the water bath at 37°C for 1-2 min. Remove from the water bath when a tiny ice crystal remains.
- b. Pipette mix the cells and transfer to a 15-ml conical tube containing 10 ml pre-warmed media (RPMI + 10% FBS).
- c. Centrifuge at 300 rcf for 5 min.
- d. Remove the supernatant without disrupting the cell pellet and resuspend in 1 ml PBS + 0.04% BSA. Transfer to a 2-ml microcentrifuge tube. Rinse the 15-ml tube with 0.5 ml PBS + 0.04% BSA and transfer the rinse to the 2-ml tube containing the cells.
- e. Centrifuge cells at 300 rcf for 5 min.
- f. Remove the supernatant without disrupting the cell pellet and resuspend in 1 ml PBS + 0.04% BSA.
- g. Pass cell suspension through a 40 µm Flowmi Cell Strainer.
- h. Determine the cell concentration using a Countess II FL Automated Cell Counter (see Appendix) or a hemocytometer.
- i. Proceed to Nuclei Isolation (step 2).
 If cell count is <100,000, nuclei may be isolated using the Low Cell Input Nuclei Isolation protocol (see Appendix).

For primary cells/fragile cells (used for PBMCs):

- a. Remove cryovial from storage, thaw in the water bath at 37°C for 1-2 min. Remove from the water bath when a tiny ice crystal remains.
- b. Transfer the thawed cells to a 50-ml conical tube. Rinse the cryovial with 1 ml pre-warmed media (RPMI + 10% FBS) and add the rinse drop-wise to the 50-ml conical tube while gently shaking the tube (300 rpm at 37°C on the thermomixer).
- c. Sequentially dilute cells in the 50-ml conical tube by incremental 1:1 volume additions of media for a total of 5 times (including dilution at step b), with ~1 min wait between additions (see Appendix). Add media (RPMI + 10% FBS) at a speed of 1 ml/3-5 sec to the tube and swirl.
- d. Centrifuge at 300 rcf for 5 min.
- e. Remove most of the supernatant, leaving ~1 ml and resuspend cell pellet in this volume.
- f. Add an additional 9 ml media (1 ml/3-5 sec) to achieve a total volume of ~10 ml.
- g. Centrifuge at 300 rcf for 5 min.
- h. Remove the supernatant without disrupting the cell pellet and resuspend in 1 ml PBS + 0.04% BSA, gently pipette mix 5x.

- i. Transfer to a 2-ml microcentrifuge tube. Rinse the 50-ml tube with 0.5 ml PBS + 0.04% BSA and transfer the rinse to the 2-ml tube containing the cells. Mix by gently inverting the tube.
 - **OPTIONAL** Primary/fragile cells may have high amounts of ambient/background DNA. Treating the cells with DNase I prior to nuclei isolation can remove the ambient DNA, improving the-quality of Single Cell ATAC libraries (see Appendix for protocol). DNase treatment is not necessary if cells are being sorted prior to nuclei isolation.
- j. Centrifuge cells at 300 rcf for 5 min.
- k. Remove the supernatant without disrupting the cell pellet and resuspend in 1 ml PBS + 0.04% BSA.
- l. Pass cell suspension through a 40 μm Flowmi Cell Strainer.
- m.Determine the cell concentration using a Countess II FL Automated Cell Counter (see Appendix) or a hemocytometer.



Cell sorting for granulocyte removal before nuclei isolation is highly recommended for granulocyte-rich samples, such as PBMCs and BMMCs (see Appendix for Cell Sorting Guidelines). Removing RNase-rich granulocytes prevents degradation of mRNA and excludes highly transposed DNA from Neutrophil Extracellular traps (NETs), resulting in cleaner data.

DO NOT sort cells if the cell count is <100.000.

 Proceed to Nuclei Isolation (step 2).
 If cell count is <100,000, nuclei may be isolated using the Low Cell Input Nuclei Isolation protocol (see Appendix).

Protocol

2. Nuclei Isolation

If using fresh cells, perform 1-2 washes with PBS + 0.04% BSA and determine cell count, before proceeding to step 2a. Nuclei may be isolated from 100,000-1,000,000 cells using this protocol.

- a. Add 100,000-1,000,000 cells to a 2-ml microcentrifuge tube. Centrifuge at 300 rcf for 5 min at 4°C.
- Remove ALL the supernatant without disrupting the cell pellet.
- c. Add 100 µl chilled Lysis Buffer. Pipette mix 10x.
- d. Incubate for 3-5 min* on ice.
 *Cryopreserved PBMCs were incubated for 3 min
 *Cryopreserved cell lines were incubated for 5 min



Optimize incubation time based on cell type. Suboptimal or prolonged lysis times can both alter assay performance. Assess lysis efficacy via the Countess II FL Automated Cell Counter/microscopy (see Results).

- e. Add 1 ml chilled Wash Buffer to the lysed cells. Pipette mix 5x.
- f. Centrifuge at 500 rcf for 5 min at 4°C.
- g. Remove the supernatant without disrupting the nuclei pellet.
- h. Repeat steps e-g two more times for a total of 3 washes.
- i. Based on cell concentration step 2a and assuming ~50% nuclei loss during cell lysis, resuspend in chilled Diluted Nuclei Buffer. See Nuclei Stock Concentration Table and Example Calculation in Appendix. Maintain on ice.




The resuspension in Diluted Nuclei Buffer is critical for optimal Single Cell ATAC assay performance. The composition of the Tris-based Diluted Nuclei Buffer, including Magnesium concentration, has been optimized for the Transposition and Barcoding steps in the Single Cell Multiome ATAC + GEX protocol. Suspension of nuclei in a different buffer may not be compatible.

- j. OPTIONAL If cell debris and large clumps are observed, pass through a cell strainer. For low volume, use a 40 μm Flowmi Cell Strainer to minimize volume loss.
- k. Determine the nuclei concentration using a Countess II FL Automated Cell Counter (see Appendix) or a hemocytometer.

Proceed immediately to Chromium Next GEM Single Cell Multiome ATAC + Gene Expression User Guide (CG000338).

Results

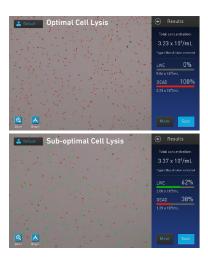

Troubleshooting

High fraction of non-viable cells in input material prior to starting nuclei isolation High fraction of non-viable cells in input material prior to starting nuclei isolation Reduce fraction of dead cells. Refer to Demonstrated Protocol Removal of Dead Cells from Single Cell Suspensions for Single Cell RNA Sequencing (Document CG000093) Sort cells using flow cytometry Gently handle cell suspensions by following best practices and reduce cell processing times High fraction of viable cells post cell lysis Difficult to count nuclei/excess debris Difficult to count for compatible cell counter or microscope Excessive debris Sample may be cleaned by extra washes/filtering/desnity centrigugation/FACS (7-AAD stain) Low nuclei recovery Use a swing-bucket rotor for centrifugation steps		
non-viable cells in input material prior to starting nuclei isolation Reduce fraction of dead cells. Refer to Demonstrated Protocol Removal of Dead Cells from Single Cell Suspensions for Single Cell RNA Sequencing (Document CG000093) Sort cells using flow cytometry Gently handle cell suspensions by following best practices and reduce cell processing times High fraction of viable cells post cell lysis Difficult to count nuclei/excess debris Difficult to count financie/excess debris Sample may be cleaned by extra washes/filtering/desnity centrigugation/FACS (7-AAD stain) Low nuclei Use a swing-bucket rotor for	Problem	Possible Solution
viable cells post cell lysis Difficult to count nuclei/excess debris Excessive debris Use a fluorescent dye (ethidium-homodimer-1) and fluorescence compatible cell counter or microscope Excessive debris Sample may be cleaned by extra washes/filtering/desnity centrigugation/FACS (7-AAD stain) Low nuclei Use a swing-bucket rotor for	non-viable cells in input material prior to starting	quality Reduce fraction of dead cells. Refer to Demonstrated Protocol Removal of Dead Cells from Single Cell Suspensions for Single Cell RNA Sequencing (Document CG000093) Sort cells using flow cytometry Gently handle cell suspensions by following best practices and reduce cell
nuclei/excess debris homodimer-1) and fluorescence compatible cell counter or microscope Excessive debris Sample may be cleaned by extra washes/ filtering/desnity centrigugation/FACS (7-AAD stain) Low nuclei Use a swing-bucket rotor for	viable cells post	
filtering/desnity centrigugation/FACS (7-AAD stain) Low nuclei Use a swing-bucket rotor for	nuclei/excess	homodimer-1) and fluorescence
9	Excessive debris	filtering/desnity centrigugation/FACS
		•

Trypan Blue Precipitate in the Countess II Slide



DO NOT use nuclei resuspended in 20X Nuclei Buffer. Repeat nuclei isolation and resuspend in Diluted Nuclei Buffer (1X).


Nuclei Quality - Representative Images (Panel A: recommended quality)



Appendix

Illustrative Overview of Incremental 1:1 Volume Additions

Incrementally add media volumes at the indicated points & swirl

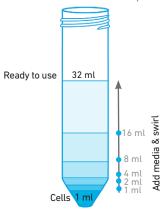

Nuclei Stock Concentration Table

Based on the Targeted Nuclei Recovery, prepare the nuclei suspension in Diluted Nuclei Buffer to achieve the corresponding Nuclei Stock concentrations.

Targeted Nuclei Recovery	Nuclei Stock Concentration (nuclei/µl)
500	160-400
1,000	320-810
2,000	650-1,610
3,000	970-2,420
4,000	1,290-3,230
5,000	1,610-4,030
6,000	1,940-4,840
7,000	2,260-5,650
8,000	2,580-6,450
9,000	2,900-7,260
10,000	3,230-8,060

Appendix

Nuclei Counting and Viability

Countess II FL Automated Cell Counter is recommended for determining nuclei concentrations. The optimal range of cell concentration for Cell Counter is 1,000-10,000 cells/µl. Refer to manufacturer's instructions for details on operations.

- a. Vortex 0.4% trypan blue stain, centrifuge briefly and aliquot 10 μl per tube.
- b. Pipette mix the nuclei suspension. Immediately add 10 μl nuclei suspension to 10 μl aliquot of 0.4% trypan blue stain. Gently pipette mix 10x.
- c. Transfer 10 μl trypan blue stained nuclei to a Countess II Cell Counting Slide chamber.
- d. Insert the slide into the Countess II FL Cell Counter, and determine the nuclei concentration and viability. <5% of input cells should be viable. Optimize focusing and light exposure.

Low Cell Input Nuclei Isolation

Nuclei may be isolated from 2,000-100,000 cells using this protocol. If cell count is <40,000, centrifuge cell suspension at 300 rcf for 5 min at 4°C and resuspend the cell pellet in 50 μ l PBS + 0.04% BSA. Transfer 50 μ l cell suspension to a 0.2-ml tube. Proceed directly to step c.

- a. Centrifuge cell suspension at 300 rcf for 5 min at 4°C. Remove supernatant and resuspend pellet in PBS + 0.04% BSA for 1,000 cells/ μ l cell suspension.
- **b.** Add 2,000–40,000 cells to a 0.2-ml tube in a total volume of **50 \mul** PBS + 0.04% BSA.

Approximately 25% of the cell input is expected to be recovered during Chromium Single Cell Multiome ATAC + GEX sequencing. Always determine nuclei counts after nuclei isolation.

Cell Input	Expected Nuclei Recovery (after cell lysis)	Expected Nuclei Recovery (ATAC + GEX)
40,000	16,000	10,000
20,000	8,000	5,000
10,000	4,000	2,500
4,000	1,600	1,000
2,000	800	500

- c. Centrifuge at 300 rcf for 5 min at 4°C.
- d. Remove 45 µl supernatant without touching the bottom of the tube to avoid dislodging the cell pellet.
- e. Add 45 µl chilled Lysis Buffer. Gently pipette mix 3x.
- f. Incubate for 3-5 min* on ice.

- *Cryopreserved PBMCs were incubated for 3 min
- *Cryopreserved cell lines were incubated for 5 min



Optimize incubation time based on cell type. Suboptimal or prolonged lysis times can both alter assay performance. Assess lysis efficacy via the Countess II FL Automated Cell Counter/microscopy. See Results for optimal cell lysis.

- g. Add 50 µl chilled Wash Buffer to the tube. DO NOT mix.
- h. Centrifuge at 500 rcf for 5 min at 4°C.
- Remove 95 μl supernatant without disrupting the nuclei pellet.
- j. Add 45 μl chilled Diluted Nuclei Buffer to the pellet. DO NOT mix.
- k. Centrifuge at 500 rcf for 5 min at 4°C.
- Remove the supernatant without touching the bottom of the tube to avoid dislodging the nuclei pellet.





The supernatant may be removed in two steps, first with a 100- μ l pipette (set to 40 μ l), followed by removal with a 10- μ l pipette (set to 10 μ l).

m. Resuspend the nuclei pellet in $7 \mu l$ chilled Diluted Nuclei Buffer (pellet may not be visible).




The use of the Diluted Nuclei Buffer for nuclei suspension is critical for optimal assay performance. The composition of the Tris-based Diluted Nuclei Buffer, including Magnesium concentration, has been optimized for the Transposition and Barcoding steps in the Single Cell Multiome ATAC + GEX protocol. Suspension of nuclei in a different buffer may not be compatible with these protocol steps.


- n. Use 2 μl nuclei suspension mixed with 8 μl Diluted Nuclei Buffer and 10 μl Trypan Blue to determine the cell concentration by a Countess II FL Automated Cell Counter (see Appendix) or a hemocytometer. A final nuclei concentration of 30 nuclei/μl is needed for Targeted Nuclei Recovery of 500.
- Proceed immediately to Chromium Next GEM Single Cell Multiome ATAC + Gene Expression User Guide (CG000338).

Appendix

Cell Sorting Guidelines

- Removal of granulocytes by cell sorting is highly recommended before nuclei isolation from granulocyte-rich samples, such as PBMCs and BMMCs.
- Cell sorting is not recommended if cell count is <100,000.
- After thawing and counting cells (step 1m), cells can be sorted using a 100 µm nozzle (BD FACSMelody or comparable).
 No stain is needed for cell sorting. Granulocytes can be identified using side scatter.
- The first collection of sorted cells is for singlets as shown in the representative plot A.
- The second collection of sorted cells separates lymphocytes and monocytes (green cells) from granulocytes (pink cells) as shown in the representative plot B.

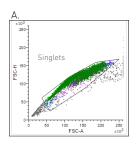

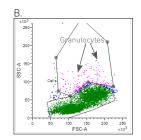

- The lymphocytes and monocytes fraction should be collected and used for isolating nuclei.
- Collect the sorted cells in a 5-ml FACS tube containing 500 µl PBS + 0.04% BSA.
- Centrifuge the collected cells at 300 rcf for 5 min at 4°C.
- Remove the supernatant without disrupting the cell pellet and resuspend in 500 µl PBS + 0.04% BSA.
- Determine the cell concentration using a Countess II FL Automated Cell Counter (see Appendix) or a hemocytometer.
- Proceed directly to Nuclei Isolation (step 2).

References

 Chromium Next GEM Single Cell Multiome ATAC + Gene Expression User Guide (CG000338)

DNase Treatment

Recommended for primary cells prior to nuclei isolation. DNase treatment is not necessary if cells are being sorted prior to nuclei isolation.

Specific Reagents

DNase I, RNase-free (includes 10x Reaction Buffer with MgCl₂) from ThermoFisher Scientific, Part Number-EN0521

Preparation - Buffers

10X TBS	Stock	Final	5 ml
Tris-HCl (pH 7.4)	1 M	200 mM	1 ml
NaCl	5 M	1.5 M	1.5 ml
Nuclease-free Water	-	3 mM	2.5 ml
DNase Solution Prepare fresh, maintain at 4°C	Stock	Final	1 ml
TBS	10X	1X	100 µl
10X Reaction Buffer with MgCl ₂	10X	1X	100 µl
DNase I	1 U/μl	0.1 U/µl	100 µl
Nuclease-free Water	-	-	700 µl

Primary cells/fragile cells may have high amounts of ambient/background DNA. Treating the cells with DNase I prior to nuclei isolation can reduce the ambient DNA, which may improve the quality of libraries.

a. Centrifuge the cells in a 2-ml microcentrifuge tube at 300 rcf for 10 min at 4°C.



Using a 2-ml microcentrifuge tube and centrifuging for a longer time (10 min) is critical in maintaining an equal proportion of all cell types.

- b. Remove supernatant without disrupting the pellet and resuspend the pellet in 300 µl DNase Solution.
- c. Pipette mix 5x and incubate on ice for 5 min.
- d. Add 1 ml PBS + 0.04% BSA.
- e. Centrifuge cells at 300 rcf for 10 min at 4°C.
- f. Remove supernatant without disrupting the pellet and resuspend the pellet in 1 ml PBS + 0.04% BSA.
- g. Repeat steps d-e for a total of 2 washes.
- h. Pass cell suspension through a 40 μm Flowmi Cell Strainer.
- Determine the nuclei concentration using a Countess II FL Automated Cell Counter (see Appendix) or a hemocytometer.
- j. Proceed directly to Nuclei Isolation (step 2).

© 2022 10x Genomics, Inc. (10x Genomics). All rights reserved. Duplication and/or reproduction of all or any portion of this document without the express written consent of 10x Genomics, is strictly forbidden. Nothing contained herein shall constitute any warranty, express or implied, as to the performance of any products described herein. Any and all warranties applicable to any products are set forth in the applicable terms and conditions of sale accompanying the purchase of such product. 10x Genomics provides no warranty and hereby disclaims any and all warranties as to the use of any third-party products or protocols described herein. The use of products described herein is subject to certain restrictions as set forth in the applicable terms and conditions of sale accompanying the purchase of such product. A non-exhaustive list of 10x Genomics' marks, many of which are registered in the United States and other countries can be viewed at: www.10xgenomics. com/trademarks. 10x Genomics may refer to the products or services offered by other companies by their brand name or company name solely for clarity, and does not claim any rights in those third-party marks or names. 10x Genomics products may be covered by one or more of the patents as indicated at: www.10xgenomics. com/patents. The use of products described herein is subject to 10x Genomics Terms and Conditions of Sale, available at www.10xgenomics.com/legal-notices, or such other terms that have been agreed to in writing between 10x Genomics and user. All products and services described herein are intended FOR RESEARCH USE ONLY and NOT FOR USE IN DIAGNOSTIC PROCEDURES.

The use of 10x Genomics products in practicing the methods set forth herein has not been validated by 10x Genomics, and such non-validated use is NOT COVERED BY 10X GENOMICS STANDARD WARRANTY, AND 10X GENOMICS HEREBY DISCLAIMS ANY AND ALL WARRANTIES FOR SUCH USE. Nothing in this document should be construed as altering, waiving or amending in any manner 10x Genomics terms and conditions of sale for the Chromium Controller or the Chromium Ch

Contact:

support@10xgenomics.com 10x Genomics 6230 Stoneridge Mall Road Pleasanton, CA 94588 USA

