

# Isolation of Cells from FFPE Tissue Sections for Chromium Fixed RNA Profiling

## Introduction

This protocol outlines how to isolate cells from formaldehyde fixed & paraffin embedded (FFPE) Tissue Sections for use with Chromium Fixed RNA Profiling workflow.

The protocol includes instructions for deparaffinizing and rehydrating the tissue sections, followed by preparation and addition of a Dissociation Enzyme Mix to the sections. The tissue sections can then be dissociated using either a pellet pestle or a gentleMACS Octo Dissociator. Once dissociated, cells are resuspended in the Tissue Resuspension Buffer or Quenching Buffer prior to counting.

The minimum recommended cell inputs for Chromium Fixed RNA Profiling hybridization following this protocol are:

- 400,000 cells for singleplex hybridization (Probe Hybridization step 1.1d in User Guide CG000477\*)



*When following User Guide CG000477, for steps 2.1g-2.1j, perform one extra 0.5 ml wash for a total of three 0.5 ml post-hybridization washes.*

- 100,000 cells per Probe Barcode for multiplex hybridization (Probe Hybridization step 1.1d in User Guide CG000527)

Performing a pilot multiplex assay run (4 rxns with 4 Probe Barcodes) is recommended prior to committing to larger studies.

*\*Using cells derived from FFPE tissue sections with Feature Barcode technology for cell surface protein detection is not supported.*

For optimal assay results, loading 10,000 cells per Probe Barcode is recommended when performing the singleplex assay or when multiplexing 4 samples. For multiplexing 16 samples, loading 8,000 cells per Probe Barcode is recommended (128,000 is the maximum cells recovery per GEM reaction).

The protocol was demonstrated using FFPE tissue blocks ranging 1-10 years in age. Each block may yield different amounts of material and data quality, depending on age, tissue type, pre-fixation tissue quality, tissue density, size/area of tissue in the scrolls, and other factors.

## Additional Guidance

This protocol was demonstrated using two or more 25 µm or 50 µm tissue sections from a diverse spectrum of tissue types. For human tissue, it is recommended to use two or more 25 µm sections and for mouse tissue, it is recommended to use two or more 50 µm sections. Some tissue types require more than two sections to yield enough cells. FFPE block properties will also impact yields. Note that poor quality blocks will likely yield data that cannot be interpreted accurately or salvaged. See [Appendix](#) for cell yields derived from tissue sections of indicated tissue types. In rare instances, a single 25 µm or 50 µm section may also yield adequate cells.

**Tissue and cells may carry potentially hazardous pathogens. Follow material supplier recommendations and local laboratory procedures and regulations for the safe handling, storage and disposal of biological materials.**

## Specific Reagents & Consumables

| Vendor                                                     | Item                                                        | Part Number |
|------------------------------------------------------------|-------------------------------------------------------------|-------------|
| <b>For Tissue Section Transfer &amp; Deparaffinization</b> |                                                             |             |
| Millipore Sigma                                            | Xylene, Reagent Grade                                       | 214736      |
|                                                            | Ethyl Alcohol, 200 Proof, anhydrous                         | E7023       |
| VWR                                                        | Ethanol absolute ≥99.5%, TechniSolv, pure (for Europe)      | 83813.360DP |
| Corning                                                    | Phosphate-Buffered Saline, 1X without Calcium and Magnesium | 21-040-CV   |
|                                                            | 15 ml PP Centrifuge Tubes                                   | 430791      |
| Thermo Fisher                                              | Nuclease-free Water (not DEPC-Treated)                      | AM9937      |
| <b>If using gentleMACS Octo Dissociator</b>                |                                                             |             |
| Miltenyi Biotec                                            | gentleMACS C Tubes                                          | 130-093-237 |

### Additional Materials

Disposable plastic Pasteur pipette  
Nuclease-free Water  
Water bath  
Centrifuge

### For Tissue Section Dissociation

|                 |             |            |
|-----------------|-------------|------------|
| Millipore Sigma | Liberase TH | 5401151001 |
| Corning         | RPMI        | 10-040-CV  |

### If using pellet pestle:

|                   |                                                                                           |            |
|-------------------|-------------------------------------------------------------------------------------------|------------|
| Fisher Scientific | RNase-Free Disposable Pellet Pestles                                                      | 12-141-364 |
| Fisher Scientific | BD Luer-Lok PrecisionGlide Disposable Syringes with Detachable Needles<br><i>OPTIONAL</i> | 14-823-37  |

### If using gentleMACS Octo Dissociator:

|                 |                                          |             |
|-----------------|------------------------------------------|-------------|
| Miltenyi Biotec | gentleMACS Octo Dissociator with Heaters | 130-096-427 |
|-----------------|------------------------------------------|-------------|

### If using Tissue Resuspension Buffer:

Refer to SDS and follow local and institutional guidelines for proper handling and disposal of all chemicals.

|                                                           |                                                             |            |
|-----------------------------------------------------------|-------------------------------------------------------------|------------|
| Millipore Sigma                                           | Protector RNase Inhibitor                                   | 3335399001 |
|                                                           | Albumin, Bovine Serum*, 10% Aqueous Solution, Nuclease-Free | 126615     |
| Thermo Fisher Scientific                                  | UltraPure Bovine Serum Albumin* (BSA, 50 mg/ml)             | AM2616     |
| *Choose either Millipore or Thermo Fisher Scientific BSA. |                                                             |            |

|     |                                      |            |
|-----|--------------------------------------|------------|
| VWR | Tris Buffer, 1M sterile sol., pH 8.0 | E199-100ML |
|-----|--------------------------------------|------------|

### If using Quenching Buffer:

|                                                                                                                                                                                  |                       |         |
|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-----------------------|---------|
| 10x Genomics                                                                                                                                                                     | Conc. Quench Buffer** | 2000516 |
| **Included in the 10x Genomics Chromium Next GEM Single Cell Fixed RNA Sample Preparation Kit, 16 rxns (PN-1000414); also includes Enhancer (PN-2000482) if storing fixed cells. |                       |         |

### For Sample Filtration

|                 |                                                                                  |                            |
|-----------------|----------------------------------------------------------------------------------|----------------------------|
| Sysmex          | Sterile Single-Pack CellTrics Filters*** (use 30 µm)                             | 04-004-2326                |
| Miltenyi Biotec | MACS SmartStrainers*** (30 µm)<br><i>OR</i><br>Pre-Separation Filters*** (30 µm) | 130-098-458<br>130-041-407 |

\*\*\*Choose either Sysmex or Miltenyi Biotec filter.

### For Cell Counting

|                          |                                                                               |              |
|--------------------------|-------------------------------------------------------------------------------|--------------|
| Nexcelom Biosciences     | ViaStain PI Staining Solution                                                 | CS2-0109-5mL |
|                          | ViaStain AOPI Staining Solution<br><i>Alternative to PI Staining Solution</i> | CS2-0106-5mL |
|                          | Cellaca MX High-throughput* Automated Cell Counter                            | MX-112-0127  |
| Thermo Fisher Scientific | Countess II FL Automated Cell Counter†                                        | AMAQAF1000   |
|                          | Countess II FL Automated Cell Counting Chamber Slides                         | C10228       |
|                          | Ethidium Homodimer-1                                                          | E1169        |

\*Choose either Countess or Cellaca.

This list may not include some standard laboratory equipment.

## Preparation - Buffers

All buffers should be prepared fresh.

Prepare Dissociation Enzyme Mix, incubate at **37°C** for **10 min** before proceeding with dissociation.

### For pestle-based protocol

| Dissociation Enzyme Mix  | Stock | Final | 1 rxn (μl)  | 4 rxn + 10% (μl) |
|--------------------------|-------|-------|-------------|------------------|
| Liberase TH (mg/ml)      | 5     | 1     | 210         | 924              |
| RPMI                     | -     | -     | 840         | 3696             |
| <b>Total Volume (μl)</b> |       |       | <b>1050</b> | <b>4620</b>      |

### For gentleMACS Octo Dissociator protocol

| Dissociation Enzyme Mix  | Stock | Final | 1 rxn (μl)  | 4 rxn + 10% (μl) |
|--------------------------|-------|-------|-------------|------------------|
| Liberase TH (mg/ml)      | 5     | 1     | 420         | 1848             |
| RPMI                     | -     | -     | 1680        | 7392             |
| <b>Total Volume (μl)</b> |       |       | <b>2100</b> | <b>9240</b>      |

Prepare either Tissue Resuspension Buffer or Quenching Buffer. They can be used interchangeably.

| Tissue Resuspension Buffer (Maintain at 4°C) | Stock | Final  | 1 rxn (μl) | 4 rxn + 10% (μl) |
|----------------------------------------------|-------|--------|------------|------------------|
| PBS                                          | 1X    | 0.496X | 248        | 1091.2           |
| Tris buffer (pH 8.0; mM)                     | 1000  | 50     | 25         | 110              |
| BSA (RNase free)                             | 10%   | 0.02%  | 1          | 4.4              |
| RNase Inhibitor (U/μl)                       | 40    | 0.24   | 3          | 13.2             |
| Nuclease-free Water                          | -     | -      | 223        | 981.2            |
| <b>Total Volume (μl)</b>                     |       |        | <b>500</b> | <b>2200</b>      |


OR

| Quenching Buffer (Maintain at 4°C)             | Stock | Final | 1 rxn (μl) | 4 rxn + 10% (μl) |
|------------------------------------------------|-------|-------|------------|------------------|
| Nuclease-free Water                            | -     | -     | 437.5      | 1925             |
| Conc. Quench Buffer* (10x Genomics PN 2000516) | 8X    | 1X    | 62.5       | 275              |
| <b>Total Volume (μl)</b>                       |       |       | <b>500</b> | <b>2200</b>      |

\*Included in the 10x Genomics Chromium Next GEM Single Cell Fixed RNA Sample Preparation Kit, 16 rxns (PN-1000414).

**Ethanol:** Prepare fresh 70% and 50% Ethanol (1 ml each/sample).

## Step Selection Overview

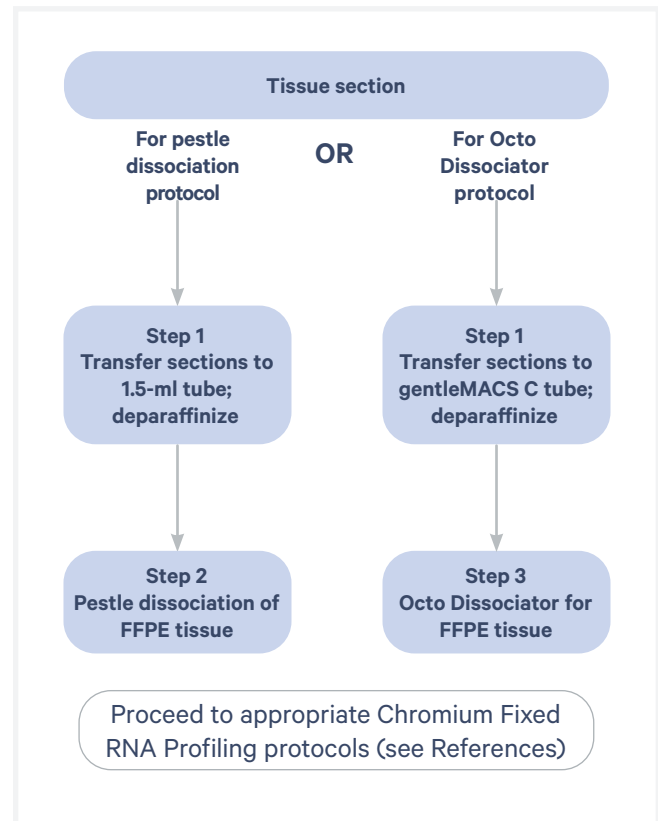


Refer to SDS and follow local and institutional guidelines for proper handling and disposal of all chemicals.


## Protocol Overview: Pestle for Isolation of Cells from FFPE Tissue Sections

### 1. Tissue Section Transfer & Deparaffinization

The scroll illustration is an abstract representation; scrolls will be curled in the tube during execution

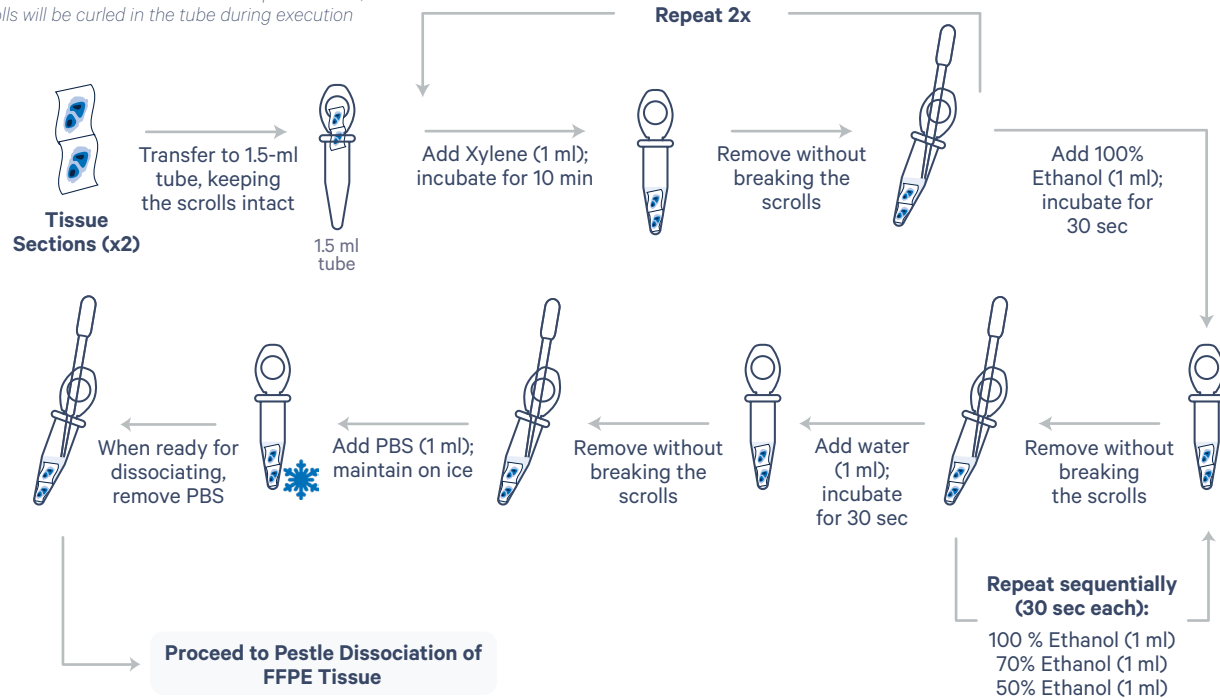


### 2. Pestle Dissociation of FFPE Tissue

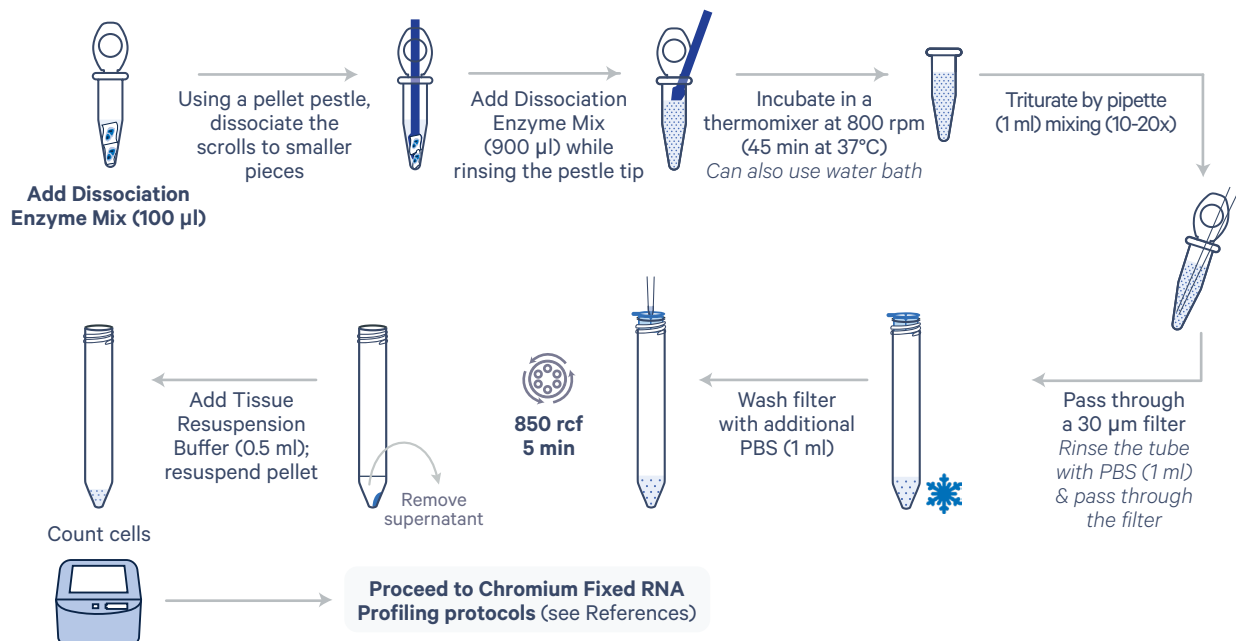


## Protocol Overview: gentleMACS Octo Dissociator for Isolation of Cells from FFPE Tissue Sections

### 1. Tissue Section Transfer & Deparaffinization

The scroll illustration is an abstract representation; scrolls will be curled in the tube during execution

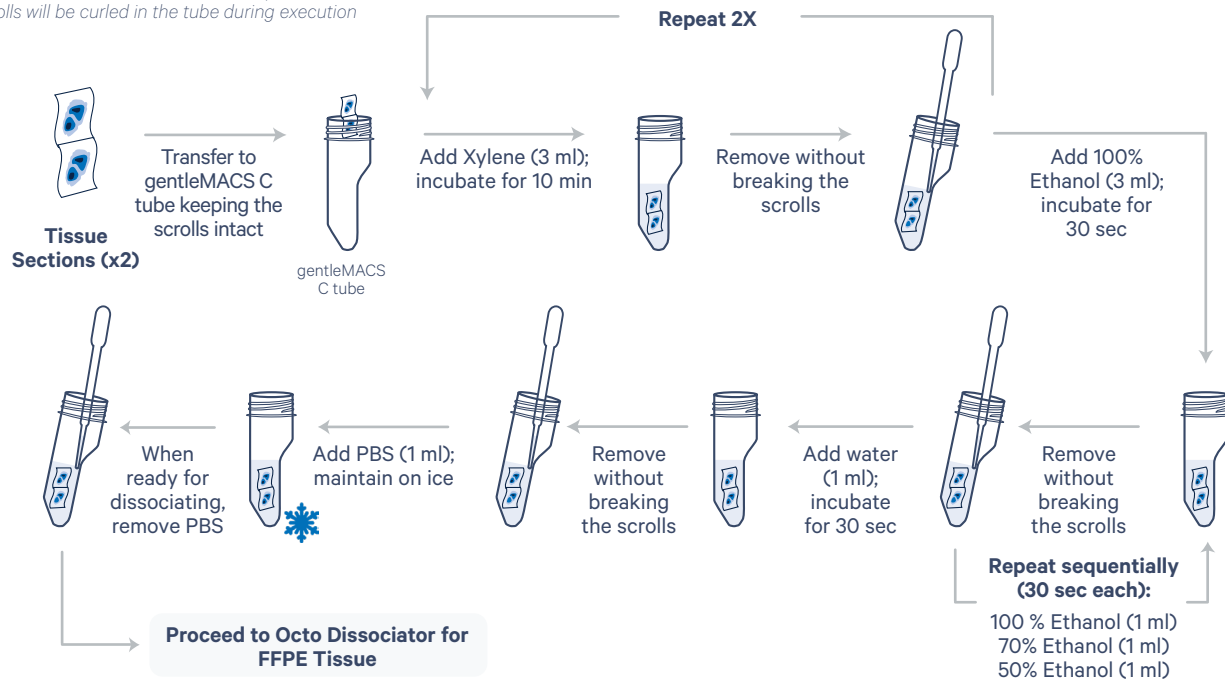

### 3. gentleMACS Octo Dissociator for FFPE Tissue

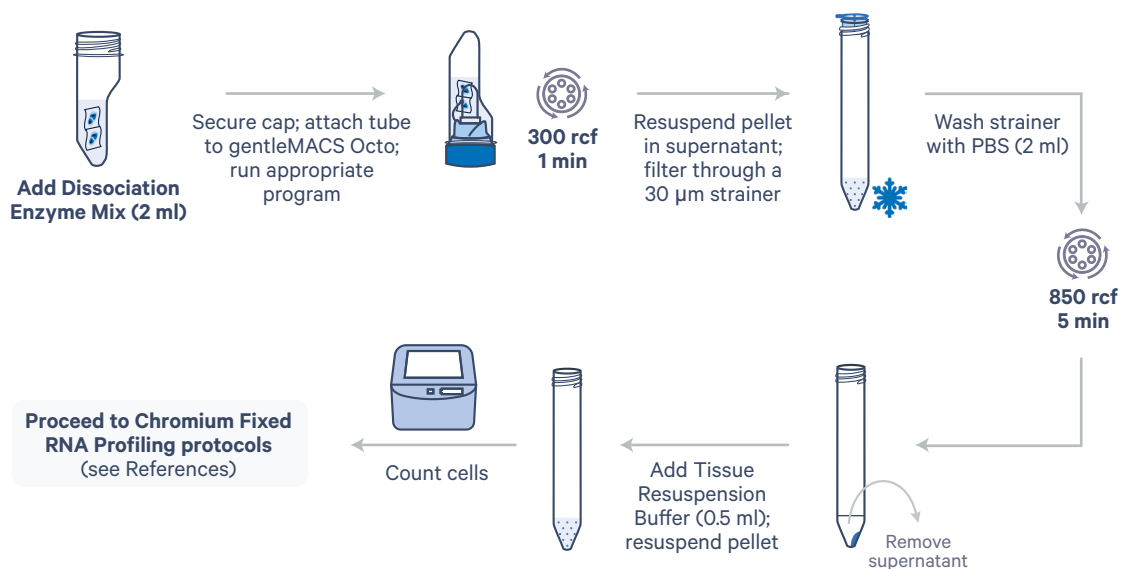

## Isolating Cells from FFPE Tissue Sections

This protocol was demonstrated using 25 µm or 50 µm FFPE tissue sections of various tissue types (see [Appendix](#) for details; maximum of six 50 µm sections were dissociated). The cells derived from the FFPE sections are compatible with the Chromium Fixed RNA Profiling assay. The recommended starting cell count for the Chromium Fixed RNA Profiling for the singleplex protocol (CG000477) is 400,000 cells and for the multiplex protocol (CG000527) is 100,000 cells.

### 1. Tissue Section Transfer & Deparaffinization

*Before starting, carefully inspect the tissue block to gauge the extent of dehydration. Face the tissue block and place for sufficient time in an ice water bath (block facing the water surface) to ensure proper rehydration of the tissue block.*

*See representative images of scrolls derived from well rehydrated and less rehydrated blocks in the [Troubleshooting](#) section.*

*Prepare the first section from the rehydrated tissue block by setting microtome to 5 µm. Discard the first section and set the microtome for 25 µm or 50 µm sections.*



*The volumes added in steps 1b, 1e, and 1g are based on the type of tube being used for the protocol.*

- a.** After discarding the first section, prepare up to two 25 or two 50 µm sections from the rehydrated tissue block and transfer into a 1.5 ml tube (or to a gentleMACS C Tube) while keeping the scrolls intact. The scrolls in the tube can be stored at 4°C for up to **one week**.



*Ensure that the scrolls stay intact during transfer and through steps 1a-1m. Add reagents along the side of the tube and not directly onto the scrolls. Intact sections enable easier liquid aspiration and minimize tissue loss. If scrolls break, see guidance in [Troubleshooting](#) section. Proceed with the protocol and based on cell yields process more sections if needed.*




*Always use fresh/freshly prepared reagents for the deparaffinization steps. Ensure the scrolls are fully submerged in reagents in the following steps.*

- b.** Add **1 ml** xylene to the 1.5-ml tube (**3 ml** to gentleMACS C Tube) and incubate for **10 min** at **room temperature**.
- c.** Using a Pasteur pipette, remove the liquid from the tube without breaking the scrolls.
- d.** Repeat **steps 1b-1c** twice.
- e.** Add **1 ml** 100% ethanol to the 1.5 ml tube (**3 ml** to gentleMACS C Tube) and incubate for **30 sec** at **room temperature**.
- f.** Using a Pasteur pipette remove the liquid from the tube without breaking the scrolls.
- g.** Add **1 ml** 100% ethanol (to both tube types) and incubate for **30 sec** at **room temperature**.
- h.** Using a Pasteur pipette remove the liquid from the tube without breaking the scrolls.
- i.** Add **1 ml** 70% ethanol (to both tube types) and incubate for **30 sec** at **room temperature**.
- j.** Using a Pasteur pipette remove the liquid from the tube without breaking the scrolls.
- k.** Add **1 ml** 50% ethanol (to both tube types) and incubate for **30 sec** at **room temperature**.
- l.** Using a Pasteur pipette remove the liquid from the tube without breaking the scrolls.
- m.** Add **1 ml** nuclease-free water to the tube and incubate for **30 sec** at **room temperature**.
- n.** Using a Pasteur pipette remove the liquid from the tube without breaking the scrolls.
- o.** Add **1 ml** PBS and maintain on ice.
- p.** Proceed to either to **Step 2 Pestle Dissociation of FFPE Tissue** (sections are in 1.5 ml tube) or to **Step 3 gentleMACS Octo Dissociator for FFPE Tissue** (sections are in gentleMACS C Tube).


## 2. Pestle Dissociation of FFPE Tissue

Refer to [Preparation - Buffers](#) section to prepare Dissociation Enzyme Mix. Incubate at **37°C** for **10 min** before proceeding with dissociation.

- a. Remove the PBS from the tube without breaking the scrolls.
- b. Add **100 µl** Dissociation Enzyme Mix to the tube.

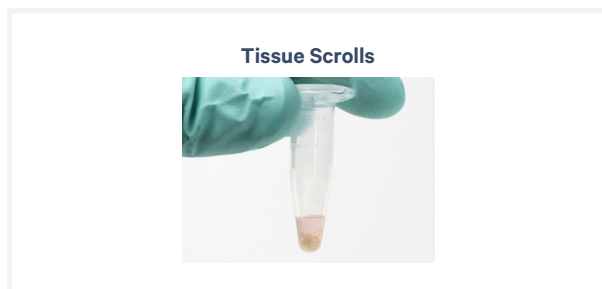

- c. Using a 1.5 ml pellet pestle, dissociate the tissue scrolls breaking them to smaller pieces.

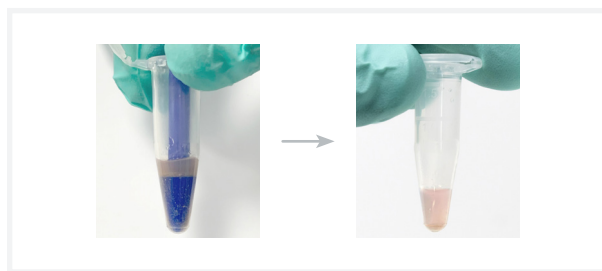

**TIPS** Grasp the pestle between the thumb and fingers and rotate the pestle (clockwise and counterclockwise 10-20X) inside the tube and go up and down with the scrolls trapped between the wall of the tube and the pestle for dissociation or until the scroll is broken into similarly sized smaller pieces.  
The pellet pestle can also be attached to a compatible cordless motor and operated inside the tube until the scroll is broken into similarly sized smaller pieces.

- d. Add **900 µl** Dissociation Enzyme Mix to the tube while rinsing the pestle tip into the tube to collect any additional tissue pieces sticking to the pestle. Pipette mix.
- e. Incubate for **45 min** at **37°C** in a thermomixer at **800 rpm**. Alternatively, incubate in a **37°C** water bath, mixing by inversion every **15 min** during the incubation.

### TIPS

Ensure no tissue pieces are stuck to the tube cap. Centrifuge at 300 rcf for ~20 sec to remove pieces stuck on the cap.

- f. Using a 1,000 µl pipette, triturate the tissue pieces in the tube by pipetting 10-20x.



### TIPS

The progress of the dissociation can be monitored by taking a 10 µl aliquot and counting. If the cell concentration is lower than recommended for hybridization, additional dissociation maybe required; recommend performing step g.

- g. **OPTIONAL:** Aspirate and push the tissue pieces and solution through a 23G needle 5x to improve cell recovery.
- h. Pass the suspension through a Pre-Separation Filter (30 µm) placed on a 5- or 15-ml tube placed on ice.
- i. Rinse the original tube (step 2f) with **1 ml** chilled PBS and rinse the 30 µm filter with additional **1 ml** chilled PBS to minimize cell loss. Collect the filtrate in the same tube as step 2h.
- j. Centrifuge the cell suspension at **850 rcf** at **4°C** for **5 min**.

### TIPS

Use of swinging-bucket rotor is recommended for higher cell recovery.

- k. Remove the supernatant without disturbing the pellet.

- l.** Resuspend the pellet in **0.5 ml** chilled Tissue Resuspension Buffer or Quenching Buffer, pipette mix 5x, and maintain on ice.
- m.** Determine cell concentration of the fixed sample using an Automated Cell Counter (Countess II/Cellaca MX) or hemocytometer. See [Appendix](#) for Fixed Cell Counting.

**TIPS**

*For accurate cell counting, it is strongly recommended that the cell suspension be stained with a fluorescent nucleic acid dye such as Ethidium Homodimer-1 or PI Staining Solution and counted using an automated fluorescent cell counter.*

- n.** Proceed **immediately** to appropriate Chromium Fixed RNA Profiling protocols – Probe Hybridization step 1.1d (see References) or store according to Fixed Sample Storage Guidance provided in [Appendix](#).

### 3. gentleMACS Octo Dissociator for FFPE Tissue

Refer to [Preparation - Buffers](#) section to prepare Dissociation Enzyme Mix. Incubate at **37°C** for **10 min** before proceeding with dissociation.

- a. Remove the PBS from the gentleMACS C Tube without breaking the scrolls.
- b. Add **2 ml** Dissociation Enzyme Mix to the gentleMACS C Tube and close securely.

**TIPS** Ensure the scrolls are fully submerged in reagents in the following steps.

- c. Place the tube in on the gentleMACS Octo Dissociator, apply Heating units and run the gentleMACS Program 37C\_FFPE\_1. Run time **~48 min**. Ensure that the Octo Dissociator blades are moving before walking away.

#### gentleMACS Program 37C\_FFPE\_1

|    |                      |
|----|----------------------|
| 1  | temp ON              |
| 2  | spin - 20 rpm, 5' 0" |
| 3  | loop 3X              |
| 4  | spin 20 rpm, 14' 0"  |
| 5  | spin 1700 rpm, 7"    |
| 6  | spin 1700 rpm, 1"    |
| 7  | spin -1700 rpm, 2"   |
| 8  | spin 1700 rpm, 1"    |
| 9  | spin 1700 rpm, 4"    |
| 10 | end loop             |
| 11 | end                  |

- d. At the end of the run, detach the tube from the gentleMACS Octo Dissociator and visually inspect to ensure the scrolls have been dissociated. If scrolls are not dissociated, see guidance in [Troubleshooting](#) section.
- e. Centrifuge at **~300 rcf** for **1 min** and resuspend the cell pellet in the supernatant.

- f. Pass the suspension through a Pre-Separation Filter (30 µm) placed on a 15-ml tube placed on ice.
- g. Rinse the original gentleMACS tube (step 3d) with **2 ml** chilled PBS and use that rinse for an additional wash of the **30 µm** filter to minimize cell loss. Collect the filtrate in the same tube as step 3f.

**TIPS** Use of swinging-bucket rotor is recommended for higher cell recovery.

- h. Centrifuge the cell suspension at **850 rcf** at **4°C** for **5 min**.
- i. Remove the supernatant without disturbing the pellet.
- j. Resuspend the pellet in **0.5 ml** chilled Tissue Resuspension Buffer or Quenching Buffer, pipette mix 5x, and maintain on ice.
- k. Determine cell concentration of the fixed sample using an Automated Cell Counter (Countess II/Cellaca MX) or hemocytometer. See [Appendix](#) for Fixed Cell Counting.

**TIPS** For accurate cell counting, it is strongly recommended that the cell suspension be stained with a fluorescent nucleic acid dye such as Ethidium Homodimer-1 or PI Staining Solution and counted using an automated fluorescent cell counter.

- l. Proceed **immediately** to appropriate Chromium Fixed RNA Profiling protocols – Probe Hybridization step 1.1d (see References) or store according to Fixed Sample Storage Guidance provided in [Appendix](#).

## Appendix

Ethidium Homodimer-1 and PI stained cells can be counted using either Countess II FL Automated Cell Counter or Cellaca Counter.

### Counting Using Ethidium Homodimer-1

This protocol provides instructions for counting samples using Ethidium Homodimer-1 and the Countess II FL Automated Cell Counter (with RFP light cube) to enable accurate quantification even in the presence of sub-cellular debris. The optimal cell concentration for the Countess is 1,000-4,000 cells/μl. Refer to manufacturer's instructions for details on operations.

- Vortex Ethidium Homodimer-1, centrifuge briefly, and dilute the concentrated stock as per manufacturer's instructions (~1:100 dilution).
- Aliquot **10 μl** diluted Ethidium Homodimer-1 in a tube.
- Gently mix the sample. Immediately add **10 μl** sample to **10 μl** diluted Ethidium Homodimer-1. Gently pipette mix 10x.
- Transfer **10 μl** sample to a Countess II Cell Counting Slide chamber.
- Insert the slide into the Countess II FL Cell Counter. Image the sample using the RFP setting for fluorescent illumination and filtering. Optimize focus and exposure settings.
- Confirm the absence of large clumps using the bright-field mode. Make sure the cell counter is circling RFP positive cells. Note the RFP-positive concentration. Multiply by dilution factor 2 to determine cell concentration.

Samples stained with Ethidium Homodimer-1 can also be counted using Cellaca Counter. Refer to manufacturer's instructions for details.

### Counting using PI Staining Solution

This protocol provides instructions for counting sample using PI staining solution and the Cellaca Counter to enable accurate quantification even in the presence of sub-cellular debris. The optimal cell concentration for the Cellaca Counter is 100-10,000 cells/μl. Refer to manufacturer's instructions for details on operations.


- Add **25 μl** PI Staining Solution into Mixing Row of Cellaca plate
- Gently mix the sample. If the sample is too concentrated, a 1:1 dilution in PBS can also be prepared. For example, add 15 μl fixed cell suspension to 15 μl PBS.
- Add **25 μl** sample to Mixing Row of plate containing PI Staining Solution. Gently pipette mix 8x.
- Transfer stained sample to Loading Row of Cellaca plate.
- For counting fixed samples, only use the PI (Propidium Iodide) channel. Refer to manufacturer's instructions for details.

Samples stained with PI staining solution can also be counted using Countess II FL Automated Cell Counter. See manufacturer's instructions for details.


When counting cells after dissociation, the suspension may include some debris. However, during the Chromium Fixed RNA Profiling workflow, the post-hybridization wash step will reduce the debris.

Representative images after dissociation and after post-hybridization wash (ready for loading on to the chip for GEM generation) are shown in Figure 1A and 1B respectively.

**Figure 1A: Representative Cells after Dissociation**

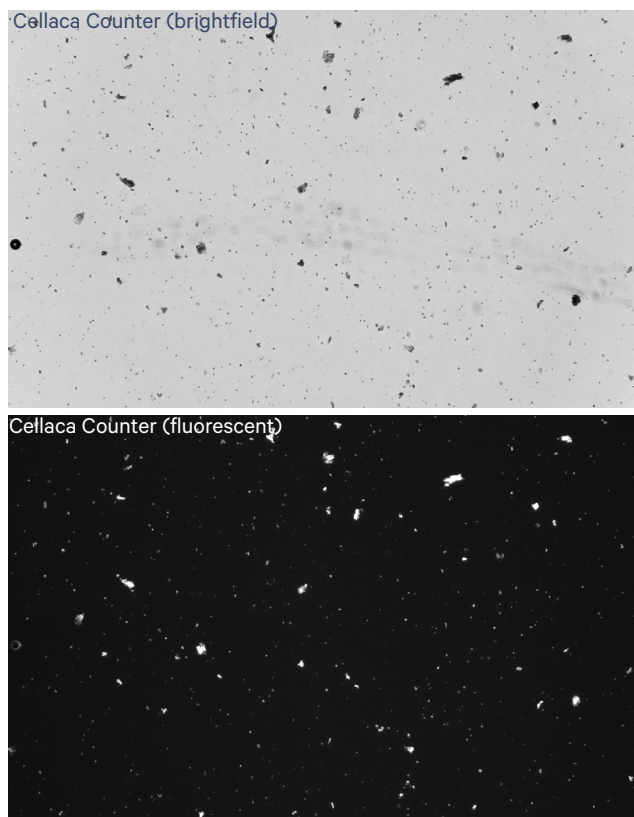

**Figure 1B: Representative Cells after Post-Hybridization Wash**

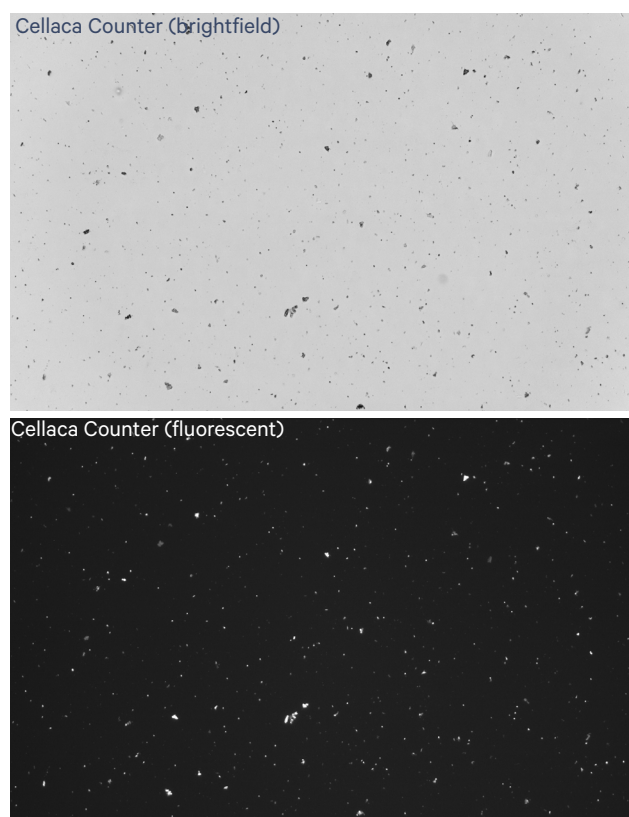

## Fixed Sample Storage Guidance

Fixed samples (dissociated cells resuspended in Quenching Buffer or in Tissue Resuspension Buffer) can be stored for short or long-term as described below.

### Short-term Storage at 4°C

- a. Thaw Enhancer (10x Genomics PN-2000482) for **10 min** at **65°C**. Vortex and centrifuge briefly. Keep warm and verify no precipitate before use.



*DO NOT keep the thawed reagent on ice, or the solution will precipitate. Once thawed, Enhancer can be kept at 42°C for up to 10 min.*

- b. Add **0.1 volume** pre-warmed Enhancer to fixed sample in Quenching Buffer or in Tissue Resuspension Buffer. For example, add **50 µl** Enhancer to **500 µl** fixed sample in Quenching Buffer or in Tissue Resuspension Buffer. Pipette mix.
- c. Store sample at **4°C** for up to **1 week**.

### Long-term Storage at –80°C

- a. Thaw Enhancer (10x Genomics PN-2000482) for **10 min** at **65°C**. Vortex and centrifuge briefly. Keep warm and verify no precipitate before use.



*DO NOT keep the thawed reagent on ice, or the solution will precipitate. Once thawed, Enhancer can be kept at 42°C for up to 10 min.*

- b. Add **0.1 volume** pre-warmed Enhancer to fixed sample in Quenching Buffer or in Tissue Resuspension Buffer. For example, add **50 µl** Enhancer to **500 µl** fixed sample in Quenching Buffer or in Tissue Resuspension Buffer. Pipette mix.
- c. Add 50% glycerol for a final concentration of 10%. For example: add **137.5 µl** 50% glycerol to **550 µl** fixed sample in Quenching Buffer or in Tissue Resuspension Buffer with Enhancer. Pipette mix.
- d. Store at **–80°C** for up to **3 months**.

## Post-Storage Processing

Samples may undergo a color change during storage (e.g. black, light gray, or green), however this will not impact assay performance.

If samples were stored at –80°C, thaw at room temperature until no ice is present.

- a. Centrifuge sample at **850 rcf** for **5 min** at **room temperature**.
- b. Remove the supernatant without disturbing the pellet.
- c. Resuspend cell pellet in **0.5 ml** 0.5X PBS + 0.02% BSA (optionally supplemented with 0.2 U/µl RNase Inhibitor) and keep on ice.

*Use RNase-free BSA at this step. See Specific Reagents & Consumables for details.*

- d. Determine cell concentration of the fixed sample using an Automated Cell Counter (Countess II/Cellaca MX) or hemocytometer. See Counting guidelines.
- e. Proceed **immediately** to appropriate Chromium Fixed RNA Profiling protocols (see References).

## Cell Yields from FFPE Tissue Sections

Listed below are the cell yields from FFPE tissue sections processed as described in this protocol. Based on this information, approximate cell yield from a specific tissue type/section may be estimated. Note that the tissue density in combination with tissue cross section area (tissue volume) will impact final cell yields.

|       | Tissue Type        | Tissue State<br>(Healthy or Diseased) | No. of Sections x<br>Thickness ( $\mu\text{m}$ ) | Tissue<br>Cross Section<br>( $x \times y \text{ mm}$ ) | Cell Yields         |                                |
|-------|--------------------|---------------------------------------|--------------------------------------------------|--------------------------------------------------------|---------------------|--------------------------------|
|       |                    |                                       |                                                  |                                                        | Manual Protocol     | gentleMACS<br>Octo Dissociator |
| MOUSE | Brain (Cerebellum) | Healthy                               | 6 x 50 $\mu\text{m}$                             | 2 x 2 mm                                               | $1.88 \times 10^6$  | $1.24 \times 10^6$             |
|       | Brain (Forebrain)  | Healthy                               | 6 x 50 $\mu\text{m}$                             | 4 x 5 mm                                               | $0.545 \times 10^6$ | $0.2 \times 10^6$              |
|       | Kidney             | Healthy                               | 3 x 50 $\mu\text{m}$                             | 5 x 9 mm                                               | $4.73 \times 10^6$  | $3.89 \times 10^6$             |
|       | Liver              | Healthy                               | 3 x 50 $\mu\text{m}$                             | 6 x 9 mm                                               | $1.7 \times 10^6$   | $1.6 \times 10^6$              |
|       | Spleen             | Healthy                               | 3 x 50 $\mu\text{m}$                             | 1 x 9 mm                                               | $5.15 \times 10^6$  | $4.865 \times 10^6$            |
| HUMAN | Brain              | Healthy                               | 2 x 25 $\mu\text{m}$                             | 5 x 10 mm                                              | $0.4 \times 10^6$   | $0.36 \times 10^6$             |
|       | Brain              | Glioblastoma                          | 2 x 25 $\mu\text{m}$                             | 10 x 13 mm                                             | $0.77 \times 10^6$  | $0.76 \times 10^6$             |
|       | Breast             | Healthy                               | 2 x 25 $\mu\text{m}$                             | 10 x 6 mm                                              | $0.04 \times 10^6$  | $0.02 \times 10^6$             |
|       | Breast             | Invasive Ductal<br>Carcinoma          | 2 x 25 $\mu\text{m}$                             | 13 x 7 mm                                              | $1.47 \times 10^6$  | $1.87 \times 10^6$             |
|       | Cervix             | Endocervical<br>Adenocarcinoma        | 2 x 25 $\mu\text{m}$                             | 10 x 7 mm                                              | $1.2 \times 10^6$   | $1.7 \times 10^6$              |
|       | Colon              | Cancer, Colorectal                    | 2 x 25 $\mu\text{m}$                             | 10 x 9 mm                                              | $0.6 \times 10^6$   | $0.76 \times 10^6$             |
|       | Kidney             | Healthy                               | 2 x 25 $\mu\text{m}$                             | 7 x 6 mm                                               | $0.36 \times 10^6$  | $0.53 \times 10^6$             |
|       | Liver              | Healthy                               | 2 x 25 $\mu\text{m}$                             | 5 x 6 mm                                               | $0.38 \times 10^6$  | $0.61 \times 10^6$             |
|       | Liver              | Hepatocellular<br>Carcinoma           | 2 x 25 $\mu\text{m}$                             | 5 x 9 mm                                               | $0.56 \times 10^6$  | $0.67 \times 10^6$             |
|       | Lung               | Healthy                               | 2 x 25 $\mu\text{m}$                             | 9 x 7 mm                                               | $0.44 \times 10^6$  | $0.79 \times 10^6$             |
|       | Lung               | Cancer                                | 2 x 25 $\mu\text{m}$                             | 9 x 8 mm                                               | $0.92 \times 10^6$  | $0.76 \times 10^6$             |
|       | Lymph Node         | Healthy                               | 2 x 25 $\mu\text{m}$                             | 7 x 5 mm                                               | $2.34 \times 10^6$  | $2.65 \times 10^6$             |
|       | Lymph Node         | Diseased, Reactive                    | 2 x 25 $\mu\text{m}$                             | 10 x 12 mm                                             | $6.75 \times 10^6$  | $7.45 \times 10^6$             |
|       | Ovary              | Cancer                                | 2 x 25 $\mu\text{m}$                             | 10 x 18 mm                                             | $2.2 \times 10^6$   | $1.92 \times 10^6$             |
|       | Pancreas           | Healthy                               | 2 x 25 $\mu\text{m}$                             | 11 x 9 mm                                              | $0.38 \times 10^6$  | $0.58 \times 10^6$             |
|       | Prostate           | Cancer                                | 2 x 25 $\mu\text{m}$                             | 9 x 9 mm                                               | $0.8 \times 10^6$   | $1.6 \times 10^6$              |
|       | Skin               | Malignant Melanoma                    | 2 x 25 $\mu\text{m}$                             | 4 x 4 mm                                               | $0.22 \times 10^6$  | $0.23 \times 10^6$             |

## Troubleshooting

| Problem                                             | Solution                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                              |
|-----------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Tissue block not adequately rehydrated              | Soak in water longer, until rehydrated                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                |
| Tissue section curl breaks during deparaffinization | <p>Ensure tissue blocks is adequately rehydrated<br/>(see adjacent representative tissue section images from well rehydrated and less rehydrated tissue blocks)</p> <p>DO NOT pipette directly on top of the curls; pipette against the tube wall</p>                                                                                                                                                                                                                                                                                 |
| Tissue section curl breaks during rehydration       | Centrifuge at 850rcf for 1 min prior to reagent exchange to pellet tissue pieces at the bottom of the tubes                                                                                                                                                                                                                                                                                                                                                                                                                           |
| Lower than expected cell yields                     | <p>Refer to Cell Yields from FFPE Tissue Sections table when planning the experiment</p> <p>If possible, repeat sample preparation protocol with more sections to increase yield</p> <p>Lower cell input recommendation for multiplex hybridization. Consider pooling cells derived from tissues of similar complexity</p> <p>If executing probe hybridization step with lower cell input, ensure that no cells are removed during post-hybridization washes by leaving behind slightly more buffer in the tube during the washes</p> |
| Debris/large chunks during cell counting            | <p>Pass sample through a 30 <math>\mu</math>m filter followed by rinsing the filter with PBS to minimize cell loss</p> <p>Remaining debris (that passes through the 30 <math>\mu</math>m filter) will be further reduced during the post-hybridization wash step in the downstream workflow (see Figure 1B)</p>                                                                                                                                                                                                                       |



### Using gentleMACS Octo Dissociator

|                                                     |                                                                                                                                     |
|-----------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------|
| Intact scroll at the end of run or run fails midway | Run a "spin only" program on the Octo Dissociator with steps 4-9 (from the protocol in step 3C)                                     |
| Octo Dissociator without heated lids                | Use a water bath and every 14 min use the "spin only" program on the Octo Dissociator with steps 4-9 (from the protocol in step 3C) |

### Representative Tissue Section Images

From well rehydrated blocks

From less rehydrated blocks

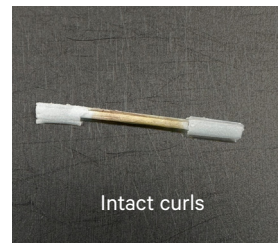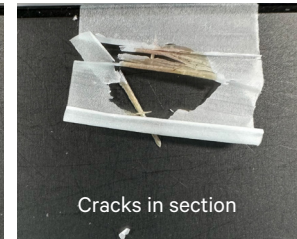

## References

1. Chromium Fixed RNA Profiling Reagent Kit for Singleplexed Samples with Feature Barcode technology for Cell Surface Protein User Guide\* (CG000477)
2. Chromium Fixed RNA Profiling Reagent Kit for Multiplexed Samples User Guide (CG000527)

\*Using cells derived from FFPE tissue sections with Feature Barcode technology for cell surface protein detection is not supported.

## Document Revision Summary

|                        |                                                                               |
|------------------------|-------------------------------------------------------------------------------|
| <b>Document Number</b> | CG000632                                                                      |
| <b>Title</b>           | Isolation of Cells from FFPE Tissue Sections for Chromium Fixed RNA Profiling |
| <b>Revision</b>        | Rev A                                                                         |
| <b>Revision Date</b>   | December 2022                                                                 |

© 2022 10x Genomics, Inc. (10x Genomics). All rights reserved. Duplication and/or reproduction of all or any portion of this document without the express written consent of 10x Genomics, is strictly forbidden. Nothing contained herein shall constitute any warranty, express or implied, as to the performance of any products described herein. Any and all warranties applicable to any products are set forth in the applicable terms and conditions of sale accompanying the purchase of such product. 10x Genomics provides no warranty and hereby disclaims any and all warranties as to the use of any third-party products or protocols described herein. The use of products described herein is subject to certain restrictions as set forth in the applicable terms and conditions of sale accompanying the purchase of such product. A non-exhaustive list of 10x Genomics' marks, many of which are registered in the United States and other countries can be viewed at: [www.10xgenomics.com/trademarks](http://www.10xgenomics.com/trademarks). 10x Genomics may refer to the products or services offered by other companies by their brand name or company name solely for clarity, and does not claim any rights in those third-party marks or names. 10x Genomics products may be covered by one or more of the patents as indicated at: [www.10xgenomics.com/patents](http://www.10xgenomics.com/patents). All products and services described herein are intended FOR RESEARCH USE ONLY and NOT FOR USE IN DIAGNOSTIC PROCEDURES.

The use of 10x Genomics products in practicing the methods set forth herein has not been validated by 10x Genomics, and such non-validated use is NOT COVERED BY 10X GENOMICS STANDARD WARRANTY, AND 10X GENOMICS HEREBY DISCLAIMS ANY AND ALL WARRANTIES FOR SUCH USE. Nothing in this document should be construed as altering, waiving or amending in any manner 10x Genomics terms and conditions of sale for the Chromium Controller or the Chromium Single Cell Controller, consumables or software, including without limitation such terms and conditions relating to certain use restrictions, limited license, warranty and limitation of liability, and nothing in this document shall be deemed to be Documentation, as that term is set forth in such terms and conditions of sale. Nothing in this document shall be construed as any representation by 10x Genomics that it currently or will at any time in the future offer or in any way support any application set forth herein.

**Contact:****[support@10xgenomics.com](mailto:support@10xgenomics.com)**

10x Genomics  
6230 Stoneridge Mall Road  
Pleasanton, CA 94588 USA

