

DEMONSTRATED PROTOCOL

Fresh Frozen Human Peripheral Blood Mononuclear Cells for Single Cell RNA Sequencing

Overview

This protocol outlines cryopreservation and thawing of human peripheral blood mononuclear cells (PBMCs) for use with 10x Genomics Single Cell protocols. While this Demonstrated Protocol is specific to PBMCs, the protocol may be used as a basis for handling other primary cells in preparation for use in the 10x Genomics Single Cell protocols.

Additional Guidance

Consult Demonstrated Protocol Cell Preparation Guide (Document CG00053) for Tips & Best Practices on handling cells and Technical Note Guidelines for Accurate Target Cell Counts using 10x Genomics Single Cell Solutions (Document CG000091) for determining accurate cell counts.

Cells carry potentially hazardous pathogens. Follow material supplier recommendations and local laboratory procedures and regulations for the safe handling, storage and disposal of biological materials.

Cell Sourcing

This protocol was demonstrated using fresh human PBMCs from AllCells.



Preparation-Buffers

Cryopreservation	
Media/Buffers	Composition
Resuspension Medium (maintain at 4°C)	40% FBS in cell culture media (e.g., IMDM/RPMI)
2X Freezing Medium (maintain at 4°C)	30% DMSO in cell culture media (e.g., IMDM/RPMI) containing 40% FBS
Thawing & Resuspension	
Media/Buffers	Composition
Complete Growth Medium (maintain at 37°C)	10% FBS in cell culture media (e.g., IMDM/RPMI)
PBS + 0.04% BSA (maintain at room temperature)	

Specific Reagents & Consumables

Vendor	Item	Part Number
Thermo Fisher Scientific	Gibco IMDM	12440-053
	UltraPure Bovine Serum Albumin (BSA, 50 mg/ml)	AM2616
	Trypan Blue Stain (0.4%)	T10282
	Nunc Biobanking & Cell Culture Cryogenic Tubes, 1.8 ml	368632
	Countess II FL Automated Cell Counter	AMAQAF1000
	Countess II FL Automated Cell Counting Chamber Slides	C10228
Fisher Scientific	Dimethyl Sulfoxide (DMSO), for molecular biology	ICN19141880
Millipore-Sigma	Bovine Serum Albumin in DPBS (10%) (alternative to Thermo Fisher product)	A1595
	Phosphate-Buffered Saline (PBS) with 10% Bovine Albumin (alternative to Thermo Fisher product)	SRE0036
Miltenyi	MACS SmartStrainers, 30 µm	130-098-458
Bel-Art	Flowmi Cell Strainer, 40 µm (alternative to Miltenyi product)	H13680-0040
Corning	Corning RPMI 1640	10-040-CM
	Phosphate-Buffered Saline without Calcium & Magnesium	21-040-CV
VWR	Seradigm Premium Grade Fetal Bovine Serum (FBS)	97068-085
	Sterile Polypropylene Centrifuge Tubes with Flat Caps, 50 ml	82018-050
Biocision	CoolCell FTS30 Cell Freezing Container	BSC-170
Eppendorf	DNA LoBind Tubes, 2.0 ml	022431048

Protocol Overview

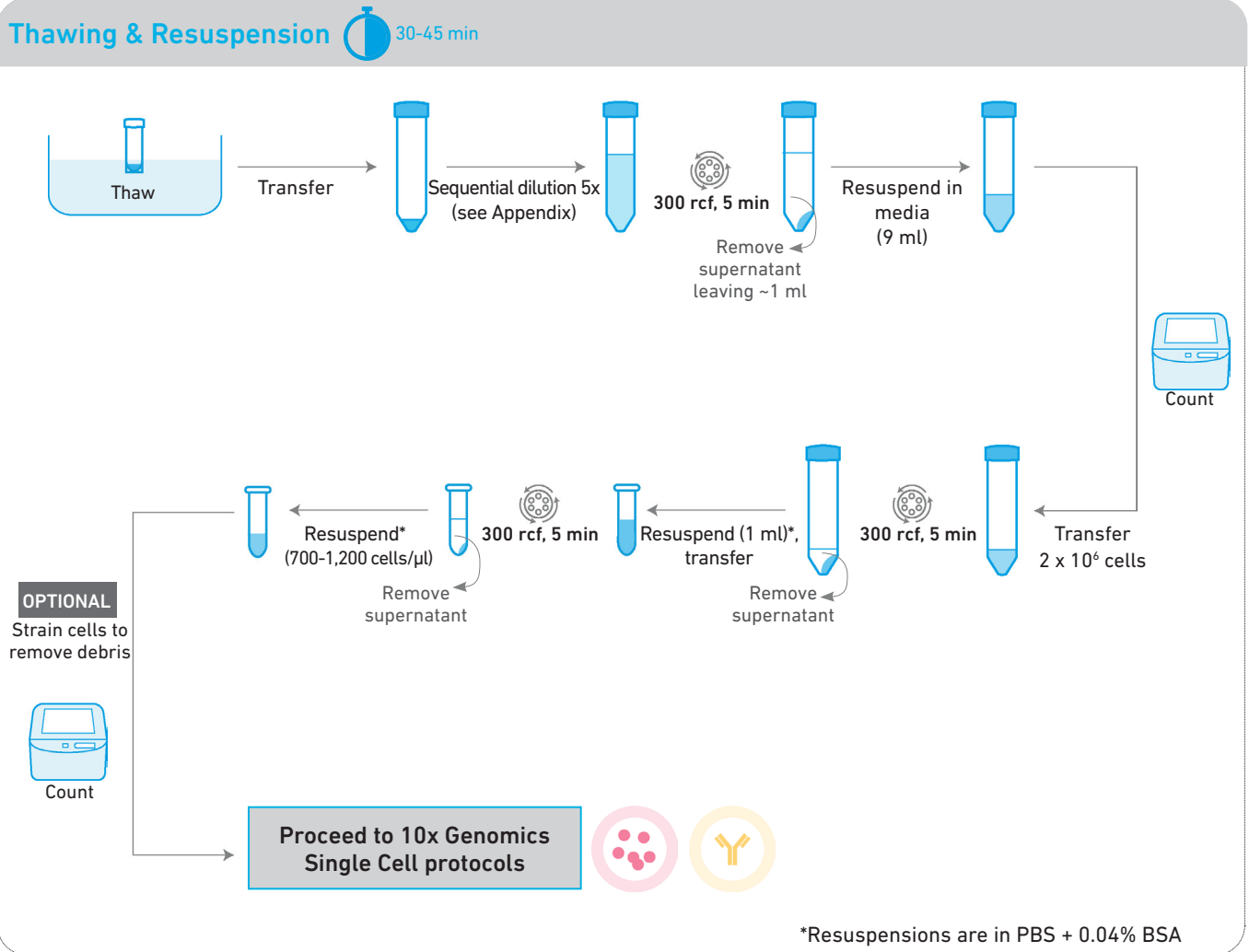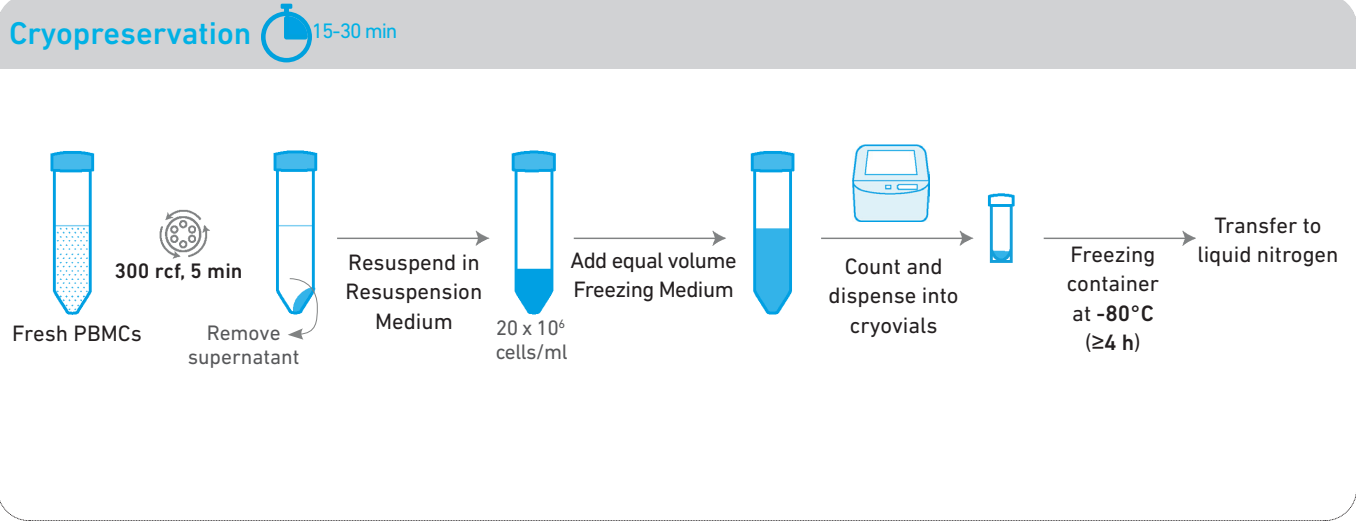

Protocol

Cryopreservation


Pre-cool a cell freezing container by placing it on ice before starting cryopreservation.

- Place PBMCs on ice.
- Gently mix the cells.
- Determine cell viability and total cell number using a Countess II Automated Cell Counter.
- Centrifuge at **300 rcf** for **5 min** at **4°C**.
- Remove the supernatant.
- Resuspend the cell pellet in an appropriate volume of chilled Resuspension Medium to achieve a concentration of 20×10^6 cells/ml. Maintain the cells on ice.
- Add an equivalent volume of chilled 2X Freezing Medium to achieve a concentration of 10×10^6 cells/ml. Gently mix the cells.
- Dispense cell suspension aliquots into the pre-cooled cryovials and place the cryovials inside a pre-cooled cell freezing container e.g., CoolCell FTS30.
- Place the cell freezing container in a -80°C freezer for ≥ 4 h. After 4 h, transfer the cryovials to liquid nitrogen for long-term storage.

Thawing & Resuspension


Set up a water bath to **37°C** before starting cell thawing. All cell washes are performed at **room temperature**.

- Remove cryovials from storage and **immediately** thaw in the water bath at **37°C** for **2-3 min**.

 **DO NOT** submerge the entire vial in the water bath. Remove from the water bath when a tiny ice crystal remains.

- In a biosafety hood, slowly transfer thawed cells to a 50-ml conical tube using a **wide-bore** pipette tip. Rinse the cryovial with **1 ml** warm complete growth medium and add the rinse **dropwise** (1 drop per 5 sec) to the 50-ml conical tube while gently shaking the tube.
- Sequentially dilute cells in the 50-ml conical tube by incremental 1:1 volume additions of media for a total of 5 times (including dilution at step b) with **~1 min** wait between additions (see Appendix). Add media at a speed of 1 ml/3-5 sec to the tube and swirl.
- Centrifuge at **300 rcf** for **5 min**.

- Remove most of the supernatant, leaving **~1 ml** and resuspend cell pellet in this volume using a **regular-bore** pipette tip.

 Cell pellet may be present on the side or on the bottom of the tube.




- Add an additional **9 ml** complete growth medium (at a speed of 1 ml/ 3-5 sec) to achieve a total volume of **~10 ml**.
- Determine the cell concentration using a Countess II Automated Cell Counter.
- Transfer $\sim 2 \times 10^6$ cells into a new 50-ml tube.
- Centrifuge at **300 rcf** for **5 min**.
- Remove the supernatant without disrupting the cell pellet.
- Using a **wide-bore** pipette tip, add **1 ml** 1X PBS + 0.04% BSA and gently pipette mix 5x.
- Transfer the cells into a 2-ml microcentrifuge tube. Rinse the 50-ml conical tube with **0.5 ml** 1X PBS + 0.04% BSA and transfer the rinse into the 2-ml tube containing cells.
- Centrifuge at **300 rcf** for **5 min**.
- Remove the supernatant without disrupting the cell pellet.
- Based on starting cell concentration at step a (Thawing & Resuspension) and assuming $\sim 50\%$ cell loss, add an appropriate volume 1X PBS + 0.04% BSA to obtain a concentration of 700-1,200 cells/ μl . Gently pipette mix using a **regular-bore** pipette tip until a single cell suspension is achieved.

OPTIONAL If cell debris and large clumps are present, pass the sample through a 40 μm Flowmi Cell Strainer.

- Determine cell concentration and viability using a Countess II Automated Cell Counter or a hemocytometer. The targeted final concentration is 700-1,200 cells/ μl .
- Once the final cell concentration is achieved, place cells on ice.
- Proceed **immediately** to the 10x Genomics Single Cell protocols.

Appendix


Sequential Dilution with Media

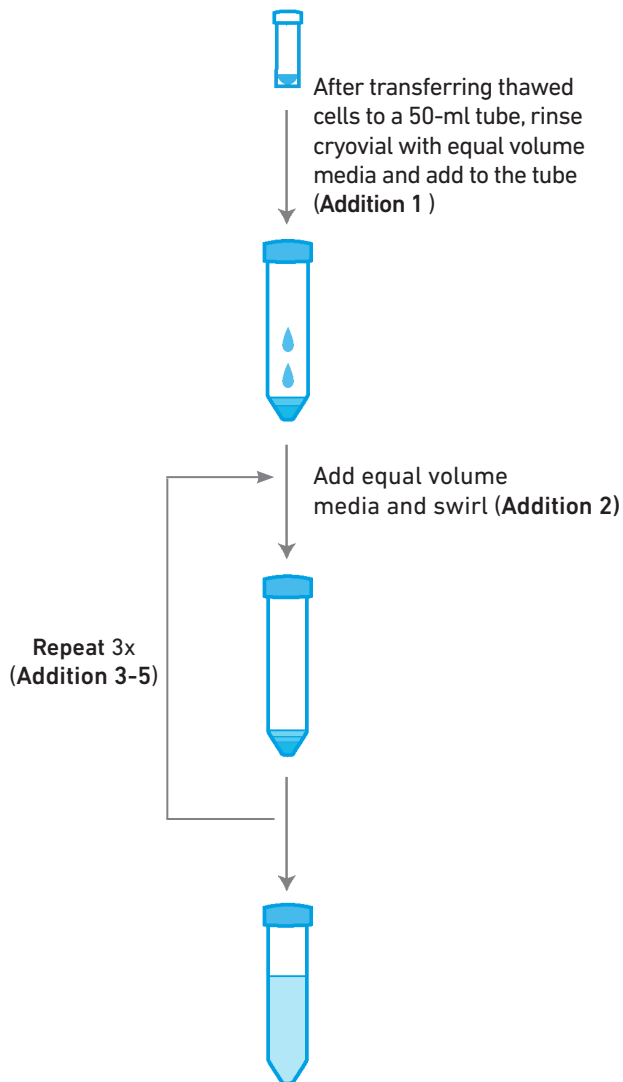

Results

The viability of PBMCs depends on both the freezing and thawing protocols, in addition to the total cell number per vial used during cryopreservation. The typical percent PBMCs viability obtained by following this protocol ranges from 88-93% based on both trypan blue and live/dead staining.

Example of Trypan Blue Staining

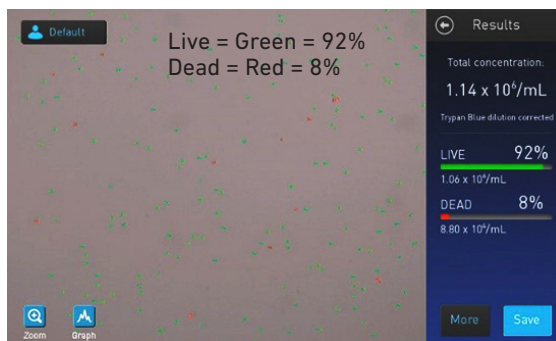

Revision Summary

Document Number	CG00039
Title	Fresh Frozen Human PBMCs for Single Cell RNA Sequencing
Revision	Rev D to Rev E
Revision Date	November 2021
Specific Changes	Added Bovine Serum Albumin in DPBS (10%) part name and number to the Specific Reagents & Consumables List
General Changes	Updated for general minor consistency of language and terms throughout.

© 2021 10x Genomics, Inc. (10x Genomics). All rights reserved. Duplication and/or reproduction of all or any portion of this document without the express written consent of 10x Genomics, is strictly forbidden. Nothing contained herein shall constitute any warranty, express or implied, as to the performance of any products described herein. Any and all warranties applicable to any products are set forth in the applicable terms and conditions of sale accompanying the purchase of such product. 10x Genomics provides no warranty and hereby disclaims any and all warranties as to the use of any third-party products or protocols described herein. The use of products described herein is subject to certain restrictions as set forth in the applicable terms and conditions of sale accompanying the purchase of such product. A non-exhaustive list of 10x Genomics' marks, many of which are registered in the United States and other countries can be viewed at: www.10xgenomics.com/ trademarks. 10x Genomics may refer to the products or services offered by other companies by their brand name or company name solely for clarity, and does not claim any rights in those third-party marks or names. 10x Genomics products may be covered by one or more of the patents as indicated at: www.10xgenomics.com/ patents. The use of products described herein is subject to 10x Genomics Terms and Conditions of Sale, available at <http://www.10xgenomics.com/legal-notices>, or such other terms that have been agreed to in writing between 10x Genomics and user. All products and services described herein are intended FOR RESEARCH USE ONLY and NOT FOR USE IN DIAGNOSTIC PROCEDURES.

The use of 10x Genomics products in practicing the methods set forth herein has not been validated by 10x Genomics, and such non-validated use is NOT COVERED BY 10x GENOMICS STANDARD WARRANTY, AND 10x GENOMICS HEREBY DISCLAIMS ANY AND ALL WARRANTIES FOR SUCH USE. Nothing in this document should be construed as altering, waiving or amending in any manner 10x Genomics terms and conditions of sale for the Chromium Controller or the Chromium Single Cell Controller, consumables or software, including without limitation such terms and conditions relating to certain use restrictions, limited license, warranty and limitation of liability, and nothing in this document shall be deemed to be Documentation, as that term is set forth in such terms and conditions of sale. Nothing in this document shall be construed as any representation by 10x Genomics that it currently or will at any time in the future offer or in any way support any application set forth herein.

Contact:
support@10xgenomics.com
 10x Genomics
 6230 Stoneridge Mall Road
 Pleasanton, CA 94588 USA

