



[User Guide](#) | [CG000477](#) | [Rev C](#)

# Chromium

## Fixed RNA Profiling Reagent Kits

**for Singleplexed Samples with Feature Barcode technology for Cell Surface Protein**

For use with:

Chromium Next GEM Single Cell Fixed RNA Sample Preparation Kit

*16 rxns PN-1000414*

Chromium Next GEM Chip Q Single Cell Kit

*48 rxns PN-1000418 | 16 rxns PN-1000422*

Chromium Fixed RNA Kit, Human Transcriptome,

*4 rxns x 1 BC PN-1000474*

Chromium Fixed RNA Kit, Mouse Transcriptome,

*4 rxns x 1 BC PN-1000495*

Fixed RNA Feature Barcode Kit *16 rxns PN-1000419*

Dual Index Kit TT Set A *96 rxns PN-1000215*

Dual Index Kit TS Set A *96 rxns PN-1000251*

# Notices

## Document Number

CG000477 | Rev C

## Legal Notices

© 2022 10x Genomics, Inc. (10x Genomics). All rights reserved. Duplication and/or reproduction of all or any portion of this document without the express written consent of 10x Genomics, is strictly forbidden. Nothing contained herein shall constitute any warranty, express or implied, as to the performance of any products described herein. Any and all warranties applicable to any products are set forth in the applicable terms and conditions of sale accompanying the purchase of such product. 10x Genomics provides no warranty and hereby disclaims any and all warranties as to the use of any third-party products or protocols described herein. The use of products described herein is subject to certain restrictions as set forth in the applicable terms and conditions of sale accompanying the purchase of such product. A non-exhaustive list of 10x Genomics' marks, many of which are registered in the United States and other countries can be viewed at: [www.10xgenomics.com/trademarks](http://www.10xgenomics.com/trademarks). 10x Genomics may refer to the products or services offered by other companies by their brand name or company name solely for clarity, and does not claim any rights in those third-party marks or names. 10x Genomics products may be covered by one or more of the patents as indicated at: [www.10xgenomics.com/patents](http://www.10xgenomics.com/patents). The use of products described herein is subject to 10x Genomics Terms and Conditions of Sale, available at [www.10xgenomics.com/legal-notices](http://www.10xgenomics.com/legal-notices), or such other terms that have been agreed to in writing between 10x Genomics and user. All products and services described herein are intended FOR RESEARCH USE ONLY and NOT FOR USE IN DIAGNOSTIC PROCEDURES.

## Instrument & Licensed Software Updates Warranties

Updates to existing Instruments and Licensed Software may be required to enable customers to use new or existing products. In the event of an Instrument failure resulting from an update, such failed Instrument will be replaced or repaired in accordance with the 10x Limited Warranty, Assurance Plan or service agreement, only if such Instrument is covered by any of the foregoing at the time of such failure. Instruments not covered under a current 10x Limited Warranty, Assurance Plan or service agreement will not be replaced or repaired.

## Support

Email: [support@10xgenomics.com](mailto:support@10xgenomics.com)  
10x Genomics  
6230 Stoneridge Mall Road  
Pleasanton, CA

# Document Revision Summary

## Document Number

CG000477 | Rev C

## Title

Chromium Fixed RNA Profiling Reagent Kits for Singleplexed Samples with Feature Barcode technology for Cell Surface Protein

## Revision

Rev B to Rev C

## Revision Date

October 25, 2022

## Specific Changes

- Updated sample storage guidelines after post-hybridization washing in step 2.1.
- Updated to include mouse transcriptome kit (Chromium Fixed RNA Kit, Mouse Transcriptome, 4 rxns x 1 BC).
- Updated Sample Index PCR cycle number for Fixed RNA Gene Expression Library for cells from fixed & dissociated tissues and cells from FFPE tissue sections in step 4.1.

## General Changes

Updated for general minor consistency of language and terms throughout.

# Table of Contents

## Introduction

Chromium Fixed RNA Profiling Reagent Kits	7
10x Genomics Accessories	11
Recommended Thermal Cyclers	11
Additional Kits, Reagents & Equipment	13
Protocol Steps & Timing	17
Stepwise Objectives	18
Sample Fixation Guidelines	25

## Tips & Best Practices

### Step 1: Probe Hybridization

1.0 Get Started	41
1.1 Probe Hybridization	42

### Step 2: GEM Generation and Barcoding

2.0 Get Started	45
2.1 Post-Hybridization Wash	47
2.2 Prepare GEM Master Mix + Sample Dilution	50
2.3 Load Chromium Next GEM Chip Q	54
2.4 Run the Chromium X/iX	56
2.5 Transfer GEMs	57
2.6 GEM Incubation	58

### Step 3: GEM Recovery and Pre-Amplification

3.0 Get Started	60
3.1 Post-GEM Incubation – Recovery	61
3.2 Pre-Amplification PCR	62
3.3 DNA Cleanup – SPRIselect	64

### Step 4: Fixed RNA – Gene Expression Library Construction

4.0 Get Started	66
4.1 Sample Index PCR	67
4.2 Post Sample Index PCR Size Selection – SPRIselect	69
4.3 Post Library Construction QC	70



# Table of Contents

<b>Step 5: Fixed RNA – Cell Surface Protein Library Construction</b>	
5.0 Get Started	72
5.1 Sample Index PCR	73
5.2 Post Sample Index PCR Size Selection – SPRIselect	75
5.3 Post Library Construction QC	76
<b>Step 6: Sequencing</b>	
Sequencing Libraries	78
Illumina Sequencer Compatibility	78
Sample Indices	79
Fixed RNA – Gene Expression Library Sequencing Parameters	79
Fixed RNA – Cell Surface Protein Library Sequencing Parameters	79
Library Loading	80
Library Pooling	80
Data Analysis and Visualization	80
<b>Troubleshooting</b>	
GEM Generation & Barcoding	83
Post Library Construction QC	84
Chromium X Series Errors	85
<b>Appendix</b>	
Chromium Fixed RNA Profiling: Chip Loading Overview	87
Post Library Construction Quantification	88
Agilent TapeStation Traces	89
LabChip Traces	89
Oligonucleotide Sequences	91


# Introduction

Chromium Fixed RNA Profiling Reagent Kits	7
10x Genomics Accessories	11
Recommended Thermal Cyclers	11
Additional Kits, Reagents & Equipment	13
Protocol Steps & Timing	17
Stepwise Objectives	18
Sample Fixation Guidelines	25

## Chromium Fixed RNA Profiling Reagent Kits

*Refer to SDS for handling and disposal information*

### Chromium Next GEM Single Cell Fixed RNA Sample Preparation Kit, 16 rxns PN-1000414

Chromium Next GEM Single Cell Fixed RNA Sample Preparation Kit 16 rxns, PN-1000414 Store at -20°C			
		#	PN
●	Conc. Fix & Perm Buffer	2	2000517
●	Conc. Quench Buffer	2	2000516
●	Enhancer	2	2000482
			

## Chromium Fixed RNA Kit, Human Transcriptome, 4 rxns x 1 BC PN-1000474

Chromium Next GEM  
**Single Cell Fixed RNA Hybridization & Library Kit**  
 4 rxns, PN-1000415  
*Store at -20°C*

		#	PN
●	Hyb Buffer B	1	2000483
●	Enhancer	1	2000482
●	Conc. Post-Hyb Buffer	1	2000533
○	Reducing Agent B	1	2000087
●	GEM Enzyme Mix	1	2000490
●	GEM Reagent Mix	1	2000491
●	Pre-Amp Primers B	1	2000529
○	Amp Mix	1	2000103

**10x**  
GENOMICS™

Chromium Next GEM  
**Single Cell Fixed RNA Human Transcriptome Probe Kit**  
 4 rxns, PN-1000423  
*Store at -20°C*

		#	PN
●	Human WTA Probes BC001	1	2000495

**10x**  
GENOMICS™









Chromium Next GEM  
**Single Cell Fixed RNA Gel Bead Kit**  
 4 rxns, PN-1000421  
*Store at -80°C*

		#	PN
	Single Cell TL v1 Gel Beads (4 rxns)	1	2000538

**10x**  
GENOMICS™


## Chromium Fixed RNA Kit, Mouse Transcriptome, 4 rxns x 1 BC PN-1000495

Chromium Next GEM  
**Single Cell Fixed RNA Hybridization & Library Kit**  
 4 rxns, PN-1000415  
*Store at -20°C*

		#	PN
	Hyb Buffer B	1	2000483
	Enhancer	1	2000482
	Conc. Post-Hyb Buffer	1	2000533
	Reducing Agent B	1	2000087
	GEM Enzyme Mix	1	2000490
	GEM Reagent Mix	1	2000491
	Pre-Amp Primers B	1	2000529
	Amp Mix	1	2000103

**10x**  
GENOMICS

Chromium Next GEM  
**Single Cell Fixed RNA Mouse Transcriptome Probe Kit**  
 4 rxns, PN-1000490  
*Store at -20°C*

		#	PN
	Mouse WTA Probes BC001	1	2000703

**10x**  
GENOMICS

Chromium Next GEM  
**Single Cell Fixed RNA Gel Bead Kit**  
 4 rxns, PN-1000421  
*Store at -80°C*

	#	PN
Single Cell TL v1 Gel Beads (4 rxns)	1	2000538


**10x**  
GENOMICS

### Chromium Next GEM Chip Q Single Cell Kit, 48 rxns PN-1000418

Chromium Partitioning Oil <i>Store at ambient temperature</i>			
	#	PN	
<input checked="" type="radio"/> Partitioning Oil	6	2000190	

Chromium Recovery Agent <i>Store at ambient temperature</i>			
	#	PN	
<input type="radio"/> Recovery Agent	6	220016	

Chromium Next GEM Chip Q & Gaskets <i>Store at ambient temperature</i>			
	#	PN	
Next GEM Chip Q	6	2000518	
Chip Gasket, 6-pack	1	370017	




### Chromium Next GEM Chip Q Single Cell Kit, 16 rxns PN-1000422



Chromium Partitioning Oil <i>Store at ambient temperature</i>			
	#	PN	
<input checked="" type="radio"/> Partitioning Oil	2	2000190	

Chromium Recovery Agent <i>Store at ambient temperature</i>			
	#	PN	
<input type="radio"/> Recovery Agent	2	220016	

Chromium Next GEM Chip Q & Gaskets <i>Store at ambient temperature</i>			
	#	PN	
Next GEM Chip Q	2	2000518	
Chip Gasket, 2-pack	1	3000072	



## Fixed RNA Feature Barcode Kit, 16 rxns PN-1000419

Fixed RNA Feature Barcode Kit 16 rxns, PN-1000419 Store at -20°C			
		#	PN
	Feature Pre-Amp Primers	1	2000515
	Amp Mix	1	2000047
<b>10x</b> GENOMICS™			

## Dual Index Kit TT Set A, 96 rxns PN-1000215

Dual Index Kit TT Set A Store at -20°C		
	#	PN
Dual Index Plate TT Set A	1	3000431

## Dual Index Kit TS Set A, 96 rxns PN-1000251

Dual Index Kit TS Set A Store at -20°C		
	#	PN
Dual Index Plate TS Set A	1	3000511

## 10x Genomics Accessories

Product	Part Number (Kit)	Part Number (Item)
10x Vortex Adapter	120251	330002
10x Magnetic Separator	120250	230003
Chromium Next GEM Secondary Holder	1000142	3000332

## Recommended Thermal Cyclers

The table below lists the thermal cyclers that have been validated by 10x Genomics. Thermal cyclers used must support uniform heating of 100 µl emulsion volumes

Supplier	Description	Part Number
BioRad	C1000 Touch Thermal Cycler with 96-Deep Well Reaction Module	1851197
Eppendorf	MasterCycler X50s	North America 6311000010
Analytik Jena	Biometra TAdvanced 96 SG with 96-well block (silver, 0.2 mL) and gradient function	846-x-070-241

If using thermal cyclers other than BioRad C1000, the ramp rates should be adjusted for all the steps as described below

- Eppendorf MasterCycler X50s: 3°C/sec heating and 2°C/sec cooling
- Analytik Jena Biometra TAdvanced: 2°C/sec heating and cooling



## Additional Kits, Reagents & Equipment

The items in the table below have been validated by 10x Genomics and are highly recommended for 10x Genomics workflows, training, and system operations. Substituting materials may adversely affect system performance. This list does not include standard laboratory equipment, such as water baths, centrifuges, vortex mixers, pH meters, freezers, etc.

Supplier	Description	Part Number	
<b>Plastics</b>			
Eppendorf	DNA LoBind Tubes, 1.5 ml	0224-31021	
	DNA LoBind Tubes, 5.0 ml	0030108310	
	PCR Tubes 0.2 ml 8-tube strips	Choose either Eppendorf, USA Scientific or Thermo Fisher Scientific PCR 8-tube strips.	
USA Scientific	TempAssure PCR 8-tube strip		1402-4700
Thermo Fisher Scientific	MicroAmp 8-Tube Strip, 0.2 ml		N8010580
	MicroAmp 8-Cap Strip, clear	N8010535	
Corning	Corning Centrifuge Tubes with CentriStar Cap (15 ml), sterile	430790	
	Self-Standing Polypropylene Centrifuge Tubes (50 ml), sterile	430921	
Sysmex	Sterile single-pack CellTrics filters	Choose either Sysmex or Miltenyi Biotec filter.	
Miltenyi Biotec	Pre-Separation Filters (30 µm) <i>Alternative to Sysmex product</i>		130-041-407
<b>Kits &amp; Reagents</b>			
Thermo Fisher Scientific	Nuclease-free Water	AM9937	
	Low TE Buffer (10 mM Tris-HCl pH 8.0, 0.1 mM EDTA)	12090-015	
	Tween 20 Surfact-Amps Detergent Solution (10% solution)	28320	
Beckman Coulter	SPRIselect Reagent Kit	B23318	
Millipore Sigma	Ethyl Alcohol, Pure (200 Proof, anhydrous)	E7023-500ML	
	Glycerol for Molecular biology, ≥99.0%	G5516	
Acros Organics	Glycerol, 99.5%, for molecular biology, DNase, RNase and Protease free <i>Alternative to Millipore Sigma product</i>	327255000	
Qiagen	Qiagen Buffer EB	19086	
Ricca Chemical Company	Glycerin (glycerol), 50% (v/v) Aqueous Solution	3290-32	
<b>Cell Counting</b>			
Thermo Fisher Scientific	Ethidium Homodimer-1	Choose Ethidium Homodimer-1, Trypan Blue, or PI Staining Solution based on the presence of debris in the sample.	E1169
	Trypan Blue Stain (0.4%)		T10282
	Countess II FL Automated Cell Counter <i>Discontinued</i>		AMAQAF1000
	Countess II FL Automated Cell Counting Chamber Slides		C10228
Nexcelom Biosciences	Cellaca MX High-throughput Automated Cell Counter	MX-112-0127	
	ViaStain PI Staining Solution	CS1-0109-5mL	
	ViaStain AO/PI Staining Solution <i>Alternative to PI Staining Solution</i>	CS2-0106-5mL	

Equipment			
VWR	VWR Mini Centrifuge <i>Or any equivalent mini centrifuge</i>	76269-064	
	Divided Polystyrene Reservoirs	41428-958	
	Vortex Mixer	10153-838	
Eppendorf	Eppendorf ThermoMixer C	5382000023	
	Eppendorf ThermoTop with condens.protect technology	5308000003	
	Eppendorf SmartBlock 1.5 ml, Thermoblock for 24 reaction vessels	5360000038	
	Eppendorf SmartBlock 2.0 ml, thermoblock for 24 reaction vessels	5362000035	
Quantification & Quality Control			
Agilent	2100 Bioanalyzer Laptop Bundle	G2953CA	
	High Sensitivity DNA Kit	Choose Bioanalyzer, TapeStation or LabChip based on availability & preference.	5067-4626
	4200 TapeStation		G2991AA
	High Sensitivity D5000: ScreenTape/ Reagents		5067-5592/ 5067-5593
	High Sensitivity D1000: ScreenTape/ Reagents		5067-5584/ 5067-5585
PerkinElmer	LabChip GX Touch HT Nucleic Acid Analyzer		CLS137031
	DNA High Sensitivity Reagent Kit	CLS760672	
KAPA Biosystems	KAPA Library Quantification Kit for Illumina Platforms	KK4824	

This list does not include reagents required for sample fixation and cell surface protein labeling. See [Sample Fixation Guidelines on page 25](#) for the key reagents needed for fixation and consult Fixation of Cells & Nuclei for Chromium Fixed RNA Profiling (Demonstrated Protocol CG000478) for a complete list of reagents for sample fixation and storage.










## Recommended Pipette Tips

10x Genomics recommends using only validated emulsion-safe pipette tips for all Single Cell protocols. Rainin pipette tips have been extensively validated by 10x Genomics and are highly recommended for all single cell assays. If Rainin tips are unavailable, any of the listed alternate pipette tips validated by 10x Genomics may be used.

Supplier	Description	Part Number (US)
<b>Recommended Pipettes &amp; Pipette tips</b>		
Rainin	<b>Pipettes</b>	
	Pipet-Lite Multi Pipette L8-50XLS+	17013804
	Pipet-Lite Multi Pipette L8-200XLS+	17013805
	Pipet-Lite Multi Pipette L8-10XLS+	17013802
	Pipet-Lite Multi Pipette L8-20XLS+	17013803
	Pipet-Lite LTS Pipette L-2XLS+	17014393
	Pipet-Lite LTS Pipette L-10XLS+	17014388
	Pipet-Lite LTS Pipette L-20XLS+	17014392
	Pipet-Lite LTS Pipette L-100XLS+	17014384
	Pipet-Lite LTS Pipette L-200XLS+	17014391
	Pipet-Lite LTS Pipette L-1000XLS+	17014382
	<b>Pipette Tips</b>	
	Tips LTS 200UL Filter RT-L200FLR	30389240
	Tips LTS 1ML Filter RT-L1000FLR	30389213
Tips LTS 20UL Filter RT-L10FLR	30389226	
<b>Alternate Recommendations</b>		
<i>(If Rainin pipette tips are unavailable, any of the listed pipette tips may be used)</i>		
Eppendorf	<b>Pipettes</b>	
	Eppendorf Research Plus, 8-channel, epT.I.P.S. Box, 0.5 – 10 µl	3125000010
	Eppendorf Research Plus, 8-channel, epT.I.P.S. Box, 10 – 100 µl	3125000036
	Eppendorf Research Plus, 8-channel, epT.I.P.S. Box, 30 – 300 µl	3125000052
	Eppendorf Research Plus, 1-channel, epT.I.P.S.® Box, 0.1 – 2.5 µl	3123000012
	Eppendorf Research Plus, 1-channel, epT.I.P.S.® Box, 0.5 – 10 µl	3123000020
	Eppendorf Research Plus, 1-channel, epT.I.P.S.® Box, 2 – 20 µl	3123000039
	Eppendorf Research Plus, 1-channel, epT.I.P.S.® Box, 2 – 200 µl	3123000055
Eppendorf Research Plus, 1-channel, epT.I.P.S.® Box, 100 – 1000 µl	3123000063	

Supplier	Description	Part Number (US)
	<b>Pipette Tips</b> (compatible with Eppendorf pipettes only)	
	ep Dualfilter T.I.P.S., 2-20 µl	0030078535
	ep Dualfilter T.I.P.S., 2-200 µl	0030078551
	ep Dualfilter T.I.P.S., 2-1,000 µl	0030078578
Labcon*	ZAP SLIK 20 µL Low Retention Aerosol Filter Pipet Tips for Rainin LTS	4-1143-965-008
	ZAP SLIK 200 µL Low Retention Aerosol Filter Pipet Tips for Rainin LTS	4-1144-965-008
	ZAP SLIK 1200 µL Low Retention Aerosol Filter Pipet Tips for Rainin LTS	4-1145-965-008
Biotix*	xTIP4 Racked Pipette Tips, Rainin LTS Pipette Compatible, 0.1-20 µl	63300931
	xTIP4 Racked Pipette Tips, Rainin LTS Pipette Compatible, 200 µl	63300001
	xTIP4 Racked Pipette Tips, Rainin LTS Pipette Compatible, 1200 µl	63300004
*Compatible with Rainin pipettes		

## Protocol Steps & Timing

Steps	Timing	Stop & Store
<b>Sample Fixation</b>	2-3 h	 4°C ≤1 week/-20°C or -80°C ≤6 months (fixed single cells/nuclei) 4°C ≤1 week (fixed tissue pieces)
<b>Step 1: Probe Hybridization (page 40)</b>		
1.1 Probe Hybridization (page 42)	16-24 h	
<b>Step 2: GEM Generation and Barcoding (page 44)</b>		
2.1 Post-Hybridization Wash (page 47)	60 min	 -20°C/-80°C ≤1 month
2.2 Prepare GEM Master Mix + Sample Dilution (page 50)	30 min	
2.3 Load Chromium Next GEM Chip Q (page 54)	10 min	
2.4 Run the Chromium X/iX (page 56)	5.5 min	
2.5 Transfer GEMs (page 57)	5 min	
2.6 GEM Incubation (page 58)	125 min	 4°C ≤1 week (GEMs) -20°C/-80°C ≤1 month (washed undiluted sample)
<b>Step 3: GEM Recovery and Pre-Amplification (page 59)</b>		
3.1 Post-GEM Incubation – Recovery (page 61)	10 min	
3.2 Pre-Amplification PCR (page 62)	55 min	 4°C ≤72 h/-20°C ≤1 week
3.3 DNA Cleanup – SPRIselect (page 64)	30 min	 4°C ≤72 h/-20°C ≤4 weeks
<b>Step 4: Fixed RNA – Gene Expression Library Construction (page 65)</b>		
4.1 Sample Index PCR (page 67)	40 min	 4°C ≤72 h
4.2 Post Sample Index PCR Size Selection – SPRIselect (page 69)	30 min	 4°C ≤72 h or -20°C long-term
4.3 Post Library Construction QC (page 70)	60 min	
<b>Step 5: Fixed RNA – Cell Surface Protein Library Construction on page 71</b>		
5.1 Sample Index PCR (page 73)	40 min	 4°C ≤72 h
5.2 Post Sample Index PCR Size Selection – SPRIselect (page 75)	30 min	 4°C ≤72 h or -20°C long-term
4.3 Post Library Construction QC (page 70)	60 min	

## Stepwise Objectives


Chromium Fixed RNA Profiling offers comprehensive, scalable solutions to measure gene expression in formaldehyde fixed samples. Gene expression is measured using probe pairs designed to hybridize to mRNA specifically. Using a microfluidic chip, the fixed and probe-hybridized single cell and nuclei suspensions are partitioned into nanoliter-scale Gel Beads-in-emulsion (GEMs). A pool of ~737,000 10x GEM Barcodes (also referred to as 10x Barcodes) is sampled separately to index the contents of each partition. Inside the GEMs, probes are ligated and the 10x GEM Barcode is added, and all ligated probes within a GEM share a common 10x GEM Barcode. Barcoded and ligated probes are then pre-amplified in bulk, after which gene expression libraries are generated and sequenced. The protocol also enables the detection of cell surface proteins in fixed single cell suspensions that are labeled with antibody-oligonucleotide conjugates prior to fixation.

## Sample Fixation for Chromium Fixed RNA Profiling

Samples are fixed in a 4% formaldehyde fixative solution, as described in the Fixation of Cells & Nuclei for Chromium Fixed RNA Profiling (Demonstrated Protocol CG000478) and using Chromium Next GEM Single Cell Fixed RNA Sample Preparation Kit (10x Genomics PN-1000414). Fixed samples are taken through a series of steps to generate Fixed RNA – Gene Expression and Cell Surface Protein libraries. A high-level overview of each step in the protocol, including gene expression and cell surface protein library construction, is provided in the following sections.

### Step 1: Probe Hybridization


The whole transcriptome probe pairs, consisting of a left hand side (LHS) and a right hand side (RHS) for each targeted gene, are added to the fixed sample. Together, probe pairs hybridize to their complementary target RNA in an overnight incubation.

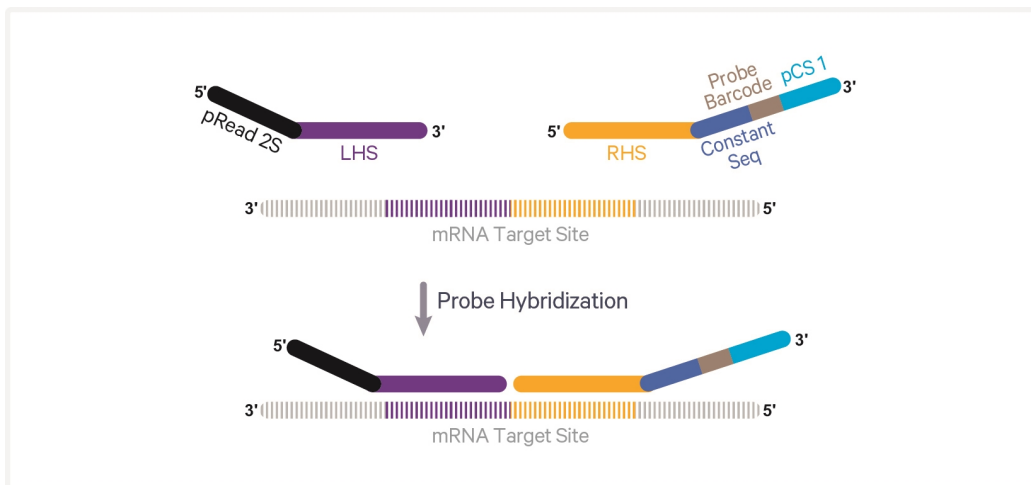

## Step 2: GEM Generation & Barcoding


After hybridization, the unbound probes are washed off. GEMs are generated by combining barcoded Gel Beads, a Master Mix containing cells, and Partitioning Oil onto Chromium Next GEM Chip Q.

### GEM Generation

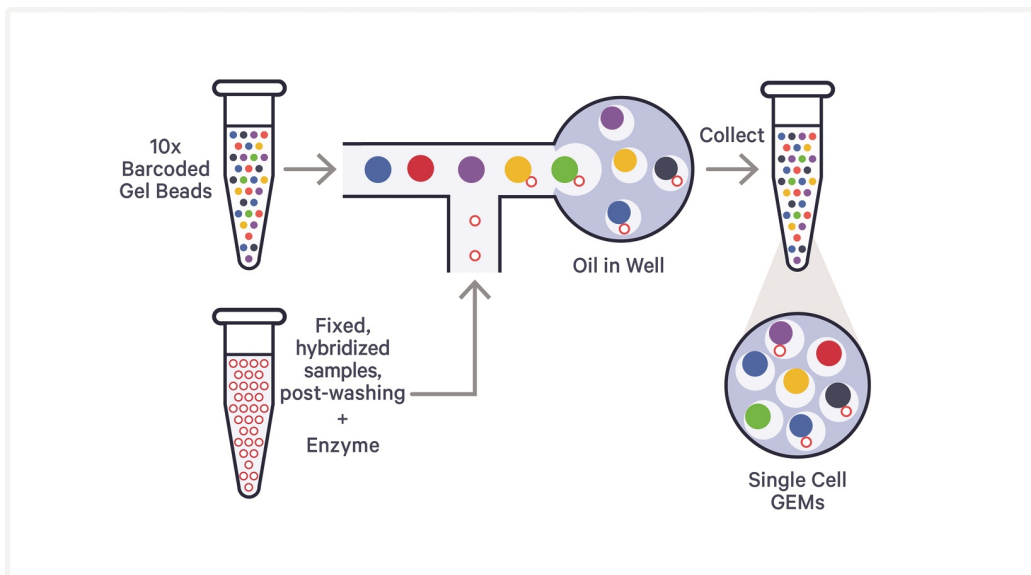

Immediately following GEM generation, the Gel Bead is dissolved, releasing the barcoded Gel Bead primers, and any co-partitioned cell is lysed. Gel Bead primers contain an Illumina TruSeq Read 1 sequence (Read 1T, read 1 sequencing primer), a 16 nt 10x GEM Barcode (or 10x Barcode), a 12 nt unique molecular identifier (UMI), and partial Capture Sequence 1 (sequence complementary to the probe and Feature Barcode oligonucleotide).

### Gel Bead

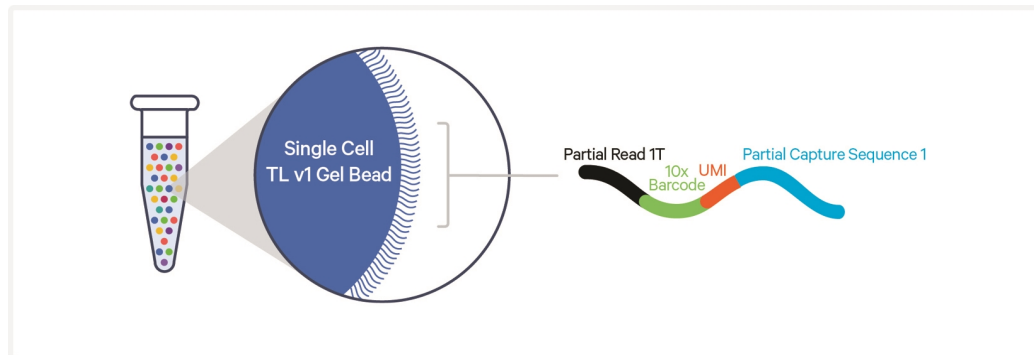


After GEM generation, the partitioned cells, Gel Beads, and Master Mix are placed in a thermal cycler and taken through several steps. First, a ligation step seals the nick between the left hand and right hand probe, while the probes remain hybridized to their target RNA. Second, the Gel Bead primer hybridizes to the capture sequence on the ligated probe pair and is extended by a polymerase to add the UMI, 10x GEM Barcode, and partial Read 1T.

If the sample was previously labeled with antibody-oligonucleotide conjugates, the Gel Bead primer captures the cell surface protein Feature Barcode conjugated to the antibody containing (i) a TruSeq Read 2 (Read 2T), (ii) a 15 nt Feature Barcode, and (iii) Capture Sequence 1. Incubation of the GEMs with the Master Mix, produces 10x Barcoded DNA from the cell surface protein Feature Barcode.


Finally, a heat denaturation step inactivates the enzymes in the GEM reaction.



**Inside Individual GEMs**

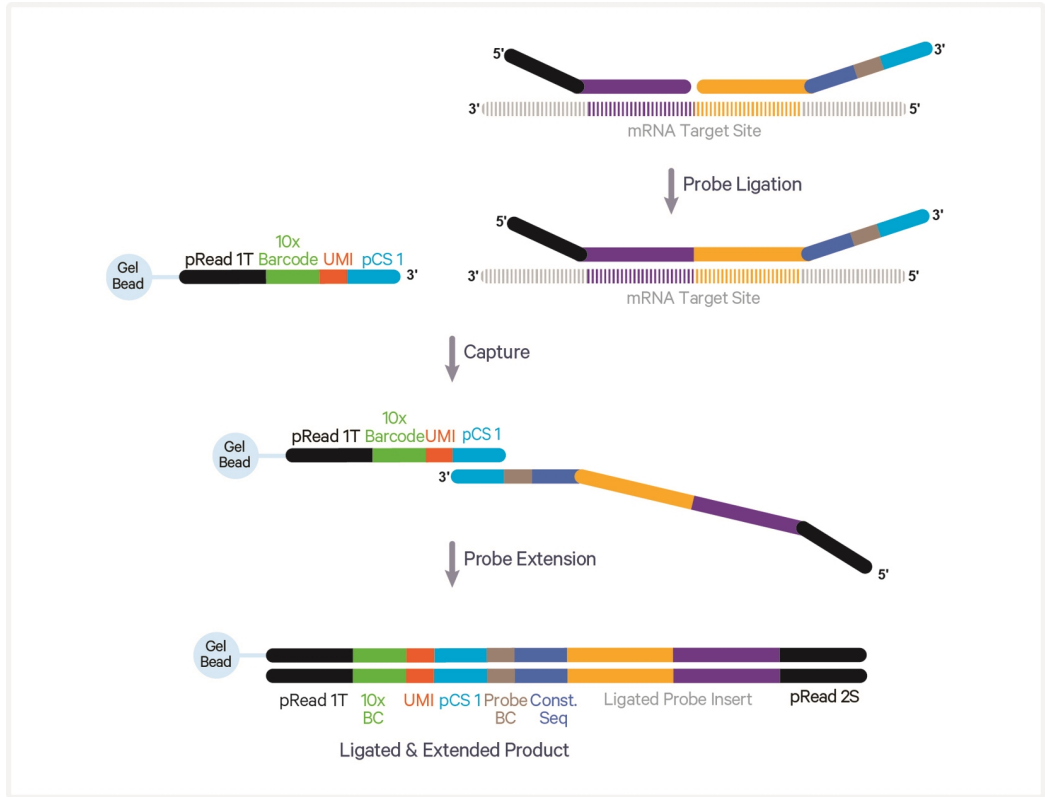

**Inside Individual GEMs - If cells were also labeled with antibody-oligonucleotide conjugates**

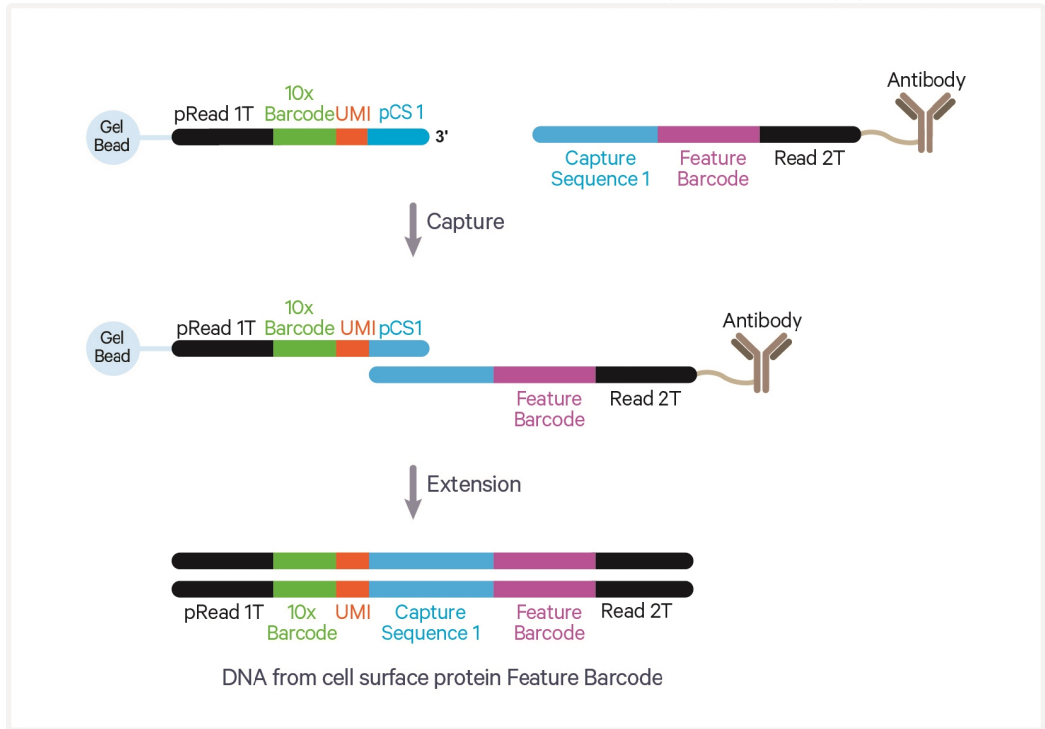

### Step 3: GEM Recovery & Pre-Amplification

Once the ligation and barcoding steps are completed, the GEMs are broken by the addition of Recovery Agent, inverting the mixture, and removing the Recovery Agent. A PCR master mix is added directly to the post-GEM aqueous phase to pre-amplify the ligated product, along with (optionally) barcoded Feature Barcodes identifying cell surface proteins. The pre-amplified products are then cleaned up by SPRIselect.

#### DNA Pre-Amplification




### Step 4: Fixed RNA – Gene Expression Library Construction

The 10x barcoded, ligated probe products undergo indexing via Sample Index PCR. This, in turn, generates final library molecules that are cleaned up by SPRIselect, assessed on a bioanalyzer or a similar instrument, quantified, and then sequenced.

P5, P7, i5 and i7 sample indexes, and Illumina Small Read 2 (Read 2S) sequences are added via Sample Index PCR. The final libraries contain the P5 and P7 priming sites used in Illumina sequencers.

**Pooled Amplified DNA Processed in Bulk**

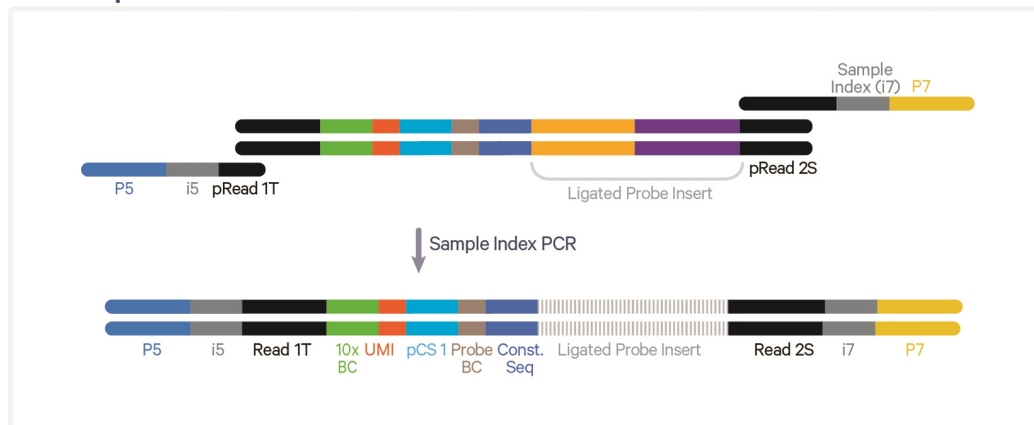


## Step 5: Fixed RNA – Cell Surface Protein Library Construction

The amplified DNA from cell surface protein Feature Barcode undergo indexing via Sample Index PCR. This, in turn, generates final library molecules that are cleaned up by SPRIselect, assessed on a bioanalyzer or a similar instrument, quantified, and then sequenced.

P5, P7, i5, and i7 sample indexes are added via Sample Index PCR. The final libraries contain the P5 and P7 priming sites used in Illumina sequencers.


**Pooled Amplified DNA Processed in Bulk**

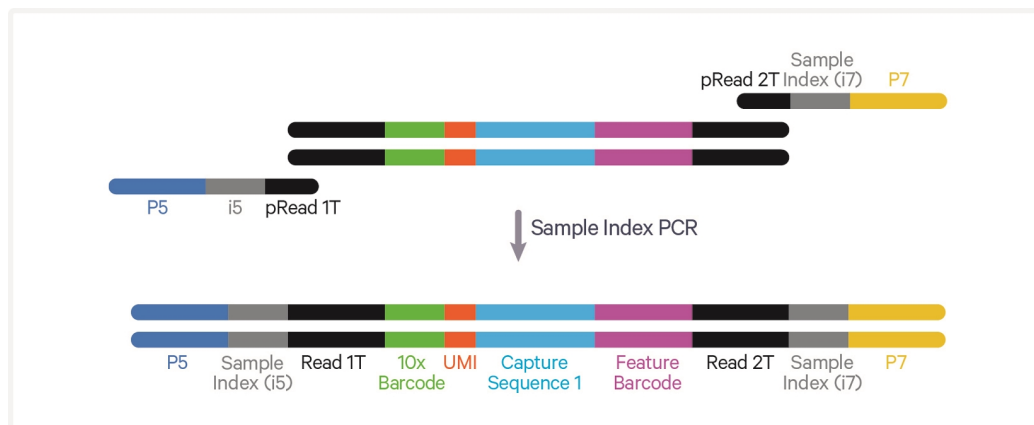

## Step 6: Sequencing


A Chromium Fixed RNA Profiling – Gene Expression library comprises standard Illumina paired-end constructs which begin and end with P5 and P7. The 16 bp 10x GEM Barcode and 12 bp UMI are encoded in Read 1T. Small RNA Read 2 (Read 2S) sequences the ligated probe insert.

### Chromium Fixed RNA Profiling – Gene Expression Library

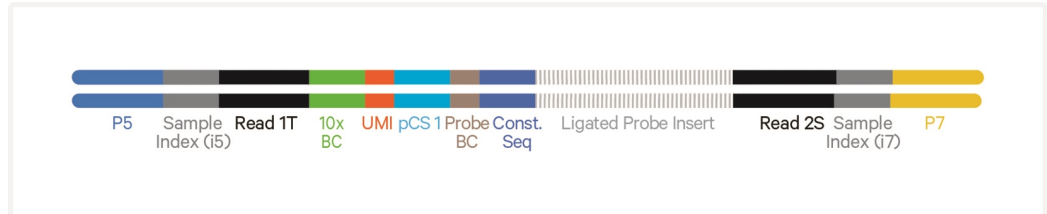

A Chromium Fixed RNA Profiling – Cell Surface Protein library comprises standard Illumina paired-end constructs which begin and end with P5 and P7. The 16 bp 10x GEM Barcode (10x Barcode) and 12 bp UMI are encoded in Read 1T. Illumina TruSeq Read 2 (Read 2T) sequences the Feature Barcode.

### Chromium Fixed RNA Profiling – Cell Surface Protein Library

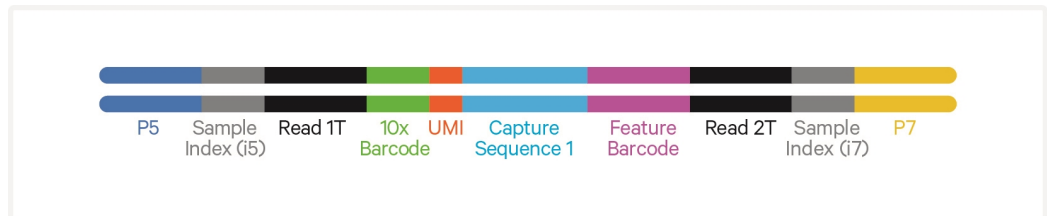

See Appendix for [Oligonucleotide Sequences on page 91](#)

## Sample Fixation Guidelines

### Overview



Chromium Fixed RNA Profiling is compatible with samples (single cell and nuclei suspensions) that are appropriately fixed. Samples are fixed in a 4% formaldehyde fixative solution and then quenched as described in the Fixation of Cells & Nuclei for Chromium Fixed RNA Profiling (Demonstrated Protocol CG000478) and Tissue Fixation & Dissociation for Chromium Fixed RNA Profiling (CG000553). Chromium Next GEM Single Cell Fixed RNA Sample Preparation Kit (10x Genomics PN-1000414) is required for fixation.

Prior to fixation, cells can also be labeled using a specific protein binding molecule, such as an antibody conjugated to a Feature Barcode oligonucleotide. Consult Demonstrated Protocol Cell Surface Protein Labeling for Chromium Fixed RNA Profiling for Singleplexed Samples with Feature Barcode technology (CG000529) for guidance.

Fixed samples can be stored short-term at 4°C or long-term (up to 6 months for fixed cells or nuclei) at -20°C or -80°C. Fixation and quenching protocols as well as the storage guidelines are described in detail in the Fixation of Cells & Nuclei for Chromium Fixed RNA Profiling (Demonstrated Protocol CG000478) and Tissue Fixation & Dissociation for Chromium Fixed RNA Profiling (CG000553).



Failure to label cell surface proteins with a Feature Barcode oligonucleotide conjugated to a specific protein binding molecule prior to fixing the cells will preclude generation of Cell Surface Protein library.

#### Key Reagents for Sample Fixation

Vendor	Description	Part Number
10x Genomics	Chromium Next GEM Single Cell Fixed RNA Sample Preparation Kit, 16 rxns	1000414
Thermo Fisher Scientific	Formaldehyde (37% by Weight/Molecular Biology), Fisher BioReagents	BP531-500

*Consult Fixation of Cells & Nuclei for Chromium Fixed RNA Profiling (Demonstrated Protocol CG000478) and Tissue Fixation & Dissociation for Chromium Fixed RNA Profiling (CG000553) for a complete list of reagents for sample fixation.*

# Tips & Best Practices



## Icons



Tips & Best Practices section includes additional guidance



Signifies critical step requiring accurate execution



Troubleshooting section includes additional guidance



Chip Q specific steps

## Emulsion-safe Plastics

Use validated emulsion-safe plastic consumables when handling GEMs as some plastics can destabilize GEMs.

## General Reagent Handling

- Fully thaw the reagents at indicated temperatures. Thoroughly mix reagents before use.
- Keep all enzymes and Master Mixes at indicated temperatures during setup and use. Promptly move reagents back to the recommended storage.
- Calculate reagent volumes with indicated % excess of 1 reaction values.
- Cover Partitioning Oil tubes and reservoirs to minimize evaporation.
- Thoroughly mix samples with the beads during bead-based cleanup steps.

## Pipette Calibration

- Follow manufacturer's calibration and maintenance schedules.
- Pipette accuracy is particularly important when using SPRIselect reagents.

## Probe Hybridization

### Cell Counts

- The minimum input is 200,000 cells or 400,000 nuclei per hybridization and the maximum input is  $2 \times 10^6$  for cells/nuclei per hybridization.
- During post-hybridization washing steps, some cell loss is expected. It is recommended to start the hybridization reaction with  $\sim 1 \times 10^6$  cells/nuclei per hybridization, if possible.
- If proceeding with <500,000 fixed cells in a hybridization reaction, use a swinging bucket rotor for centrifugation and carefully remove the


supernatant without disturbing the pellet. In such cases, complete removal of the supernatant is not required. Up to 15  $\mu\text{l}$  of supernatant may be left behind prior to resuspending the cell pellet in the Hyb Mix to optimize cell recovery without significantly impacting assay performance.

### Incubation Time

- Recommended incubation time for probe hybridization is 16-24 h.
- Incubation time should be of same length for all samples. DO NOT mix samples with different hybridization times in one experiment.

## Sample Washing & Recovery

- Swinging bucket centrifuge can increase cell recovery during washing.
- When performing post-hybridization washing with low cell numbers (i.e. <500,000 cells), complete removal of the supernatant is not required. Up to 30  $\mu\text{l}$  of supernatant may be left behind to optimize cell recovery without significantly impacting assay performance.

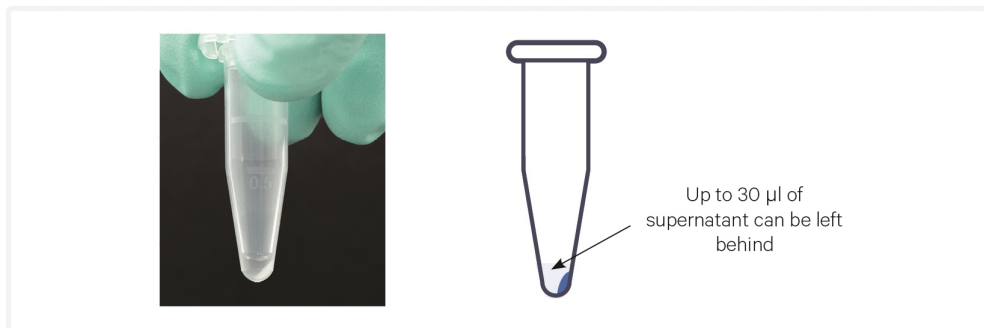

## Cell Counts for Chip Loading

- The Chromium Fixed RNA Profiling is designed to target 500-10,000 cells per sample with a per sample undetected multiplet rate of 0.4% to 8.0%.
- Recommended starting point is to target ~4,000 cells, and a multiplet rate of ~3.2%.
- The minimum cell input concentration to get maximum cell recovery is 413 cells/ $\mu\text{l}$ .



Undetectable Multiplet Rate (%)	# of Cells Loaded	# of Cells Recovered
~0.4	825	500
~0.8	1,650	1,000
~1.6	3,300	2,000
~2.4	4,950	3,000
~3.2	6,600	4,000
~4.0	8,250	5,000
~4.8	9,900	6,000
~5.6	11,550	7,000
~6.4	13,200	8,000
~7.2	14,850	9,000
~8.0	16,500	10,000


## Sample Filtration

- After post-hybridization wash, pass the sample through a 30 µm filter (Sysmex CellTrics or Miltenyi Biotec Pre-Separation) into a new 1.5-ml microcentrifuge tube.



DO NOT use 40 µm Flowmi Tip Strainer for filtration.

- Hold the pipette tip at an angle and touch the filter membrane where the filter meets the wall. Slowly pipette through the filter. Tap gently or centrifuge briefly if liquid remains at the end of the filter.
- To maximize recovery, residual volume can be pipetted from underneath the filter.

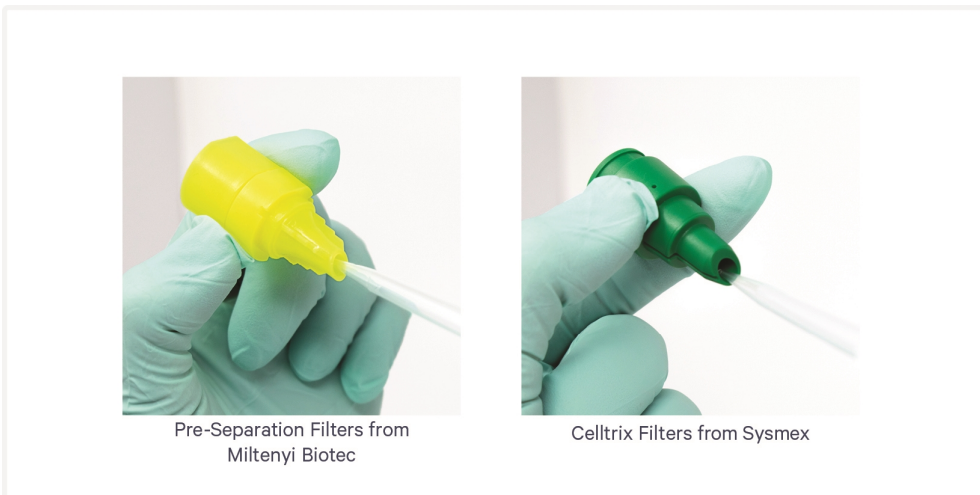


## Cell Counting

- Accurate counting is critical for optimal assay performance.
- It is strongly recommended that the sample be stained with a fluorescent dye such as Ethidium Homodimer-1 or PI staining solution and counted using an automated fluorescent cell counter (Countess II Automated Cell Counter or a Cellaca Counter).
- If using an automated cell counter, ensure that the cells are being circled correctly. The settings of the automated cell counters may need to be adjusted for optimal cell-detection accuracy

### Counting using Ethidium Homodimer-1

This protocol provides instructions for counting sample using Ethidium Homodimer-1 and the Countess II FL Automated Cell Counter (with RFP light cube) to enable accurate quantification even in the presence of sub-cellular debris. The optimal cell concentration for the Countess is 1,000-4,000 cells/ $\mu$ l. Refer to manufacturer's instructions for details on operations.

- Vortex Ethidium Homodimer-1, centrifuge briefly, and dilute the concentrated stock as per manufacturer's instructions (~1:100 dilution).
- Aliquot 10  $\mu$ l diluted Ethidium Homodimer-1 in each tube.
- Gently mix the sample. Immediately add 10  $\mu$ l sample to 10  $\mu$ l diluted Ethidium Homodimer-1. Gently pipette mix 10x.
- Transfer 10  $\mu$ l sample to a Countess II Cell Counting Slide chamber.
- Insert the slide into the Countess II FL Cell Counter. Image the sample using the RFP setting for fluorescent illumination and filtering. Optimize focus and exposure settings and confirm the absence of large clumps using the bright-field mode. Make sure the Countess is circling RFP positive cells. Note the RFP-positive concentration. Multiply by dilution factor 2 to determine cell concentration.

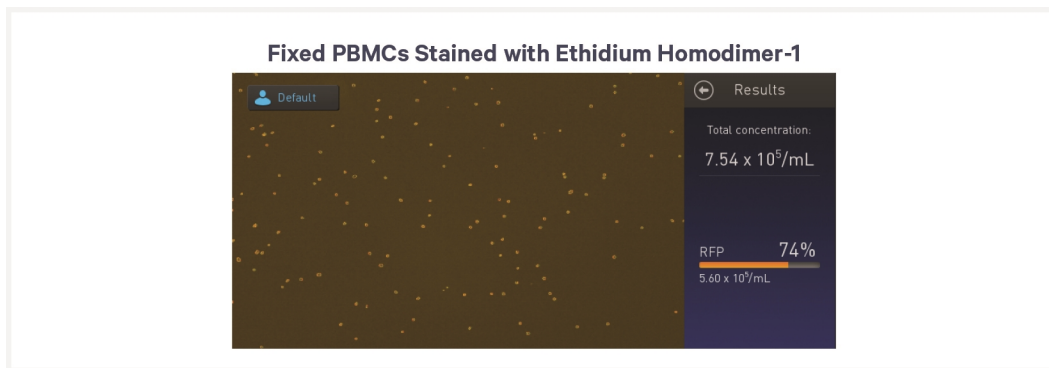

Samples stained with Ethidium Homodimer-1 can also be counted using Cellca counter. See manufacturer's instructions for details.

## Counting using PI Staining Solution

This protocol provides instructions for counting sample using PI staining solution and the Cellca Counter to enable accurate quantification even in the presence of sub-cellular debris. The optimal cell concentration for the Cellca Counter is 100-10,000 cells/ $\mu$ L. Refer to manufacturer's instructions for details on operations.

- Add 25  $\mu$ L PI Staining Solution into Mixing Row of Cellca plate.
- Gently mix the sample. If the sample is too concentrated, a 1:1 dilution in PBS can also be prepared. For example, add 15  $\mu$ L fixed cell suspension to 15  $\mu$ L PBS.
- Add 25  $\mu$ L sample to Mixing Row of plate containing PI Staining Solution. Gently pipette mix 8x.
- Transfer stained sample to Loading Row of Cellca plate.
- For counting fixed samples, only use the PI (Propidium Iodide) channel. Refer to manufacturer's instructions for details.

Samples stained with PI Staining Solution can also be counted using Countess II FL Automated Cell Counter. See manufacturer's instructions for details.


## Counting using Trypan Blue (Only for Debris-Free Samples)

Debris-free samples (cells or nuclei suspensions) can also be counted using trypan blue. This protocol provides instructions for counting sample using trypan blue and a hemocytometer or Countess II Automated Cell Counter.

- Mix 1 part 0.4% trypan blue and 1 part sample.

- Transfer 10  $\mu$ l sample to a Countess II Cell Counting Slide chamber or a hemocytometer.
- Insert the slide into the Countess II Cell Counter and determine the cell concentration. Or if using hemocytometer, count fixed cells by placing hemocytometer under the microscope.
- The majority of fixed cells or nuclei suspensions will be stained with trypan blue stain and appear non-viable.

#### Brightfield image from Automated Cell Counter – Samples Stained with Trypan Blue

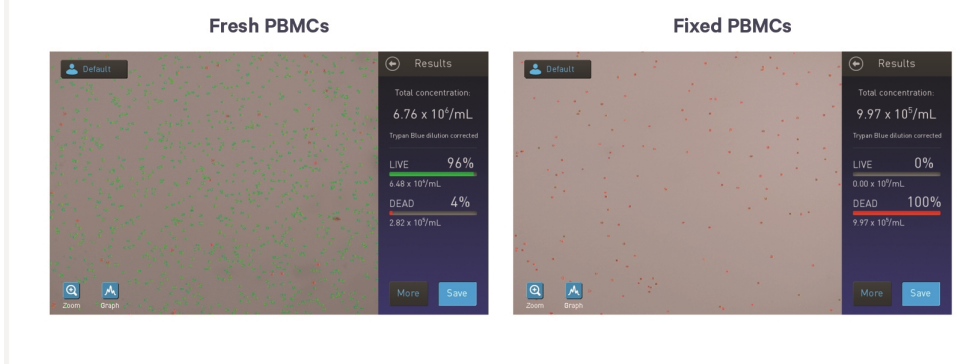



## Chromium Next GEM Chip Handling

- Chromium Fixed RNA Profiling uses Chromium Chip Q.
- Minimize exposure of reagents, chips, and gaskets to sources of particles and fibers, laboratory wipes, frequently opened flip-cap tubes, clothing that sheds fibers, and dusty surfaces.
- After removing the chip from the sealed bag, use in **≤24 h**.
- Execute steps without pause or delay, unless indicated. When using multiple chips, load, run, and collect the content from one chip before loading the next.
- Fill all unused input wells in rows labeled 1, 2, and 3 on a chip with an appropriate volume of 50% glycerol solution before loading the used wells. **DO NOT** add glycerol to the wells in the bottom **NO FILL** row.
- Avoid contacting the bottom surface of the chip with gloved hands and other surfaces. Frictional charging can lead to inadequate priming of the channels, potentially leading to either clogs or wetting failures.
- Minimize the distance that a loaded chip is moved to reach the Chromium X/iX.
- Keep the chip horizontal to prevent wetting the gasket with oil, which depletes the input volume and may adversely affect the quality of the resulting emulsion.


## Chromium Next GEM Secondary Holders



- Chromium Next GEM Secondary Chip Holders encase Chromium Next GEM Chips.
- The holder lid flips over to become a stand, holding the chip at 45 degrees for optimal recovery from each well.
- Squeeze the black sliders on the back side of the holder together to unlock the lid and return the holder to a flat position.

## Chromium Next GEM Chip & Holder Assembly with Gasket

- Chromium Next GEM Chip Q is only compatible with Chromium Next GEM Secondary Holder (PN-3000332). DO NOT use any other holder
- Close the holder lid. Attach the gasket by holding the tongue (curved end, to the right) and hook the gasket on the left-hand tabs of the holder. Gently pull the gasket toward the right and hook it on the two right-hand tabs.
- DO NOT touch the smooth side of the gasket.
- Open the chip holder.
- Align notch on the chip (upper left corner) and the open holder with the gasket attached.
- Slide the chip to the left until the chip is inserted under the guide on the holder. Depress the right hand side of the chip until the spring-loaded clip engages.
- Keep the assembled unit with the attached gasket until ready for dispensing reagents into the wells.

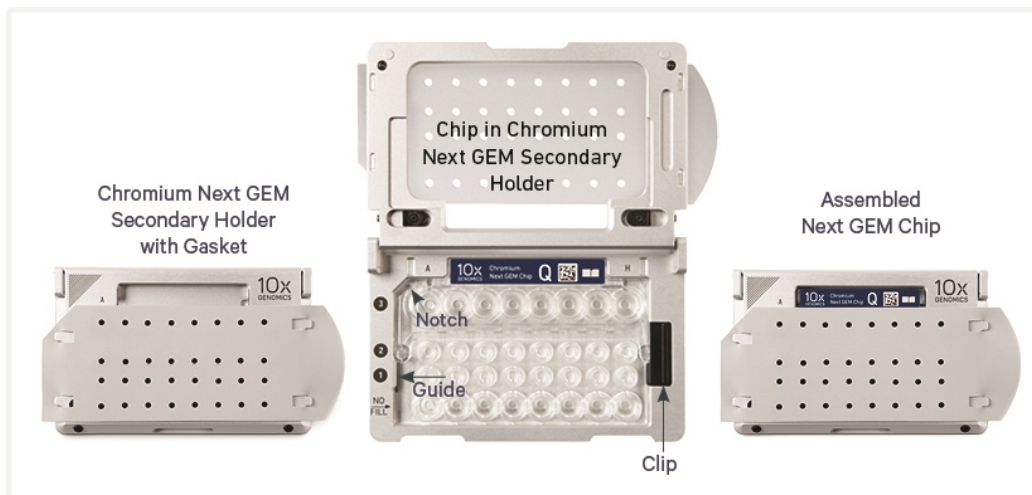



## Chromium Next GEM Chip Loading

- Place the assembled chip and holder flat (gasket attached) on the bench with the lid open.
- Dispense at the bottom of the wells without introducing bubbles.
- When dispensing Gel Beads into the chip, wait for the remainder to drain into the bottom of the pipette tips and dispense again to ensure complete transfer.



- Refer to [2.3 Load Chromium Next GEM Chip Q](#) on page 54 for specific instructions.

## 50% Glycerol Solution for Addition to Unused Chip Wells


- Purchase 50% glycerol solution from Ricca Chemical Company, Glycerin (glycerol), 50% (v/v) Aqueous Solution, PN-3290-32.  
OR
- Prepare 50% glycerol solution:
  - Mix an equal volume of water and  $\geq 99\%$  Glycerol, Molecular Biology Grade.
  - Filter through a 0.2  $\mu\text{m}$  filter.
  - Store at  $-20^{\circ}\text{C}$  in 1-ml LoBind tubes. 50% glycerol solution should be equilibrated to room temperature before use.
- Adding glycerol to non-sample chip wells is essential to avoid chip failure.



## 50% Glycerol Solution for Sample Storage


- Use nuclease-free water and molecular biology grade Glycerol from Millipore Sigma, PN-G5516, to prepare 50% glycerol solution as described previously. DO NOT use 50% glycerol solution from Ricca Chemical Company, Glycerin (glycerol), 50% (v/v) Aqueous Solution, PN-3290-32.

## Gel Bead Handling

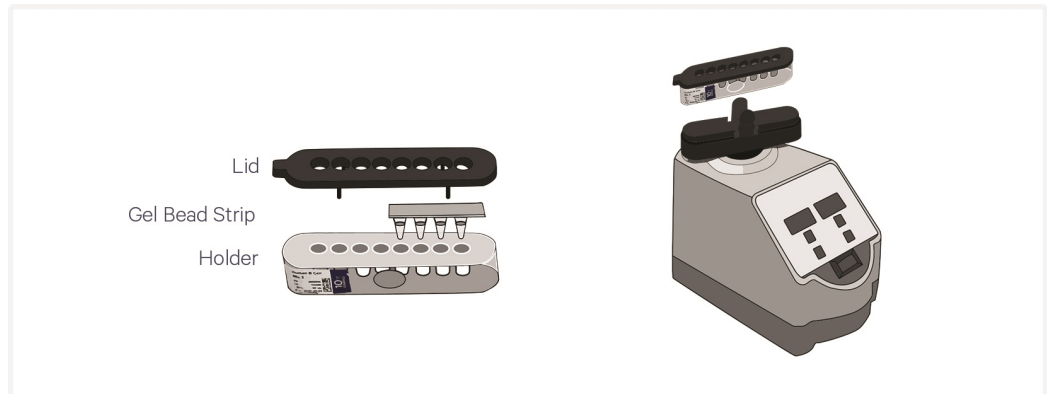

- Use one tube of Gel Beads per sample. DO NOT puncture the foil seals of tubes not used at the time.
- After removing the Gel Bead strip from the packaging, equilibrate the Gel Bead strip to **room temperature** for at least **30 min** before use.
- Store unused Gel Beads at **-80°C** and avoid more than 12 freeze-thaw cycles. DO NOT store Gel Beads at **-20°C**.
- Snap the tube strip holder with the Gel Bead strip into a 10x Vortex Adapter. Vortex **30 sec**.
- Centrifuge the Gel Bead strip for **~5 sec** after removing from the holder. Confirm there are no bubbles at the bottom of tubes and the liquid levels look even. Place Gel Bead strip back in the holder and secure the holder lid.


## 10x Magnetic Separator

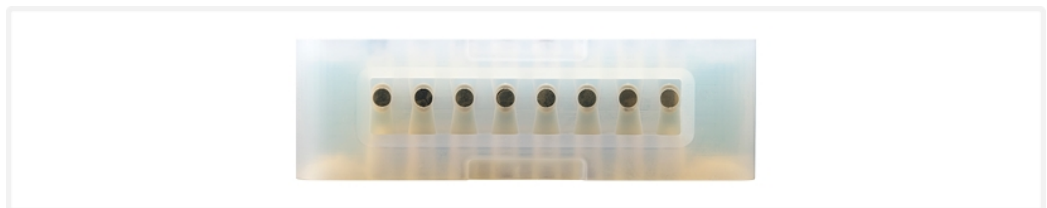

- Offers two positions of the magnets (high and low) relative to a tube, depending on its orientation. Flip the magnetic separator over to switch between high (magnet•**High**) or low (magnet•**Low**) positions.
- If using MicroAmp 8-Tube Strips, use the high position (magnet•**High**) only throughout the protocol.

## Magnetic Bead Cleanup Steps

- During magnetic bead based cleanup steps that specify waiting “until the solution clears”, visually confirm clearing of solution before proceeding to the next step. See panel below for an example.
- The time needed for the solution to clear may vary based on specific step, reagents, volume of reagents etc.
- **Visually Confirm Clearing of Magnetic Bead Solution**

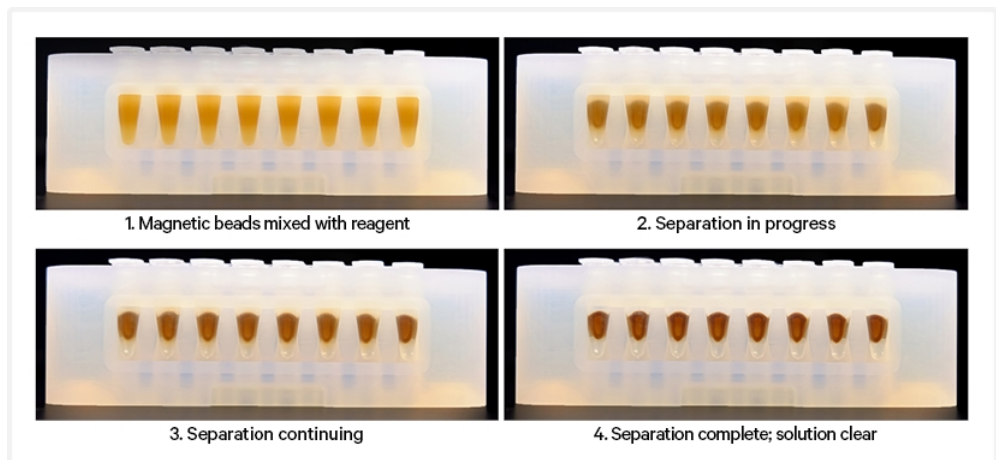

## SPRIselect Cleanup & Size Selection

- After aspirating the desired volume of SPRIselect reagent, examine the pipette tips before dispensing to ensure the correct volume is transferred.
- Pipette mix thoroughly as insufficient mixing of sample and SPRIselect reagent will lead to inconsistent results.
- Use fresh preparations of 80% Ethanol.

### Tutorial — SPRIselect Reagent:DNA Sample Ratios

SPRI beads selectively bind DNA according to the ratio of SPRIselect reagent (beads).

$$\text{Example Ratio: } = \frac{\text{Volume of SPRIselect reagent added to the sample}}{\text{Volume of DNA sample}} = \frac{50 \mu\text{l}}{100 \mu\text{l}} = \mathbf{0.5X}$$

## Sample Indices in Sample Index PCR

- Choose the appropriate sample index sets to ensure that no sample indices overlap in a multiplexed sequencing run.
- Verify and use the specified index plate only. DO NOT use the plates interchangeably.
- Each well in the Dual Index Plate contains a unique i7 and a unique i5 oligonucleotide.







# Step 1:

## Probe Hybridization

1.0 Get Started	41
1.1 Probe Hybridization	42

1

## 1.0 Get Started

Action	Item	10x PN	Preparation & Handling	Storage
<b>Thaw &amp; Keep Warm</b>				
<input type="checkbox"/>	 <b>Hyb Buffer B</b>	2000483	Thaw at 42°C. Vortex and centrifuge briefly. Keep warm and verify no precipitate before use.  DO NOT keep the thawed buffer on ice, or the solution will precipitate. Thawed Hyb Buffer B can be kept at 42°C for up to 1 h.	-20°C
<input type="checkbox"/>	 <b>Enhancer</b>	2000482	Thaw for 10 min at 65°C. Vortex and centrifuge briefly. Keep warm and verify no precipitate before use.  DO NOT keep the thawed reagent on ice, or the solution will precipitate. Once thawed, Enhancer can be kept at 42°C for up to 10 min.	-20°C
<b>Place on Ice</b>				
<input type="checkbox"/>	 <b>Fixed Cell Suspension</b>	—	Consult Demonstrated Protocols Fixation of Cells & Nuclei for Chromium Fixed RNA Profiling (CG000478) and Tissue Fixation & Dissociation for Chromium Fixed RNA Profiling (CG000553).	—
<input type="checkbox"/>	 <b>Human WTA Probes BC001</b> OR <b>Mouse WTA Probes BC001</b>	2000495 OR 2000703	Thaw on ice. Vortex and centrifuge briefly.	-20°C

## 1.1 Probe Hybridization



Before starting this protocol, ensure that samples have been appropriately fixed and quenched. Consult Fixation of Cells & Nuclei for Chromium Fixed RNA Profiling (Demonstrated Protocol CG000478) for details.

- a. Set a thermomixer with heated lid to 42°C or prepare a thermal cycler with the following incubation protocol and start the program.

Lid Temperature	Reaction Volume	Run Time
42°C	100 µl	Overnight
Step	Temperature	Time
Pre-equilibrate	42°C	Hold
Probe Hybridization	42°C	16-24 h

- b. Prepare Hyb Mix at **room temperature**. Pipette mix 10x.

Hyb Mix	PN	1X* (µl)	1X* + 20% (µl)	4X* + 20% (µl)
<i>Add reagents in the order listed</i>				
● <b>Hyb Buffer B</b> <i>Thaw at 42°C. Add warm to the mix and if appears milky keep it back on 42°C.</i>	2000483	70.0	84.0	336.0
● <b>Enhancer</b> <i>Heat at 65°C for 10 min. Vortex and verify no precipitate. Add warm to the mix.</i>	2000482	10.0	12.0	48.0
<b>Total</b>	-	<b>80.0</b>	<b>96.0</b>	<b>384.0</b>

\*1X = 1 fixed sample, 4X = 4 fixed samples

- c. Incubate Hyb Mix at **42°C** for **5 min**.
- d. Centrifuge fixed cells/nuclei resuspended in Quenching Buffer/post-storage processing buffer (0.5X PBS + 0.02% BSA) at **850 rcf** for **5 min** at **4°C**. The following table provides guidelines for number of cells/nuclei recommended during hybridization.

Sample Type	Recommended #
Fixed cells	<b>200,000-2 x 10<sup>6</sup></b>
Fixed nuclei	<b>400,000-2 x 10<sup>6</sup></b>
FFPE dissociated suspension	<b>400,000-2 x 10<sup>6</sup></b>



*DO NOT exceed  $2 \times 10^6$  cells in one hybridization reaction.*

- e.** Remove the supernatant.

*If proceeding with <math>500,000</math> fixed cells, use a swinging bucket rotor for centrifugation and carefully remove the supernatant without disturbing the pellet. In such cases, complete removal of the supernatant is not required. Up to 15  $\mu\text{l}$  of supernatant may be left behind to optimize cell recovery without significantly impacting assay performance.*



- f.** Resuspend each pellet in **80  $\mu\text{l}$**  Hyb Mix. Keep sample at **room temperature**. DO NOT place on ice.

- g.** Add **20  $\mu\text{l}$**  Human/Mouse WTA Probes BC001 to the **80  $\mu\text{l}$**  mixture of Hyb Mix and fixed sample and gently pipette mix 10x with pipette set at 80  $\mu\text{l}$ .

- h.** Incubate sample for **16-24 h** at **42°C** in a thermomixer with heated lid and no shaking. If a thermomixer with heated lid is not available, samples can be transferred into 8-tube strips and incubated in a thermocycler.



*Incubation for less than 16 h is not recommended. Incubation time should be consistent across all samples in an experiment.*

# Step 2:

## GEM Generation and Barcoding

2.0 Get Started	45
2.1 Post-Hybridization Wash	47
2.2 Prepare GEM Master Mix + Sample Dilution	50
2.3 Load Chromium Next GEM Chip Q	54
2.4 Run the Chromium X/iX	56
2.5 Transfer GEMs	57
2.6 GEM Incubation	58












## 2.0 Get Started



Firmware Version 1.1.0 or higher is required in the Chromium X/iX used for this Chromium Fixed RNA Profiling protocol.



Action	Item	10x PN	Preparation & Handling	Storage
<b>Equilibrate to Room Temperature</b>				
<input type="checkbox"/>	<b>Single Cell TL v1 Gel Beads</b>	2000538	Equilibrate to room temperature 30 min before loading the chip.	-80°C
<input type="checkbox"/>	 <b>Reducing Agent B</b>	2000087	Vortex, verify no precipitate, centrifuge briefly.	-20°C
<b>Thaw &amp; Keep Warm</b>				
<input type="checkbox"/>	 <b>Enhancer</b>	2000482	Thaw for 10 min at 65°C. Vortex and centrifuge briefly. Keep warm and verify no precipitate before use. <div style="margin-top: 10px;">                      DO NOT keep the thawed reagent on ice, or the solution will precipitate. Once thawed, Enhancer can be kept at 42°C for up to 10 min.                 </div>	-20°C
<b>Place on Ice</b>				
<input type="checkbox"/>	 <b>Conc. Post-Hyb Buffer</b>	2000533	Thaw at room temperature and keep on ice.	-20°C
<input type="checkbox"/>	 <b>GEM Enzyme Mix</b>	2000490	Centrifuge briefly before adding to the mix.	-20°C
<input type="checkbox"/>	 <b>GEM Reagent Mix</b>	2000491	Thaw at room temperature. Vortex, verify no precipitate, centrifuge briefly. Keep on ice.	-20°C
<b>Obtain</b>				
<input type="checkbox"/>	 <b>Partitioning Oil</b>	2000190	—	Ambient
<input type="checkbox"/>	<b>Next GEM Chip Q</b>	2000518	See Tips & Best Practices.	Ambient
<input type="checkbox"/>	<b>Chromium Next GEM Secondary Chip Holder</b>	3000332	See Tips & Best Practices.	Ambient
<input type="checkbox"/>	<b>10x Gasket</b>	370017/ 3000072	See Tips & Best Practices.	Ambient

Action	Item	10x PN	Preparation & Handling	Storage
<input type="checkbox"/>	<b>Sample Filters</b> <b>Sysmex Sterile Single-pack CellTrics Filters/Miltenyi Biotec Pre-Separation Filters (30 µm)</b>	—	Manufacturer's recommendations.	Ambient
<input type="checkbox"/>	<b>10x Vortex Adapter</b>	330002	See Tips & Best Practices.	Ambient
<input type="checkbox"/>	<b>50% glycerol solution</b> for adding to unused wells	—	See Tips & Best Practices.	
<input type="checkbox"/>	<b>Glycerol for molecular biology, ≥99%</b> Prepare 50% glycerol solution for sample storage	—	See Tips & Best Practices.	—

## 2.1 Post-Hybridization Wash



- a. Prepare Post-Hyb Wash Buffer. Vortex briefly and keep at **room temperature**. DO NOT keep at 4°C.

Post-Hyb Wash Buffer	PN	1X + 10% (ml)*	4X + 10% (ml)*
<i>Add reagents in the order listed</i>			
<b>Nuclease-free Water</b>	-	1.98	7.92
 <b>Conc. Post-Hyb Buffer</b>	2000533	0.11	0.44
 <b>Enhancer</b> <i>Heat at 65°C for 10 min. Vortex and verify no precipitate. Add warm to the mix.</i>	2000482	0.11	0.44
<b>Total</b>	-	<b>2.2</b>	<b>8.8</b>

*\*Volumes are in ml*

- b. Remove tubes from thermomixer (1.5-ml microcentrifuge tubes) after overnight incubation. If hybridization was performed in 8-tube strips, remove tubes from thermal cycler.
- c. Add **900 µl** Post-Hyb Wash Buffer to the sample in 1.5-ml microcentrifuge tube. Pipette mix 5x.

*If the hybridization was performed in 8-tube strips, add 175 µl Post-Hyb Wash Buffer to the sample, gently pipette mix, and transfer to a 1.5-ml microcentrifuge tube. Wash the tube strips with additional Post-Hyb Wash Buffer, transfer to the microcentrifuge tube, and add the remaining volume of Post-Hyb Wash Buffer for a total of 900 µl Post-Hyb Wash Buffer to the sample.*


- d. Incubate at **42°C** for **10 min** in a thermomixer or a heat block.
- e. Centrifuge at **850 rcf** for **5 min** at **room temperature**.
- f. Remove the supernatant without disturbing the pellet.



*See Tips & Best Practices for [Sample Washing & Recovery on page 28](#).*

- g. Resuspend cell pellet in **0.5 ml** Post-Hyb Wash Buffer. Pipette mix 5x.
- h. Incubate at **42°C** for **10 min** in a thermomixer or a heat block.
- i. Centrifuge at **850 rcf** for **5 min** at **room temperature**.
- j. Remove the supernatant without disturbing the pellet.
- k. Resuspend cell pellet in **0.5 ml** Post-Hyb Wash Buffer. Pipette mix 5x.

- l.** Incubate sample at **42°C** for **10 min** in a thermomixer or a heat block.
- m.** Prepare Post-Hyb Resuspension Buffer. Pipette mix 10x and maintain at **4°C**.

Post-Hyb Resuspension Buffer		PN	1X + 10% ( $\mu$ l)	4X + 10% ( $\mu$ l)
<i>Add reagents in the order listed</i>				
	<b>Nuclease-free Water</b>	-	522.5	2090.0
	<b>Conc. Post-Hyb Buffer</b>	2000533	27.5	110.0
	<b>Total</b>	-	<b>550.0</b>	<b>2200.0</b>

- n.** Centrifuge the sample at **850 rcf** for **5 min** at **room temperature**.
- o.** Remove the supernatant without disturbing the pellet.
- p.** Resuspend cell pellet in **500  $\mu$ l chilled** Post-Hyb Resuspension Buffer. Pipette mix 20x to resuspend and breakup any cell clumps and maintain on ice.

*500  $\mu$ l is the recommended volume for most samples. The volume might need adjustment if the starting cell number <200,000 and if >50% cell loss is observed during washing.*

- q.** Pass the sample through a 30  $\mu$ m filter (Sysmex CellTrics or Miltenyi Biotec Pre-Separation Filters) into a new 1.5-ml microcentrifuge tube and place on ice.

DO NOT use 40  $\mu$ m Flowmi Tip Strainer for filtration.



*Filtration is essential for optimal microfluidic performance. Hold the pipette tip at an angle and touch the filter membrane where the filter meets the wall. Slowly pipette through the filter. Tap gently or centrifuge briefly if liquid remains at the end of the filter. To maximize recovery, residual volume can be pipetted from underneath the filter.*



See [Sample Filtration on page 30](#) for details.

- r.** Determine cell concentration of the sample using a Countess II Automated Cell Counter, a Cellaca counter, or a hemocytometer.



See Tips & Best Practices for [Cell Counting on page 31](#). A serial dilution may be needed to accurately determine cell concentration.

*If the sample concentration is not sufficient to achieve the desired targeted cell recovery, concentrate the sample as follows:*

- *Centrifuge a known volume of sample at 850 rcf for 5 min at room temperature.*

- Carefully remove only a fraction of the supernatant, and pipette thoroughly to resuspend the cell pellet in the remaining volume. The amount of supernatant removed should be proportional to the desired increase in concentration.

*For example, to increase the concentration 4-fold from a starting volume of 400  $\mu$ l, centrifuge, then remove 300  $\mu$ l supernatant, and finally resuspend the cell pellet in the remaining 100  $\mu$ l ( $400/100 = 4$ ).*

- Recount to confirm final concentration.



- S.** Store the sample (see Sample Storage below) at **-20°C** or **-80°C** for up to **1 month** or proceed immediately to the next step. If directly proceeding with next step, the undiluted samples can be placed on ice and then stored later after GEM Incubation. See [Prepare GEM Master Mix + Sample Dilution on page 50](#) for details.

### Sample Storage

- Add 0.1 volume Enhancer to sample in Post-Hyb Resuspension Buffer. For example, add 50  $\mu$ l Enhancer to 500  $\mu$ l of sample in Post-Hyb Resuspension Buffer.
- Add 50% glycerol for a final concentration of 10%. For example, add 137.5  $\mu$ l 50% glycerol to 550  $\mu$ l sample in Post-Hyb Resuspension Buffer and Enhancer.
- Store at **-20°C** or **-80°C** for up to **1 month**. For best results, storage at **-80°C** is strongly recommended.

### Using Stored Samples

- When ready to use samples stored at **-20°C** or **-80°C** from this step, thaw at room temperature until no ice remains and then continue from **2.1m** (Prepare Post-Hyb Resuspension Buffer) to wash the sample once before proceeding to the step 2.2. Samples may undergo a color change during storage (e.g. black, light gray, or green), however this will not impact assay performance.



## 2.2 Prepare GEM Master Mix + Sample Dilution

- a. Prepare Master Mix on ice. Pipette mix 15x and centrifuge briefly.

GEM Master Mix		PN	1X ( $\mu$ l)	4X + 10% ( $\mu$ l)
<i>Add reagents in the order listed</i>				
	GEM Reagent Mix	2000491	20.9	92.1
	Reducing Agent B	2000087	1.7	7.3
	GEM Enzyme Mix	2000490	12.4	54.6
<b>Total</b>		-	<b>35.0</b>	<b>154.0</b>

- b. Add the appropriate volume of Post-Hyb Resuspension Buffer to the appropriate volume of sample into each tube of a PCR 8-tube strip on ice. Refer to the Cell Suspension Volume Calculator on the next page for the volumes.

Use the Post-Hyb Resuspension Buffer prepared at the previous step (2.1) for sample dilution. Additional buffer can be prepared using the buffer preparation table in step 2.1.



*Place remaining undiluted sample on ice. These samples can be stored later after GEM incubation. Guidelines for storage of remaining samples are provided in step 2.6 [GEM Incubation on page 58](#).*

- c. Add **35  $\mu$ l** of prepared GEM Master Mix into each tube containing diluted sample and **immediately** proceed to the next step.

*Pipette mixing at this step is not required, and will be performed prior to loading into the chip.*


## Cell Suspension Volume Calculator


Volume of Cell Suspension Stock per reaction (µl) | Volume of Post-Hyb Resuspension Buffer per reaction (µl)

Cell Stock Concentration (Cells/µl)	Targeted Cell Recovery										
	500	1000	2000	3000	4000	5000	6000	7000	8000	9000	10000
400	2.1 37.9	4.1 35.9	8.3 31.8	12.4 27.6	16.5 23.5	20.6 19.4	24.8 15.3	28.9 11.1	33.0 7.0	37.1 2.9	n/a
600	1.4 38.6	2.8 37.3	5.5 34.5	8.3 31.8	11.0 29.0	13.8 26.3	16.5 23.5	19.3 20.8	22.0 18.0	24.8 15.3	27.5 12.5
800	1.0 39.0	2.1 37.9	4.1 35.9	6.2 33.8	8.3 31.8	10.3 29.7	12.4 27.6	14.4 25.6	16.5 23.5	18.6 21.4	20.6 19.4
1000	0.8 39.2	1.7 38.4	3.3 36.7	5.0 35.1	6.6 33.4	8.3 31.8	9.9 30.1	11.6 28.5	13.2 26.8	14.9 25.2	16.5 23.5
1200	0.7 39.3	1.4 38.6	2.8 37.3	4.1 35.9	5.5 34.5	6.9 33.1	8.3 31.8	9.6 30.4	11.0 29.0	12.4 27.6	13.8 26.3
1400	0.6 39.4	1.2 38.8	2.4 37.6	3.5 36.5	4.7 35.3	5.9 34.1	7.1 32.9	8.3 31.8	9.4 30.6	10.6 29.4	11.8 28.2
1600	0.5 39.5	1.0 39.0	2.1 37.9	3.1 36.9	4.1 35.9	5.2 34.8	6.2 33.8	7.2 32.8	8.3 31.8	9.3 30.7	10.3 29.7
1800	0.5 39.5	0.9 39.1	1.8 38.2	2.8 37.3	3.7 36.3	4.6 35.4	5.5 34.5	6.4 33.6	7.3 32.7	8.3 31.8	9.2 30.8
2000	0.4 39.6	0.8 39.2	1.7 38.4	2.5 37.5	3.3 36.7	4.1 35.9	5.0 35.1	5.8 34.2	6.6 33.4	7.4 32.6	8.3 31.8
2200	0.4 39.6	0.8 39.3	1.5 38.5	2.3 37.8	3.0 37.0	3.8 36.3	4.5 35.5	5.3 34.8	6.0 34.0	6.8 33.3	7.5 32.5
2400	0.3 39.7	0.7 39.3	1.4 38.6	2.1 37.9	2.8 37.3	3.4 36.6	4.1 35.9	4.8 35.2	5.5 34.5	6.2 33.8	6.9 33.1
2600	0.3 39.7	0.6 39.4	1.3 38.7	1.9 38.1	2.5 37.5	3.2 36.8	3.8 36.2	4.4 35.6	5.1 34.9	5.7 34.3	6.3 33.7
2800	0.3 39.7	0.6 39.4	1.2 38.8	1.8 38.2	2.4 37.6	2.9 37.1	3.5 36.5	4.1 35.9	4.7 35.3	5.3 34.7	5.9 34.1
3000	0.3 39.7	0.6 39.5	1.1 38.9	1.7 38.4	2.2 37.8	2.8 37.3	3.3 36.7	3.9 36.2	4.4 35.6	5.0 35.1	5.5 34.5
3200	0.3 39.7	0.5 39.5	1.0 39.0	1.5 38.5	2.1 37.9	2.6 37.4	3.1 36.9	3.6 36.4	4.1 35.9	4.6 35.4	5.2 34.8
3400	0.2 39.8	0.5 39.5	1.0 39.0	1.5 38.5	1.9 38.1	2.4 37.6	2.9 37.1	3.4 36.6	3.9 36.1	4.4 35.6	4.9 35.1
3600	0.2 39.8	0.5 39.5	0.9 39.1	1.4 38.6	1.8 38.2	2.3 37.7	2.8 37.3	3.2 36.8	3.7 36.3	4.1 35.9	4.6 35.4
3800	0.2 39.8	0.4 39.6	0.9 39.1	1.3 38.7	1.7 38.3	2.2 37.8	2.6 37.4	3.0 37.0	3.5 36.5	3.9 36.1	4.3 35.7
4000	0.2 39.8	0.4 39.6	0.8 39.2	1.2 38.8	1.7 38.4	2.1 37.9	2.5 37.5	2.9 37.1	3.3 36.7	3.7 36.3	4.1 35.9


Yellow boxes Indicate a low transfer volume that may result in higher cell load variability

## Assemble Chromium Next GEM Chip Q

 Chromium Next GEM Chip Q is only compatible with Chromium Next GEM Secondary Holder (PN-3000332). DO NOT use any other holder.

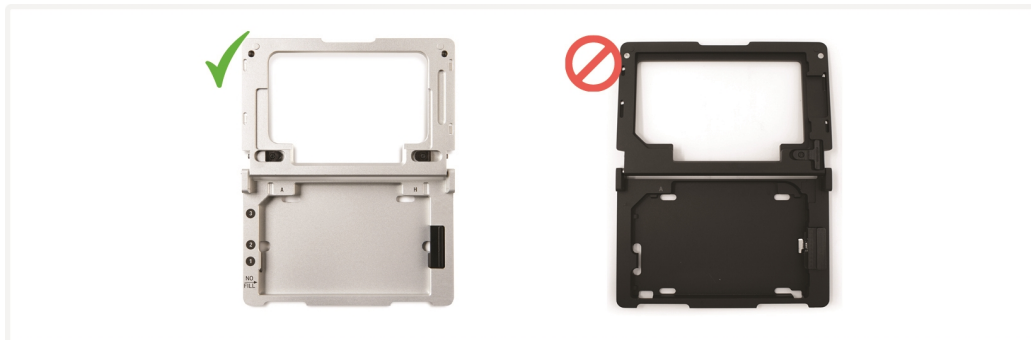

 See [Tips & Best Practices on page 26](#) for chip handling instructions.

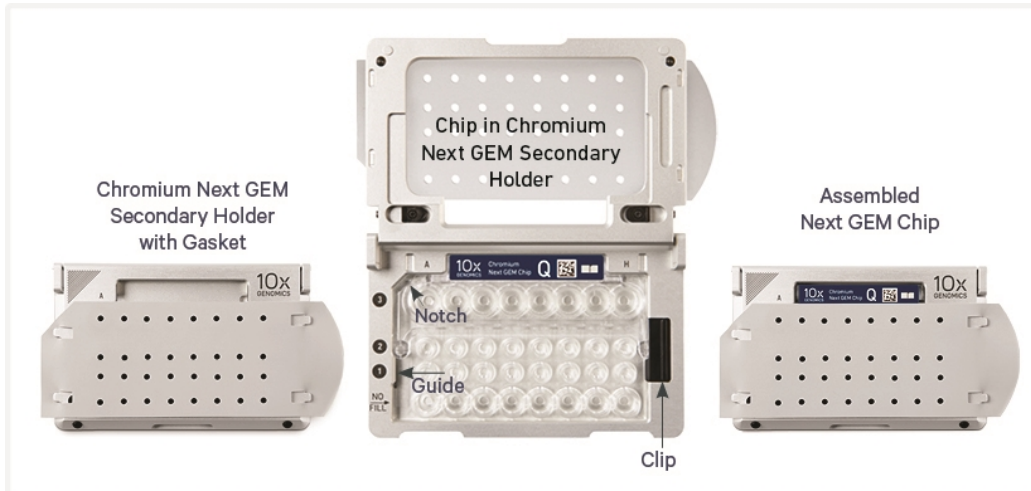

- a. Close the holder lid. Attach the gasket by holding the tongue (curved end, to the right) and hook the gasket on the left-hand tabs of the holder. Gently pull the gasket toward the right and hook it on the two right-hand tabs.
- b. DO NOT touch the smooth side of the gasket.
- c. Open the chip holder.
- d. Remove the chip from the sealed bag. Use the chip within  $\leq 24$  h.
- e. Align notch on the chip (upper left corner) and the open holder with the gasket attached.
- f. Slide the chip to the left until the chip is inserted under the guide on the holder. Depress the right hand side of the chip until the spring-loaded clip



engages.

- g.** Keep the assembled unit with the attached gasket open until ready for and while dispensing reagents into the wells.
- h.** DO NOT touch the smooth side of the gasket.
- i.** The assembled chip is ready for loading the indicated reagents. Refer to [2.3 Load Chromium Next GEM Chip Q on the next page](#) for reagent volumes and loading order.
- j.** After loading reagents, close the chip holder. DO NOT press down on the top of the gasket.



For GEM generation, load the indicated reagents only in the specified rows, starting from row labeled 1, followed by rows labeled 2 & 3. DO NOT load reagents in the bottom row labeled NO FILL.

## Chip Q 2.3 Load Chromium Next GEM Chip Q

Chip loading instructions are unique to Chip Q.

### TIPS

- After removing chip from the sealed bag, use in **≤24 h**.
- Open the lid (gasket attached) of the assembled chip and lay flat for loading.
- When loading the chip, raising and depressing the pipette plunger should each take **~5 sec**. When dispensing, raise the pipette tips at the same rate as the liquid is rising, keeping the tips slightly submerged.




### a. Add 50% glycerol solution to each unused well

- **70 µl** in each unused well in row labeled 1
- **50 µl** in each unused well in row labeled 2
- **45 µl** in each unused well in row labeled 3

**!** DO NOT add 50% glycerol solution to the bottom row of NO FILL wells. DO NOT use any substitute for 50% glycerol solution.


### Glycerol

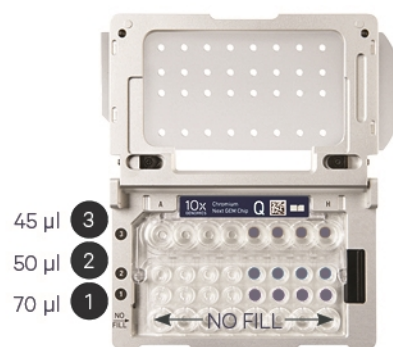

### b. Prepare Gel Beads

- Snap the tube strip holder with the Gel Bead strip into a 10x Vortex Adapter. Vortex **30 sec**.
- Centrifuge the Gel Bead strip for **~5 sec**. Confirm there are no bubbles at the bottom of the tubes & the liquid levels are even.
- Place the Gel Bead strip back in the holder. Secure the holder lid.

### Prep Gel Beads

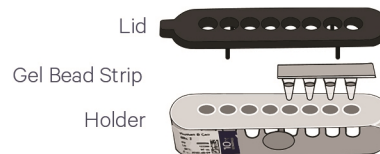

### c. Load Row 1


- With pipette set to 70 µl, gently **pipette mix** the GEM Master Mix + Sample 15x.
- Using the same pipette tips, dispense **70 µl** GEM Master Mix + Sample into the bottom center of wells in **row labeled 1** without introducing bubbles.

### Master Mix + Sample + Buffer



**d. Load Row Labeled 2**


- Puncture the foil seal of the Gel Bead tubes. Slowly aspirate **50 µl** Gel Beads.
- Dispense into the wells in **row labeled 2** without introducing bubbles.
- Wait **60 sec.**

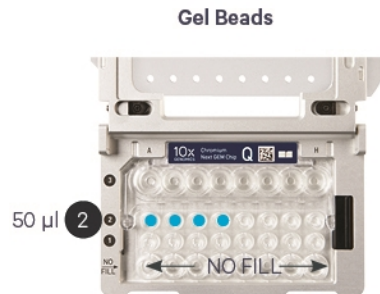

**e. Load Row Labeled 3**

- Dispense **45 µl** Partitioning Oil into the wells in **row labeled 3** from a reagent reservoir.


*Failure to add Partitioning Oil to the row labeled 3 will prevent GEM generation and can damage the Chromium X/iX.*

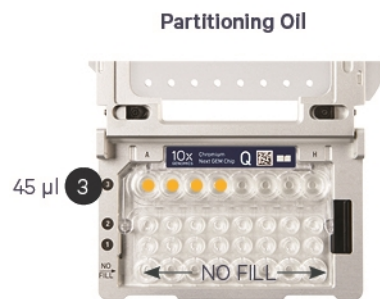

**f. Prepare for Run**

- Close the lid (gasket already attached). DO NOT touch the smooth side of the gasket. DO NOT press down on the top of the gasket.

*Run the chip in the Chromium X/iX immediately after loading the Partitioning Oil.*

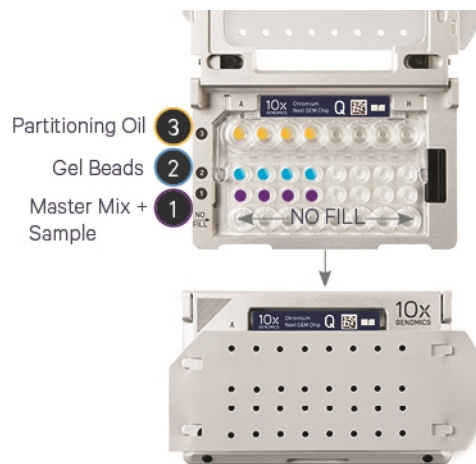



## 2.4 Run the Chromium X/iX

Consult the Chromium X Series (X/iX) User Guide (CG000396) for detailed instrument operation instructions and follow the instrument touchscreen prompts for execution. Run time for Chip Q is ~5.5 min.

- a. Press the eject button on the Chromium X to eject the tray.  
If the eject button is not touched within **1 min**, tray will close automatically. System requires a few seconds before the tray can be ejected again.
- b. Place the assembled chip with the gasket in the tray, ensuring that the chip stays horizontal. Press the button to retract the tray.
- c. Press the play button.
- d. At completion of the run (**~5.5 min**), Chromium X/iX will chime. **Immediately** proceed to the next step.



### Run Chromium X/iX



## 2.5 Transfer GEMs

- a. Place a tube strip on ice.
- b. Press the eject button of the Chromium X/iX and remove the chip.
- c. Discard the gasket. Open the chip holder. Fold the lid back until it clicks to expose the wells at 45 degrees.



- d. Check the volume in rows labeled 1-2. Abnormally high volume in any well indicates a clog.
- e. Slowly aspirate **100 µl** GEMs from the lowest points of the recovery wells in the top row labeled 3 without creating a seal between the tips and the bottom of the wells.

*In some cases, minor clogs may result in recovery of >90 µl but <100 µl of GEMs. Though the cell recovery efficiency might be slightly reduced, it's recommended to carry forward with library preparation and sequencing to recover information from the rest of the sample.*



- f. Withdraw pipette tips from the wells. GEMs should appear opaque and uniform across all channels.



- g. Over the course of **~20 sec**, dispense GEMs into the tube strip on ice with the pipette tips against the sidewalls of the tubes.

## 2.6 GEM Incubation

Use a thermal cycler that can accommodate at least 100  $\mu$ l volume. A volume of 100  $\mu$ l is the preferred setting on Bio-Rad C1000 Touch. In alternate thermal cyclers, use highest reaction volume setting.

- a. Incubate in a thermal cycler with the following protocol.

Lid Temperature	Reaction Volume	Run Time
80°C	100 $\mu$ l	~125 min
Step	Temperature	Time <i>hh:mm:ss</i>
1	25°C	00:60:00
2	60°C	00:45:00
3	80°C	00:20:00
Hold	4°C	Hold

- b. Store at 4°C for up to **a week**, or proceed to the next step.



**DO NOT store the GEMs at -20°C.**

- c. Sample placed on ice at step [2.2 Prepare GEM Master Mix + Sample Dilution on page 50](#) can either be discarded or stored at -20°C or -80°C for up to **1 month**. See Sample Storage below:

### Sample Storage

- Add 0.1 volume Enhancer to sample in Post-Hyb Resuspension Buffer. For example, add 50  $\mu$ l Enhancer to 500  $\mu$ l of sample in Post-Hyb Resuspension Buffer.
- Add 50% glycerol for a final concentration of 10%. For example, add 137.5  $\mu$ l 50% glycerol to 550  $\mu$ l sample in Post-Hyb Resuspension Buffer and Enhancer.
- Store at -20°C or -80°C for up to **1 month**. For best results, storage at -80°C is strongly recommended.

### Using Stored Samples

- When ready to use samples stored at -20°C or -80°C from this step, thaw at room temperature until no ice remains and then continue from step **2.1m** of [2.1 Post-Hybridization Wash on page 47](#) step to wash the sample once before proceeding to the step 2.2. Samples may undergo a color change during storage (e.g. black, light gray, or green), however this will not impact assay performance.

# Step 3:

## GEM Recovery and Pre-Amplification

3.0 Get Started	60
3.1 Post-GEM Incubation - Recovery	61
3.2 Pre-Amplification PCR	62
3.3 DNA Cleanup - SPRIselect	64



## 3.0 Get Started

Action	Item	10x PN	Preparation & Handling	Storage	
<b>Equilibrate to Room Temperature</b>					
<input type="checkbox"/>	<input type="radio"/>	<b>Reducing Agent B</b>	2000087	Thaw, vortex, verify no precipitate, centrifuge briefly.	-20°C
<input type="checkbox"/>	<input checked="" type="radio"/>	<b>Pre-Amp Primers B*</b>	2000529	Thaw, vortex, centrifuge briefly.	-20°C
<input type="checkbox"/>	<input checked="" type="radio"/>	<b>Feature Pre-Amp Primers*</b> Verify name & PN	2000515	Thaw, vortex, centrifuge briefly.	-20°C
<input type="checkbox"/>		<b>Beckman Coulter SPRIselect Reagent</b>	—	Manufacturer's recommendations.	Ambient
<b>Place on Ice</b>					
<input type="checkbox"/>	<input type="radio"/>	<b>Amp Mix</b>	2000047/ 2000103	Vortex and centrifuge briefly.	-20°C
<b>Obtain</b>					
<input type="checkbox"/>	<input type="radio"/>	<b>Recovery Agent</b>	220016	—	Ambient
<input type="checkbox"/>		<b>Qiagen Buffer EB</b>	—	Manufacturer's recommendations.	Ambient
<input type="checkbox"/>		<b>10% Tween 20</b>	—	Manufacturer's recommendations.	Ambient
<input type="checkbox"/>		<b>10x Magnetic Separator</b>	230003	—	Ambient
<input type="checkbox"/>		<b>Prepare 80% Ethanol</b> Prepare 2.5 ml for 4 GEM reactions.	—	Prepare fresh.	—

*\*Use Feature Pre-Amp Primers, PN-2000515, when samples were labeled with antibody-oligonucleotide conjugates prior to fixation. Use Pre-Amp Primers B, PN-2000529 when samples were not labeled.*





## 3.1 Post-GEM Incubation – Recovery



- a. Add **125 µl** Recovery Agent to each sample at **room temperature**. DO NOT pipette mix or vortex the biphasic mixture.
- b. Firmly secure the cap on the tube strip, ensuring that no liquid is trapped between the cap and the tube rim. Mix by inverting the capped tube strip 5x.



DO NOT invert without firmly securing the caps.

- c. Wait **2 min**.

*The resulting biphasic mixture contains Recovery Agent/Partitioning Oil (pink) and aqueous phase (translucent/opaque).*

*A smaller aqueous phase volume indicates a clog during GEM generation.*

- d. Centrifuge briefly.



- e. Slowly remove and discard **125 µl** Recovery Agent/Partitioning Oil (pink) from the bottom of the tube. DO NOT aspirate any aqueous sample.
- f. Proceed directly to Pre-Amplification PCR. No cleanup step is required.

### 3.2 Pre-Amplification PCR

- a. Prepare appropriate Pre-Amplification Mix on ice. Vortex and centrifuge briefly.

**For Samples Labeled with Antibody-oligonucleotide Conjugates:**

Pre-Amplification Mix		PN	1X (μl)	4X + 10% (μl)
<i>Add reagents in the order listed</i>				
○	<b>Amp Mix</b>	2000047	25.0	110.0
●	<b>Feature Pre-Amp Primers</b>	2000515	10.0	44.0
<b>Total</b>			<b>35.0</b>	<b>154.0</b>

**For Samples NOT Labeled with Antibody-oligonucleotide Conjugates:**

Pre-Amplification Mix		PN	1X (μl)	4X + 10% (μl)
<i>Add reagents in the order listed</i>				
○	<b>Amp Mix</b>	2000103	25.0	110.0
●	<b>Pre-Amp Primers B</b>	2000529	10.0	44.0
<b>Total</b>			<b>35.0</b>	<b>154.0</b>

- b. Add **35 μl** Pre-Amplification Mix to aqueous sample from step 3.1f.
- c. Cap firmly and invert 8x to mix. Centrifuge briefly.
- d. Incubate in a thermal cycler with the following protocol.

Lid Temperature	Reaction Volume	Run Time
105°C	100 μl	~30-45 min
Step	Temperature	Time hh:mm:ss
1	98°C	00:03:00
2	98°C	00:00:15
3	67°C	00:00:20
4	72°C	00:01:00
5	Go to Step 2, 7x (total 8 cycles)	
6	72°C	00:01:00
7	4°C	Hold




- e. Store at **4°C** for up to **72 h** or **-20°C** for **≤1 week**, or proceed to the next step.

### 3.3 DNA Cleanup – SPRIselect

- a. Prepare Elution Solution. Vortex and centrifuge briefly.

Elution Solution		PN	1000 $\mu$ l
<i>Add reagents in the order listed</i>			
	<b>Buffer EB</b>		980
	<b>10% Tween 20</b>	-	10
○	<b>Reducing Agent B</b>	2000087	10
	<b>Total</b>		<b>1000</b>

- b. Centrifuge the sample (PCR product) for 30 sec in a microcentrifuge and transfer **70  $\mu$ l** of the upper layer to a new tube.
- Presence of a cloudy precipitate at the interface between phases is normal. Avoid transferring the precipitate when transferring 70  $\mu$ l at this step.*
- c. Vortex to resuspend the SPRIselect reagent. Add **126  $\mu$ l** SPRIselect reagent (**1.8X**) to each sample and pipette mix 15x (pipette set to 180  $\mu$ l).
- d. Incubate **5 min** at **room temperature**.
- e. Place on the magnet•**High** until the solution clears.
- f. Remove the supernatant. DO NOT discard any beads.
- g. With the tube still in the magnet, add **200  $\mu$ l** 80% ethanol to the pellet. Wait **30 sec**.
- h. Remove the ethanol.
- i. **Repeat** steps g and h for a total of 2 washes.
- j. Centrifuge briefly and place on the magnet•**Low**.
- k. Remove any remaining ethanol. DO NOT let the sample dry to ensure maximum elution efficiency.
- l. Remove from the magnet. Add **101  $\mu$ l** Elution Solution. Wait **1 min** before resuspending. Pipette mix 15x.
- m. Incubate **2 min** at **room temperature**.
- n. Place the tube strip on the magnet•**High** until the solution clears.
- o. Transfer **100  $\mu$ l** sample to a new tube strip.
-  p. Store at **4°C** for **≤72 h** or at **-20°C** for **≤4 weeks**, or proceed to the next step.



# Step 4:

## Fixed RNA – Gene Expression Library Construction

4.0 Get Started	66
4.1 Sample Index PCR	67
4.2 Post Sample Index PCR Size Selection - SPRIselect	69
4.3 Post Library Construction QC	70



## 4.0 Get Started

Action	Item	10x PN	Preparation & Handling	Storage
<b>Equilibrate to Room Temperature</b>				
<input type="checkbox"/>	 <b>Dual Index Plate TS Set A</b> Verify name & PN. Use indicated plate only	3000511	Vortex and centrifuge briefly.	-20°C
<input type="checkbox"/>	<b>Beckman Coulter SPRIselect Reagent</b>	—	Manufacturer's recommendations.	—
<input type="checkbox"/>	<b>Agilent Bioanalyzer High Sensitivity Kit</b> If used for QC	—	Manufacturer's recommendations.	—
<input type="checkbox"/>	<b>Agilent TapeStation ScreenTape &amp; Reagents</b> If used for QC	—	Manufacturer's recommendations.	—
<b>Place on Ice</b>				
<input type="checkbox"/>	 <b>Amp Mix</b>	2000103/2000047	Vortex and centrifuge briefly.	-20°C
<input type="checkbox"/>	<b>KAPA Library Quantification Kit for Illumina Platforms</b>	—	Manufacturer's recommendations.	—
<b>Obtain</b>				
<input type="checkbox"/>	<b>Qiagen Buffer EB</b>	—	Manufacturer's recommendations.	Ambient
<input type="checkbox"/>	<b>10x Magnetic Separator</b>	230003	See Tips & Best Practices.	Ambient
<input type="checkbox"/>	<b>Prepare 80% Ethanol</b> Prepare 2.5 ml for 4 GEM reactions.	—	Prepare fresh.	Ambient

## 4.1 Sample Index PCR

- a. Choose the appropriate sample index sets to ensure that no sample indices overlap in a multiplexed sequencing run. Record the 10x sample index name (PN-3000511 Dual Index Plate TS Set A well ID) used.
- b. Prepare Sample Index PCR Mix on ice.

Sample Index PCR Mix					
<i>Add reagents in the order listed</i>		PN	1X (µl)	1X + 10% (µl)	4X + 10% (µl)
○	<b>Amp Mix</b>	2000103/2000047	50.0	55.0	220.0
	<b>Nuclease-free Water</b>	—	10.0	11.0	44.0
<b>Total</b>			<b>60.0</b>	<b>66.0</b>	<b>264.0</b>

- c. Transfer **ONLY 20 µl** sample from the step [DNA Cleanup – SPRiSelect on page 64](#) to a new tube strip.
- d. Add **60 µl** Sample Index PCR Mix to **20 µl** sample.
- e. Add **20 µl** of an individual Dual Index TS Set A to each sample. Pipette mix 5x (pipette set to 90 µl). Centrifuge briefly.
- f. Incubate in a thermal cycler with the following protocol.

Lid Temperature	Reaction Volume	Run Time
105°C	100 µl	~25-40 min
Step	Temperature	Time hh:mm:ss
1	98°C	00:00:45
2	98°C	00:00:20
3	54°C	00:00:30
4	72°C	00:00:20
5	Go to step 2, see table below for total # of cycles	
6	72°C	00:01:00
7	4°C	Hold

Targeted Cell Recovery	Total Cycles*			
	for Cell Lines	for PBMCs	for Cells from Fixed & Dissociated Tissues**	for Cells from FFPE Tissue Sections
500-2,000	12	16	15-16	17
2,000-4,000	11	15	14-15	16
4,000-7,000	10	14	13-14	15
7,000-10,000	9	13	12-13	14

\*Optimization of cycle number may be needed based on the total RNA content of the sample. Target concentration is between 100 nM and 200 nM. If optimization is needed, additional Amp Mix can be obtained using the Library Amplification Kit (PN-1000249). For dissociated tumor cells, cycle numbers for cell lines can be used as a starting point. For dissociated primary cells, cycle numbers for PBMCs can be used as a starting point.


\*\*For cells derived from the fixed and dissociated tissue samples, the cycle number will depend on the RNA expression level of the tissue and on overall quality of the tissue prior to fixation. Additional optimization may be required.



**g.** Store at 4°C for ≤72 h, or proceed to the next step.




## 4.2 Post Sample Index PCR Size Selection – SPRIselect

- a. Vortex to resuspend the SPRIselect reagent. Add **100 µl** SPRIselect Reagent (**1.0X**) to each sample. Pipette mix 15x (pipette set to 180 µl).
- b. Incubate **5 min** at **room temperature**.
- c. Place on the magnet•**High** until the solution clears.
- d. Remove the supernatant. DO NOT discard any beads.
- e. With the tube still in the magnet, add **200 µl** 80% ethanol to the pellet. Wait **30 sec**.
- f. Remove the ethanol.
- g. **Repeat** steps e and f for a total of 2 washes.
- h. Centrifuge briefly and place on the magnet•**Low**.
- i. Remove any remaining ethanol. DO NOT let the sample dry to ensure maximum elution efficiency.
- j. Remove from the magnet. Add **41 µl** Buffer EB. Pipette mix 15x.
- k. Incubate **2 min** at **room temperature**.
- l. Place on the magnet•**Low** until the solution clears.
- m. Transfer **40 µl** to a new tube strip.
-  n. Store at **4°C** for up to **72 h** or at **-20°C** for **long-term** storage.

## 4.3 Post Library Construction QC

Run **1  $\mu$ l** sample at **1:80** dilution on an Agilent Bioanalyzer High Sensitivity chip.

### Representative Trace

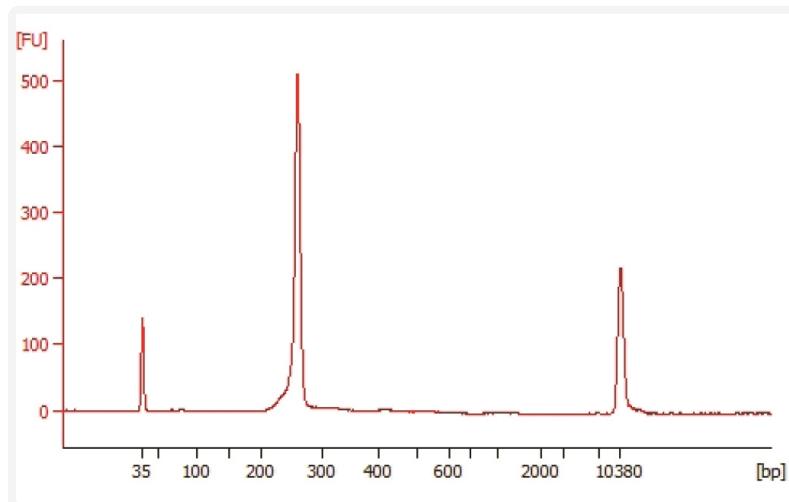

Determine the average fragment size from the Bioanalyzer trace.

### Alternate QC Method

Agilent TapeStation

LabChip

See Appendix for:

- [Post Library Construction Quantification on page 88](#)
- [Agilent TapeStation Traces on page 89](#)
- [LabChip Traces on page 89](#)

# Step 5:



## Fixed RNA – Cell Surface Protein Library Construction

5.0 Get Started	72
5.1 Sample Index PCR	73
5.2 Post Sample Index PCR Size Selection - SPRIselect	75
5.3 Post Library Construction QC	76



## 5.0 Get Started

### Overview

Action	Item	10x PN	Preparation & Handling	Storage
<b>Equilibrate to Room Temperature</b>				
<input type="checkbox"/>	 <b>Dual Index Plate TT Set A</b> <i>Verify name &amp; PN. Use indicated plate only</i>	3000431	Vortex and centrifuge briefly.	-20°C
<input type="checkbox"/>	<b>Beckman Coulter SPRiSelect Reagent</b>	—	Manufacturer's recommendations.	—
<input type="checkbox"/>	<b>Agilent Bioanalyzer High Sensitivity Kit</b> <i>If used for QC</i>	—	Manufacturer's recommendations.	—
<input type="checkbox"/>	<b>Agilent TapeStation ScreenTape &amp; Reagents</b> <i>If used for QC</i>	—	Manufacturer's recommendations.	—
<b>Place on Ice</b>				
<input type="checkbox"/>	 <b>Amp Mix</b>	2000047/ 2000103	Vortex and centrifuge briefly.	-20°C
<input type="checkbox"/>	<b>KAPA Library Quantification Kit for Illumina Platforms</b>	—	Manufacturer's recommendations.	—
<b>Obtain</b>				
<input type="checkbox"/>	<b>Qiagen Buffer EB</b>	—	Manufacturer's recommendations.	Ambient
<input type="checkbox"/>	<b>10x Magnetic Separator</b>	230003	See Tips & Best Practices.	Ambient
<input type="checkbox"/>	<b>Prepare 80% Ethanol</b> <i>Prepare 2.5 ml for 4 GEM reactions.</i>	—	Prepare fresh.	Ambient

## 5.1 Sample Index PCR

- a. Choose the appropriate sample index sets to ensure that no sample indices overlap in a multiplexed sequencing run. Record the 10x sample index name (PN-3000431 Dual Index Plate TT Set A well ID) used.
- b. Prepare Sample Index PCR Mix.

Sample Index PCR Mix		PN	1X ( $\mu$ l)	1X + 10% ( $\mu$ l)	4X + 10% ( $\mu$ l)
<i>Add reagents in the order listed</i>					
○	<b>Amp Mix</b>	2000047/ 2000103	50.0	55.0	220.0
	<b>Nuclease-free Water</b>	—	10.0	11.0	44.0
<b>Total</b>			<b>60.0</b>	<b>66.0</b>	<b>264.0</b>

- c. Transfer **ONLY 20  $\mu$ l** sample (pre-amplified DNA) from the [DNA Cleanup – SPRiSelect on page 64](#) to a new tube strip.
- d. Add **60  $\mu$ l** Sample Index PCR Mix to **20  $\mu$ l** sample.
- e. Add **20  $\mu$ l** of an individual Dual Index TT Set A to each sample and record the well ID used. Pipette mix 5x (pipette set to 90  $\mu$ l). Centrifuge briefly.
- f. Incubate in a thermal cycler with the following protocol.

Lid Temperature	Reaction Volume	Run Time
105°C	100 $\mu$ l	~25-40 min
Step	Temperature	Time hh:mm:ss
1	98°C	00:00:45
2	98°C	00:00:20
3	54°C	00:00:30
4	72°C	00:00:20
5	Go to step 2, see table below for total # of cycles	
6	72°C	00:01:00
7	4°C	Hold


*\*Optimization of cycle number may be needed based on target protein expression levels and number of antibodies used for labeling.*

Targeted Cell Recovery	Total Cycles
500-2,000	14-15
2,000-4,000	13-14
4,000-10,000	12-13



**g.** Store at **4°C** for **≤72 h**, or proceed to the next step.


## 5.2 Post Sample Index PCR Size Selection – SPRIselect

- a. Vortex to resuspend the SPRIselect reagent. Add **100 µl** SPRIselect Reagent (**1.0X**) to each sample. Pipette mix 15x (pipette set to 180 µl).
- b. Incubate **5 min** at **room temperature**.
- c. Place on the magnet•**High** until the solution clears.
- d. Remove the supernatant. DO NOT discard any beads.
- e. With the tube still in the magnet, add **200 µl** 80% ethanol to the pellet. Wait **30 sec**.
- f. Remove the ethanol.
- g. **Repeat** steps e and f for a total of 2 washes.
- h. Centrifuge briefly and place on the magnet•**Low**.
- i. Remove any remaining ethanol. DO NOT let the sample dry to ensure maximum elution efficiency.
- j. Remove from the magnet. Add **41 µl** Buffer EB. Pipette mix 15x.
- k. Incubate **2 min** at **room temperature**.
- l. Place on the magnet•**Low** until the solution clears.
- m. Transfer **40 µl** to a new tube strip.
-  n. Store at **4°C** for up to **72 h** or at **-20°C** for **long-term** storage.

## 5.3 Post Library Construction QC

Run **1  $\mu$ l** sample at **1:80** dilution on an Agilent Bioanalyzer High Sensitivity chip.

### Representative Trace

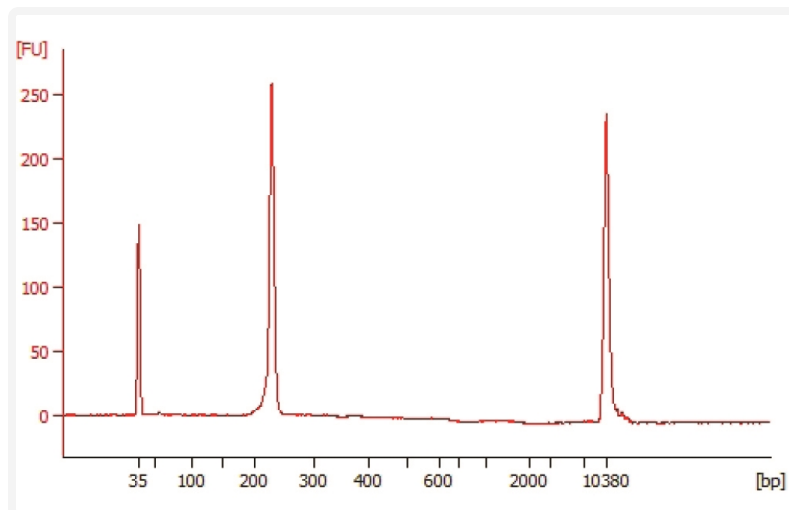

Determine the average fragment size from the Bioanalyzer trace.

### Alternate QC Method

Agilent TapeStation

LabChip

See Appendix for:

- [Post Library Construction Quantification on page 88](#)
- [Agilent TapeStation Traces on page 89](#)
- [LabChip Traces on page 89](#)



# Step 6:

## Sequencing


Sequencing Libraries	78
Illumina Sequencer Compatibility	78
Sample Indices	79
Fixed RNA - Gene Expression Library Sequencing Parameters	79
Fixed RNA - Cell Surface Protein Library Sequencing Parameters	79
Library Loading	80
Library Pooling	80
Data Analysis and Visualization	80




## Sequencing Libraries

Chromium Fixed RNA Profiling – Gene Expression and Cell Surface Protein libraries comprise standard Illumina paired-end constructs which begin with P5 and end with P7. These libraries include 16 bp 10x GEM Barcodes (10x Barcode) encoded at the start of TruSeq Read 1 (Read 1T). Sample index sequences are incorporated as the i5 and i7 index reads. TruSeq Read 1 (Read 1T) and Small RNA Read 2 (Read 2S) are used in paired-end sequencing of Fixed RNA – Gene Expression libraries. TruSeq Read 1 (Read 1T) and TruSeq Read 2 (Read 2T) are used for paired-end sequencing of Fixed RNA – Cell Surface Protein library. Sequencing these libraries produces a standard Illumina BCL data output folder.

### Chromium Fixed RNA Profiling – Gene Expression Library

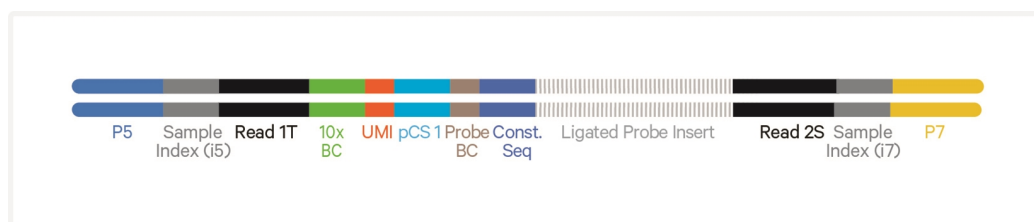

### Chromium Fixed RNA Profiling – Cell Surface Protein Library

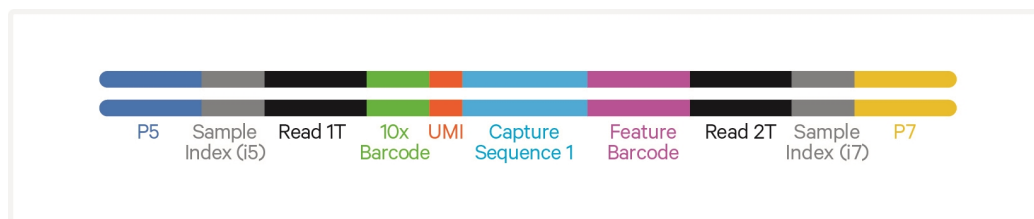

## Illumina Sequencer Compatibility

The compatibility of the listed sequencers has been verified by 10x Genomics. Some variation in assay performance is expected based on sequencer choice. For more information about performance variation, visit the 10x Genomics Support website.

- iSeq
- MiSeq
- NextSeq 500/550
- NextSeq 1000/2000
- NovaSeq

## Sample Indices

Each sample index in the Dual Index Kit TT Set A (PN-1000215) or Dual Index Kit TS Set A (PN-1000251) is a mix of one unique i7 and one unique i5 sample index. If multiple samples are pooled in a flow cell lane, the sample index name (i.e. the Dual Index TS Set A plate well ID, SI-TS-) is needed in the sample sheet used for generating FASTQs with “cellranger mkfastq”. Samples utilizing the same sample index should not be pooled together, or run on the same flow cell lane, as this would not enable correct sample demultiplexing.

## Fixed RNA – Gene Expression Library Sequencing Parameters

Parameter	Description
Sequencing Depth	Minimum 10,000 read pairs per cell
Sequencing Type	Paired-end, dual indexing
Sequencing Read	Recommended Number of Cycles
Read 1	28 cycles
i7 Index	10 cycles
i5 Index	10 cycles
Read 2	90 cycles <i>(Minimum required Read 2 length is 50 bp)</i>

## Fixed RNA – Cell Surface Protein Library Sequencing Parameters

Pooling Fixed RNA – Gene Expression & Cell Surface Protein libraries is recommended for sequencing to maintain nucleotide diversity

Parameter	Description
Sequencing Depth	Minimum 5,000 read pairs per cell
Sequencing Type	Paired-end, dual indexing
Sequencing Read	Recommended Number of Cycles
Read 1	28 cycles
i7 Index	10 cycles
i5 Index	10 cycles
Read 2	90 cycles <i>(Minimum required Read 2 length is 25 bp)</i>

## Library Loading

Once quantified and normalized, the libraries should be denatured and diluted as recommended for Illumina sequencing platforms. Refer to Illumina documentation for denaturing and diluting libraries. Refer to the 10x Genomics Support website for more information.

**Fixed RNA – Gene Expression libraries alone or in combination with Cell Surface Protein libraries**

### Library Loading

Instrument	Loading Concentration (pM)	PhiX (%)
MiSeq	12	1
NextSeq 500/550	1.6	1
NextSeq 1000/2000	650	1
NovaSeq	150*/300	1

\* Use 150 pM loading concentration for Illumina XP workflow.

## Library Pooling

Fixed RNA – Gene Expression and Cell Surface Protein libraries may be pooled for sequencing, taking into account the differences in cell number and per-cell read depth requirements between each library. Samples utilizing the same sample index should not be pooled together or run on the same flow cell lane, as this would not enable correct sample demultiplexing.

### Library Pooling Example

Libraries	Sequencing Depth (read pairs per cell)	Library Pooling Ratio
<b>Example</b>		
Fixed RNA – Gene Expression library	10,000	2
Fixed RNA – Cell Surface Protein library	5,000	1

## Data Analysis and Visualization

Sequencing data may be analyzed using Cell Ranger or 10x Genomics Cloud Analysis and visualized using Loupe Browser. Key features for these tools are listed below. For detailed product-specific information, visit the 10x Genomics Support website.

## Cell Ranger

Cell Ranger is a set of analysis pipelines that processes Chromium Single Gene Expression data to align reads, generate Feature Barcode matrices and perform clustering and gene expression analysis.

- Input: Base call (BCL) and FASTQ
- Output: BAM, MEX, CSV, HDF5, Web Summary, .cloupe/.loupe
- Operating System: Linux

## Cloud Analysis

*Cloud Analysis is currently only available for US customers.*

Cloud Analysis allows users to run Cell Ranger analysis pipelines from a web browser while computation is handled in the cloud.

- Key features: scalable, highly secure, simple to set up and run
- Input: FASTQ
- Output: BAM, MEX, CSV, HDF5, Web Summary, .cloupe/.loupe

## Loupe Browser

Loupe Browser is an interactive data visualization tool that requires no prior programming knowledge.

- Input: .cloupe
- Output: Data visualization, including t-SNE and UMAP projections, custom clusters, differentially expressed genes
- Operating System: MacOS, Windows

# Troubleshooting


GEM Generation & Barcoding	83
Post Library Construction QC	84
Chromium X Series Errors	85




## GEM Generation & Barcoding

Step	Normal	Abnormal
------	--------	----------

After loading the Chip

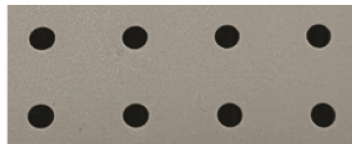

The gasket holes are aligned with the wells.

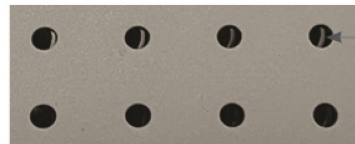

Misaligned gasket holes & chip wells

Gasket holes are misaligned with the sample and gel bead wells. Open and close the chip holder slowly once.

After Chip Q is removed from X/iX and the wells are exposed




All recovery wells are similar in volume and opacity.




Recovery well G indicates a reagent clog. Recovery well C and E indicates a wetting failure. Recovery wells B, D, and F are normal. Wells A and H contain 50% Glycerol Solution.

Transfer GEMs

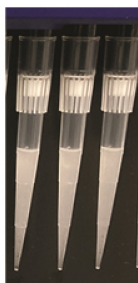

All liquid levels are similar in volume and opacity without air trapped in the pipette tips.

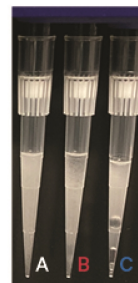

- Adequate emulsion volume (no clog or wetting failure)
- Wetting failure
- Low emulsion volume (clog)


Pipette tip A shows normal GEM generation, pipette tip B indicates a wetting failure, and pipette tip C shows a clog and wetting failure.

Consult Best Practices to Minimize Chromium Next GEM Chip Clogs and Wetting Failures (Technical Note CG000479) for more information. If a channel clogs or wetting failure occurs during GEM generation, it is recommended that the sample be remade. If any of the listed issues occur, take a picture and send it to [support@10xgenomics.com](mailto:support@10xgenomics.com) for further assistance.

## Post Library Construction QC

### Step


Fixed RNA – Gene Expression Library  
 Correct Sample Index PCR cycling

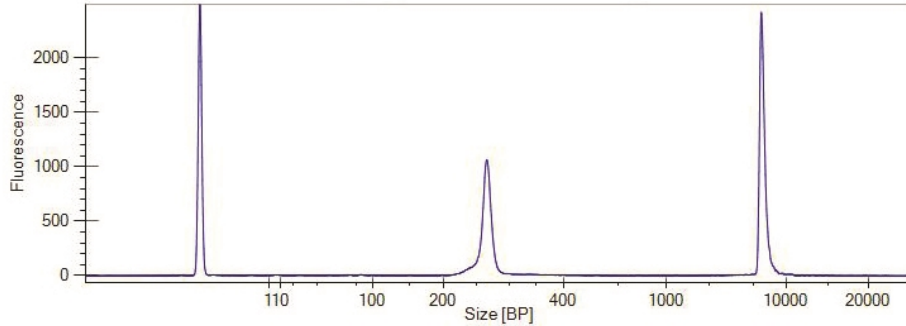

Fixed RNA – Gene Expression libraries should fall within the range of 100-200 nM. If the concentration falls below or above this range, sequencing data quality may be impacted, sample index PCR should be repeated with greater or fewer number of PCR cycles, respectively. If concentration is >200 nM and does not contain a high proportion of higher molecular weight peaks, proceed with sequencing.


### Over cycling

Additional higher molecular weight peaks present in the library trace indicate over cycling.

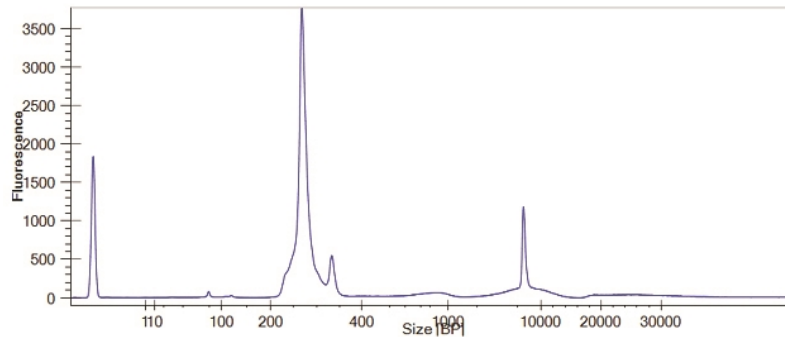

### Under cycling

Higher proportion of low molecular weight peaks present in the library trace indicate under cycling.

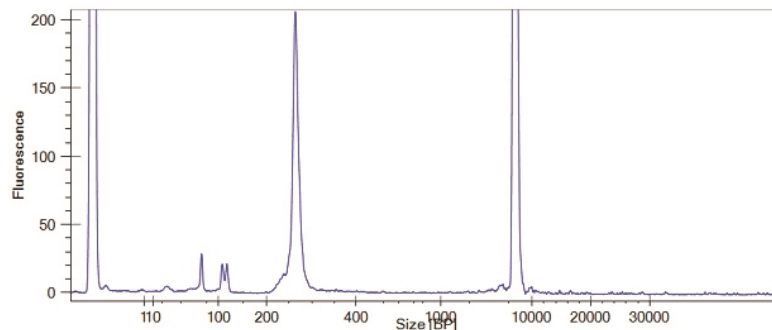



## Chromium X Series Errors

The Chromium X touchscreen will guide the user through recoverable errors. If the error continues, or if the instrument has seen critical or intermediate errors, email [support@10xgenomics.com](mailto:support@10xgenomics.com) with the displayed error code. Support will request a troubleshooting package. Upload pertinent logs to 10x Genomics by navigating to the Logs menu option on screen.

### There are two types of errors:

**Critical Errors** – When the instrument has seen a critical error, the run will immediately abort. Do not proceed with any further runs. Contact [support@10xgenomics.com](mailto:support@10xgenomics.com) with the error code.

- a. System Error
- b. Pressure Error
- c. Chip Error
- d. Run Error
- e. Temperature Error
- f. Software Error

**User Recoverable Errors** – Follow error handling instructions through the touchscreen and continue the run.

- a. Gasket Error
- b. Tray Error
- c. Chip Error
- d. Unsupported Chip Error
- e. Network Error
- f. Update Error



*Consult the Chromium X Series (X/iX) User Guide (CG000396) for additional information and follow the Chromium X touchscreen prompts for execution. The Chromium X touchscreen will guide the user through recoverable errors.*

# Appendix

Chromium Fixed RNA Profiling: Chip Loading Overview	87
Post Library Construction Quantification	88
Agilent TapeStation Traces	89
LabChip Traces	89
Oligonucleotide Sequences	91



## Chip Q Chromium Fixed RNA Profiling: Chip Loading Overview

This section provides a quick overview to the Chip Q loading and does not include detailed instructions. Refer to [Load Chromium Next GEM Chip Q on page 54](#) for details.

### Steps

**a. Add 50% glycerol solution to each unused well**

- Load 70  $\mu$ l to row labeled 1
- Load 50  $\mu$ l to row labeled 2
- Load 45  $\mu$ l to row labeled 3

**b. Prepare Gel Beads**

- Vortex for 30 sec
- Centrifuge for 5 sec

**c. Load Row Labeled 1**

- Mix GEM Master Mix + Sample
- Load 70  $\mu$ l to row labeled 1


**d. Load Row Labeled 2**

- Aspirate Gel Beads
- Load 50  $\mu$ l to row labeled 2
- Wait 60 sec

**e. Load Row Labeled 3**

- Load 45  $\mu$ l Partitioning Oil to row labeled 3

**f. Close the lid and prepare for run.**

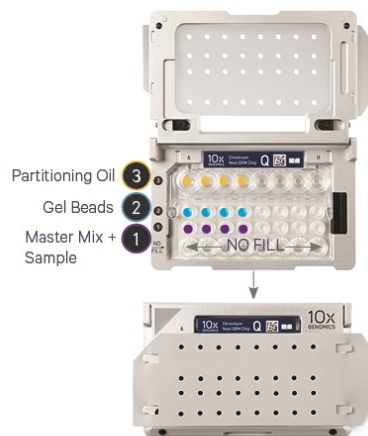

## Post Library Construction Quantification

- a. Thaw KAPA Library Quantification Kit for Illumina Platforms.
- b. Dilute **2 µl** sample with deionized water to appropriate dilutions that fall within the linear detection range of the KAPA Library Quantification Kit for Illumina Platforms. (For more accurate quantification, make the dilution(s) in duplicate).
- c. Make enough Quantification Master Mix for the DNA dilutions per sample and the DNA Standards (plus 10% excess) using the guidance for 1 reaction volume below.

Quantification Master Mix	1X (µl)
SYBR Fast Master Mix + Primer	12
Water	4
<b>Total</b>	<b>16</b>



- d. Dispense **16 µl** Quantification Master Mix for sample dilutions and DNA Standards into a 96 well PCR plate.
- e. Add **4 µl** sample dilutions and **4 µl** DNA Standards to appropriate wells. Centrifuge briefly.
- f. Incubate in a thermal cycler with the following protocol.

Step	Temperature	Run Time
1	95°C	00:03:00
2	95°C	00:00:05
3	67°C Read Signal	00:00:30
4	Go to Step 2, 29X (Total 30 cycles)	

- g. Follow the manufacturer's recommendations for qPCR-based quantification. For library quantification for sequencer clustering, determine the concentration based on insert size derived from the Bioanalyzer/TapeStation trace.

## Agilent TapeStation Traces


Agilent TapeStation High Sensitivity D5000 ScreenTape was used.

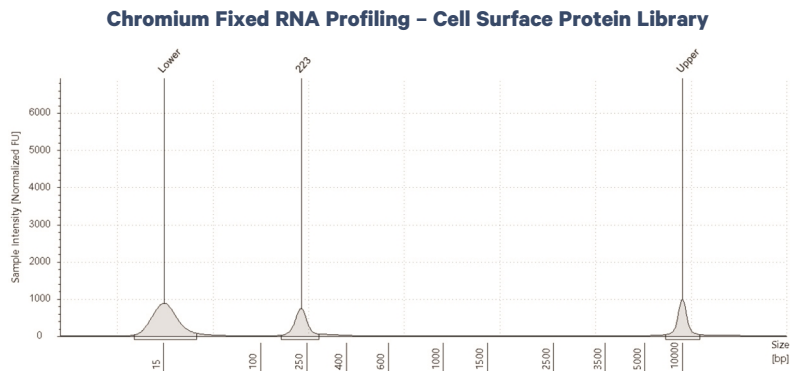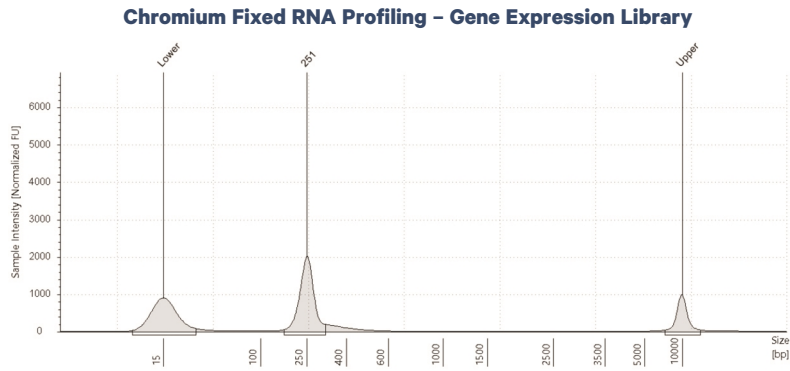


All traces are representative. Samples were run at 1:80 dilution.

## LabChip Traces

DNA High Sensitivity Reagent Kit was used.

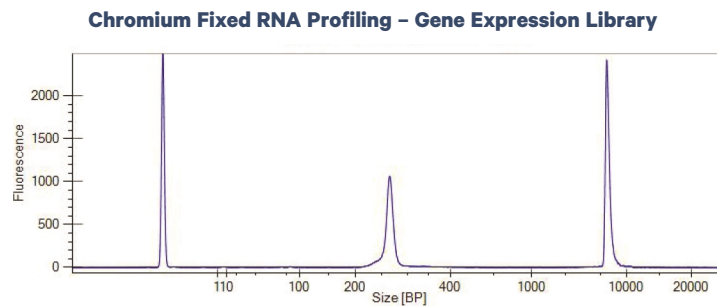

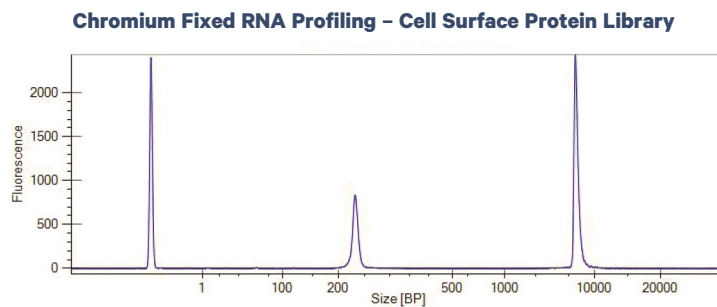

*All traces are representative. Samples were run at 1:80 dilution.*

# Oligonucleotide Sequences


## Gel Bead Primer

5'-CTACACGACGCTCTCCGATCT-N16-N12-TTGCTAGGACCG-3'



## Chromium Fixed RNA Profiling – Gene Expression Library

5'-AATGATACGGCGACCCAGAGATCTACAC-N10-ACACTCTTCCCTACACGACGCTCTCCGATCT-N16-N12-TTGCTAGGACCG-BC8-NN-TACGTGCTAACCGCGT-Ligated\_Probe\_Insert-TGGAATTCGCGGTGCCAAGGAAGCTCCAGTCAC-N10-ATCTCGTATGCCGCTCTCTCGCTTG-3'  
 3'-TACATATGCCCGCTGGTGGCTCTAGATGTG-N10-TGTGAGAAAGGGATGTGCTGCGAGAAGGCTAGA-N16-N12-AACGATCCTGCGC-BC8-NN-ATGCACGATTGGCGCA- Ligated\_Probe\_Insert-ACCCTAAGAGCCACGGTTCTCTGAGGTCAGTG-N10-TAGAGCATACGGCAGAAGACGAAC-5'

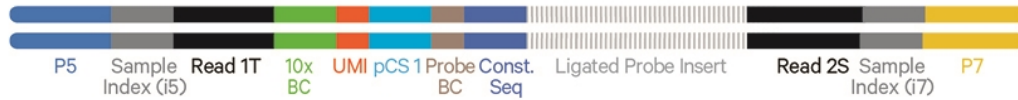

## Cell Surface Protein Feature Barcode

3'-AACGATCCTGGCCGGAATTCG-N9-N15-N10-TCTAGCCTTCTCGTGTGCAGACTTGAGGTCAGTG-5'



## Chromium Fixed RNA Profiling – Cell Surface Protein Library

5'-AATGATACGGCGACCCAGAGATCTACAC-N10-ACACTCTTCCCTACACGACGCTCTCCGATCT-N16-N12-TTGCTAGGACCGCCTTAAAGC-N9-N15-N10-AGATCGGAAGAGCACAGCTCTGAATCCAGTCAC-N10-ATCTCGTATGCCGCTCTCTCGCTTG-3'  
 3'-TACATATGCCCGCTGGTGGCTCTAGATGTG-N10-TGTGAGAAAGGGATGTGCTGCGAGAAGGCTAGA-N16-N12-AACGATCCTGCGC-BC8-NN-ATGCACGATTGGCGCA- Ligated\_Probe\_Insert-ACCCTAAGAGCCACGGTTCTCTGAGGTCAGTG-N10-TAGAGCATACGGCAGAAGACGAAC-5'

