

Xenium In Situ Gene Expression - Post-Xenium Analyzer H&E Staining

Introduction

Xenium In Situ Gene Expression measures mRNA in formalin fixed & paraffin embedded (FFPE) and fresh frozen (FF) tissue sections that are placed onto Xenium Slides. Following a Xenium Analyzer instrument run, tissue sections on Xenium slides may be optionally processed for Hematoxylin & Eosin (H&E) staining. This document provides guidance for disassembly of the Xenium Cassette, removal of Autofluorescence Solution, and H&E staining steps for both FFPE and FF tissue sections. After Xenium Cassette Removal and Quencher Removal, follow the guidelines provided for Post-Xenium Analyzer H&E Staining or proceed to any alternate H&E staining protocol, if desired.

Additional Guidance

Proceed with the Post-Xenium Analyzer H&E Staining protocol following the Xenium Analyzer User Guide (CG000584) to obtain histological data that can be combined with gene expression data from the same tissue section. H&E Staining post-Xenium Analyzer instrument run is optional.

Contents

Workflow Overview

Specific Reagents & Consumables

Tips & Best Practices

Icons	6
Pipette Calibration	6
Xenium Slide Handling	6
Slide Storage	7
Xenium Cassette Removal	7
Slide Incubation Guidance	8
Quencher Removal	9
H&E Staining	9
Protocol Steps & Timing	10

1. H&E Staining


1.0 Overview	11
Protocol Overview	12
Get Started - Quencher Removal	14
1.1 Quencher Removal Preparation	15
1.2 Cassette Removal	17
1.3 Quencher Removal	18
Get Started - H&E Staining	19
1.4 H&E Staining Preparation	20
1.5 H&E Staining	21
1.6 Coverslipping	23

2. Troubleshooting

Dry Tissue	24
------------	----

Document Revision Summary

Workflow Overview

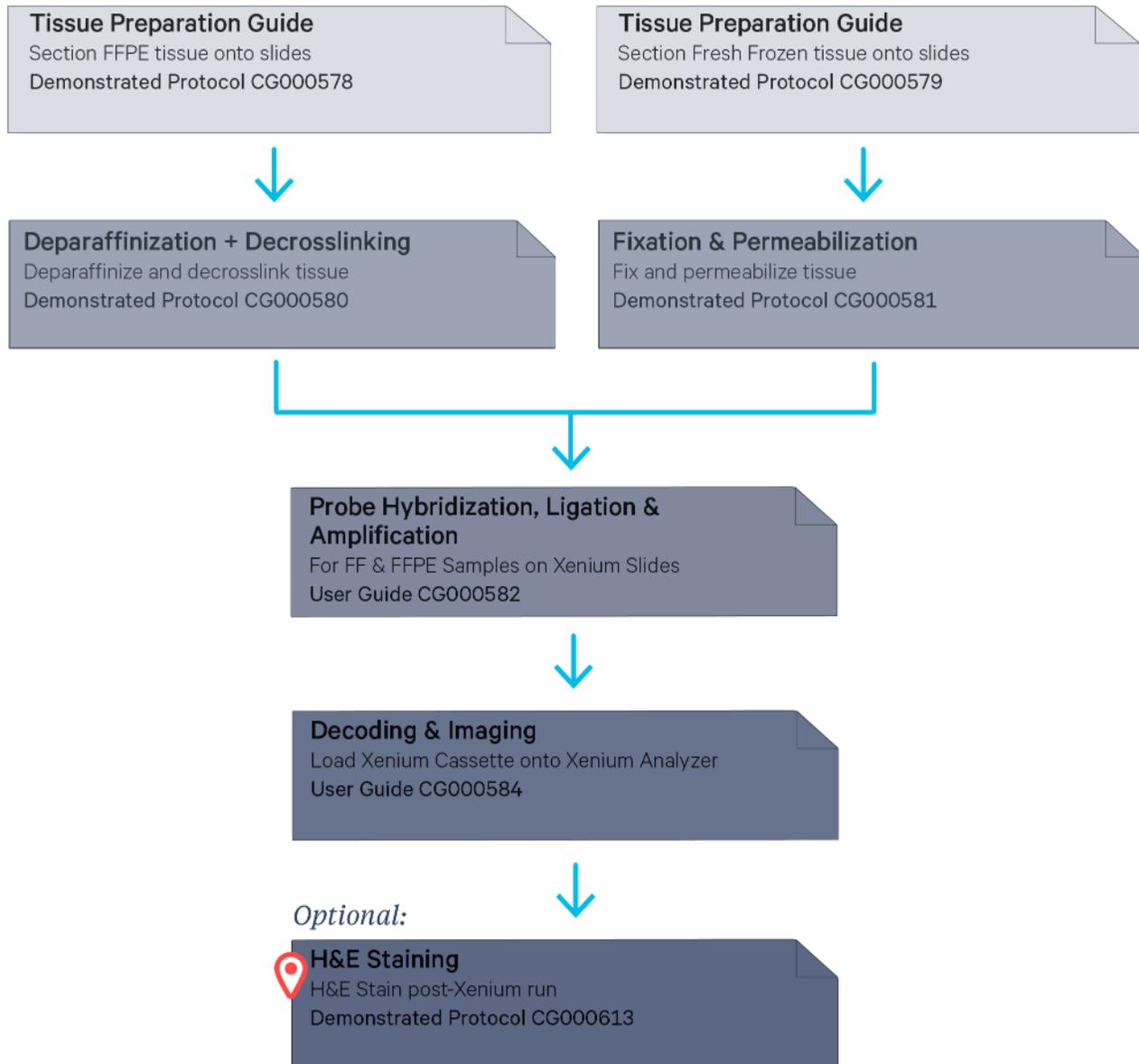

Visit the 10x Genomics Support website for the most current documentation.

Specific Reagents & Consumables

H&E Staining

The listed items have been tested by 10x Genomics and perform optimally with the assay. **Substituting materials may adversely affect system performance.** For items with multiple options listed, choose option based on availability and preference. Refer to the manufacturer's website for regional part numbers.

For Quencher Removal				
Item		Description	Vendor	Part Number
<input type="checkbox"/>	Sodium hydrosulfite	Sodium hydrosulfite, technical grade (or equivalent)	Sigma Aldrich	157953-5G
<input type="checkbox"/>	Forceps	Tweezers, 4" Wafer Handling	Excelta Corp	491P-SA-PI
<input type="checkbox"/>	PBS (optional)	PBS - Phosphate Buffered Saline (10X) pH 7.4, RNase-free	Thermo Fisher Scientific	AM9624
<input type="checkbox"/>	Slide Mailers	Simport Scientific LockMailer Tamper Evident Slide Mailer	Fisher Scientific	22-038-399
For H&E Staining				
Item		Description	Vendor	Part Number
<input type="checkbox"/>	Hematoxylin	Hematoxylin Solution, Mayer's	Sigma Aldrich	MHS16
<input type="checkbox"/>	Eosin	Eosin Y Solution, Alcoholic	Leica	3801615
<input type="checkbox"/>	Bluing Reagent	Bluing Solution	Dako	CS702
<input type="checkbox"/>	Mounting Media	Surgipath SUB-X Mounting Media (discontinued)	Leica	3801741
		Cytoseal Mountant; 60, 280 and XYL	Fisher Scientific	22-050-262
<input type="checkbox"/>	Ethanol	Ethyl Alcohol, 200 Proof, anhydrous	Millipore Sigma	E7023
		Ethanol absolute ≥99.5%, TechniSolv, pure (Europe)	VWR	83813.360DP
<input type="checkbox"/>	Xylene	Xylene, Reagent Grade	Millipore Sigma	214736
		Xylene, Histological Grade	Millipore Sigma	534056
<input type="checkbox"/>	Forceps	Tweezers, 4" Water Handling	Excelta Corp	491P-SA-PI
<input type="checkbox"/>	Filter Paper	Fisherbrand Qualitative Grade Plain Filter Paper Circles	Fisher Scientific	09-795-H
<input type="checkbox"/>	Coverslips	Fisherbrand Cover Glasses: Rectangles	Fisher Scientific	12-544-EP
		Cover Glasses, Rectangles	VWR	16004-322

For H&E Staining	
Additional Materials	
<input type="checkbox"/>	Vortex
<input type="checkbox"/>	Staining jar/dishes
<input type="checkbox"/>	Wide-bore pipette tips
<input type="checkbox"/>	Ultrapure/Milli-Q Water <i>from Milli-Q Integral Ultrapure Water System or equivalent</i>

This list may not include some standard laboratory equipment.

Tips & Best Practices

Icons



Tips & Best Practices section includes additional guidance



Signifies critical step requiring accurate execution




Troubleshooting section includes additional guidance

Pipette Calibration

- Follow manufacturer’s calibration and maintenance schedules.

Xenium Slide Handling

- Always wear gloves when handling slides.
- The bottom of the slide is indicated by the etched label, which should be readable when in the proper position.
- Hold the slide on the label. DO NOT touch the tissue sections or near fiducials.
- Minimize exposure of the slides to sources of particles and fibers.

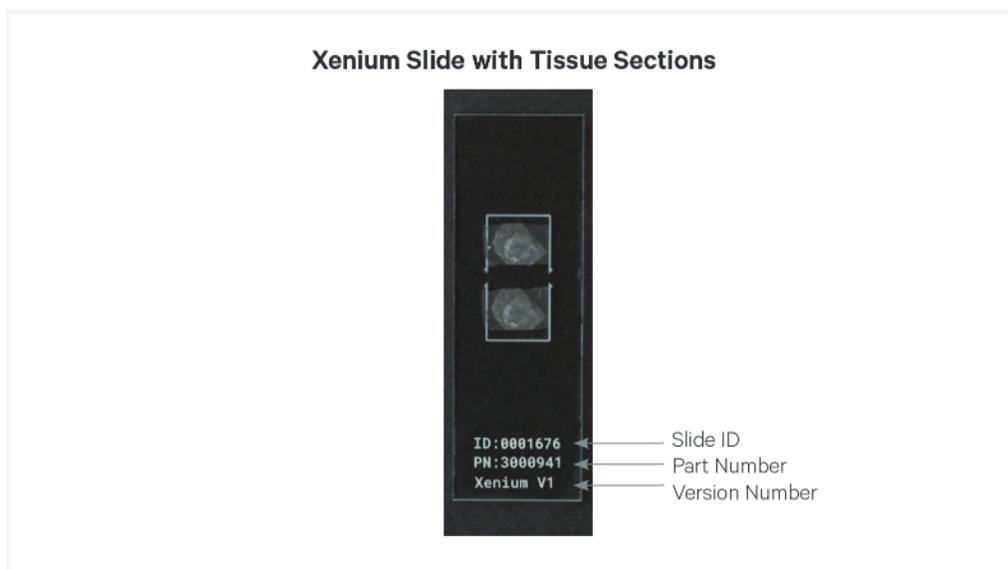

Slide Storage

- Store analyzed Xenium Cassettes in PBS-T at 4°C following a Xenium Analyzer instrument run. Consult the Xenium Analyzer User Guide (CG000584) for more information about post-Xenium Analyzer storage conditions.
- Proceed to H&E Staining within 3 days following an instrument run.
- For long-term storage after H&E staining and imaging, maintain at room temperature in the dark.

Xenium Cassette Removal

Pull inner clips from inner tabs to detach top and bottom halves of cassette



Open cassette by continuing to lift inner clips upward



Hold slide by the label and lift slide out from bottom half



Slides in images are representative.

Slide Incubation Guidance

Incubation at room temperature

- Place Xenium slide(s) with label toward the top of the slide mailer or coplin jar for incubations at room temperature.
- Separate multiple slides by at least one slotted channel inside the slide mailer.
- Avoid placing slides in the first or last slotted channel of the slide mailer or coplin jar. Slides with tissues in the first or last position may get scratched if facing the jar wall.
- Up to two slides can be processed per slide mailer and up to four slides can be processed per coplin jar. Note that buffer volume in slide mailers and coplin jars does not change if processing more than one slide.

Quencher Removal

- Quencher Removal is critical for accurate pathological assessment of tissue sections.
- Quencher Removal Solution should be prepared fresh before use for optimal results.
- The performance of the solution declines after 10 minutes. The solution will turn from clear to white in color, indicating it should be made again.
- DO NOT use alternative chemicals for Quencher Removal other than those specifically recommended in this Demonstrated Protocol.
- Customers may purchase Sodium Hydrosulfite from any supplier and in small quantities (e.g. 1 g), if preferred.
- 10x Genomics cannot provide handling instructions for Sodium Hydrosulfite outside of the following:
 - Handle in a fume hood
 - Follow institutional guidelines

H&E Staining

- Test out H&E staining protocol on a serial section placed on a blank slide before staining experimental Xenium slides.
 - If staining is too light or dark, check reagent concentration and consider optimizing staining time of Hematoxylin and Eosin solutions.
- Use fresh staining reagents as the pH of reagents (e.g. Hematoxylin) can change over time, leading to suboptimal staining.
- Filter Hematoxylin and Eosin solutions using filter paper before starting H&E staining.
- Replace ethanol frequently as any dye leftover in ethanol can stain tissues.

Protocol Steps & Timing

~2.5 h

Steps	Timing	Stop & Store
H&E Staining		
1.1 Buffer Preparation - Quencher Removal	15 min	
1.2 Cassette Removal	10 min	
1.3 Quencher Removal	15 min	 4°C ≤ 2 days
1.4 Buffer Preparation - H&E Staining	15	
1.5 H&E Staining	1 h	
1.6 Coverslipping	35 min	 4°C ≤ 3 days or room temperature (in the dark) long-term


1. H&E Staining

1.0 Overview


This chapter provides guidance on H&E staining of Xenium slides containing FFPE or fresh frozen tissue sections. First, Autofluorescence Solution is removed from tissue sections in a destaining step. Tissue sections are then stained in Hematoxylin and Eosin staining solutions and coverslipped to prepare for imaging. Proceed to imaging following the Post-Xenium Analyzer H&E Staining workflow.



Protocol Overview

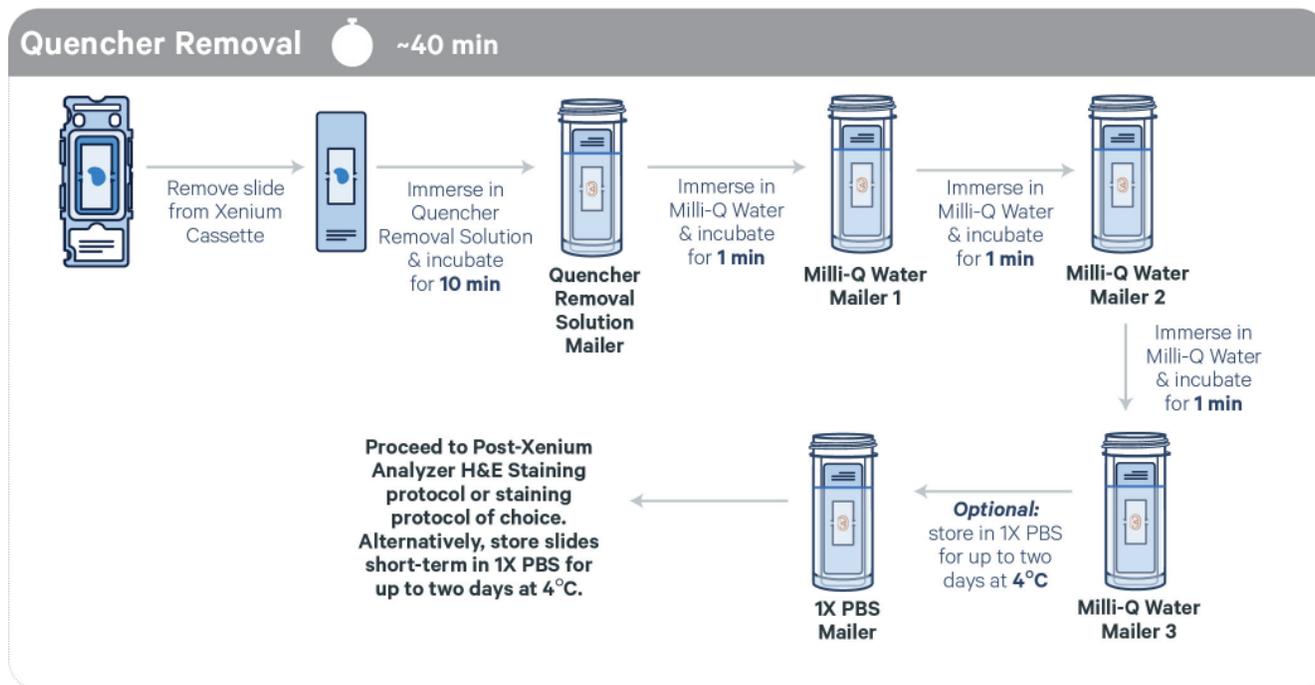

! *Quencher Removal is critical for accurate pathological assessment. For assessing tissue integrity and general tissue morphology, skipping Quencher Removal may be acceptable.*

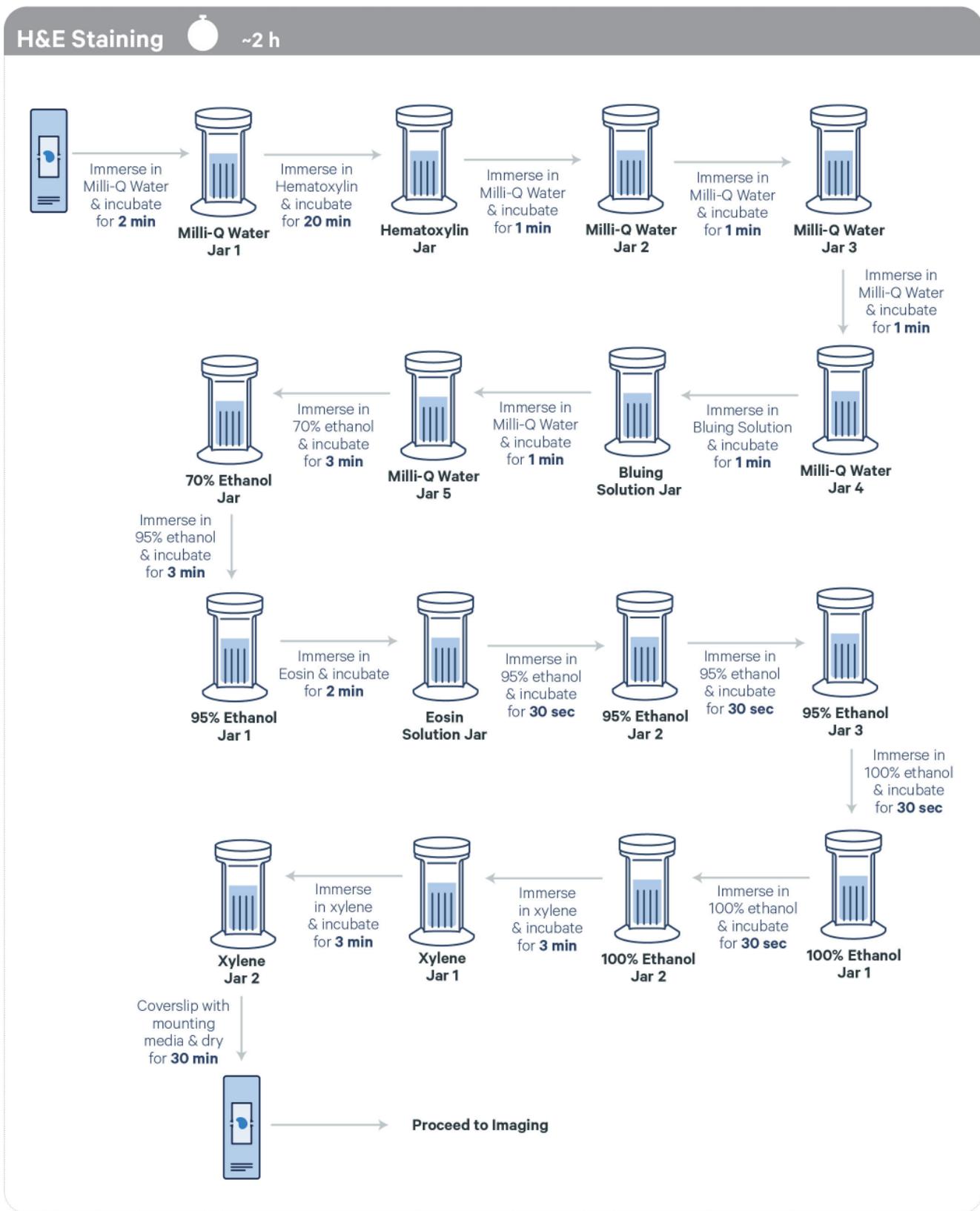

Get Started - Quencher Removal

Quencher Removal Items	10x PN	Preparation & Handling	Storage	
Obtain				
<input type="checkbox"/>	Sodium hydrosulfite	-	-	Ambient
<input type="checkbox"/>	Forceps	-	-	Ambient
<input type="checkbox"/>	Slide Mailers	-	-	Ambient
<input type="checkbox"/>	Milli-Q Water	-	-	Ambient
<input type="checkbox"/>	PBS <i>(optional)</i>	-	-	Ambient
<input type="checkbox"/>	Xenium Slides (2 pack) with FFPE or FF tissue sections	3000941	Retrieve from Xenium Analyzer and prepared according to Xenium Analyzer User Guide (CG000584).	4°C post Xenium-Analyzer run

1.1 Quencher Removal Preparation



Prepare buffers in appropriate sized conical tube or bottle and transfer carefully to corresponding slide mailer. If slide mailer internal volume is more than 12 ml, scale Quencher Removal Solution to fit mailer volume.

- a. Prepare Quencher Removal Solution **immediately** before use according to the table below.



Quencher Removal Solution must be prepared fresh before each use for optimal results, as performance of the solution declines after 10 minutes. Prepare solution in a fume hood, according to the instructions below, due to the hazardous nature of sodium hydrosulfite.

Weigh 17.4 mg of sodium hydrosulfite in a fume hood and add to conical tube containing 10 ml Milli-Q Water. Cap tube and vortex 10 sec to dissolve. Vortexing can be performed outside of fume hood if tube is tightly capped. Maintain at room temperature.

Quencher Removal Solution				
Items		Molecular Weight	Final	Total Amount (X units)
<input type="checkbox"/>	Milli-Q Water	-	-	10.0 ml
<input type="checkbox"/>	Sodium Hydrosulfite	174.104 g/mol	10 mM	17.4 mg
<input type="checkbox"/>	Total	-	-	10.0 ml

- b. Prepare 1X PBS (optional) for storage of slides post-Quencher Removal. Add reagents in the order listed. Invert gently to mix. Maintain at room temperature.

1X PBS (optional)				
Items		Molecular Weight	Final	Total Amount (ml)
<input type="checkbox"/>	Milli-Q Water	-	-	9.0
<input type="checkbox"/>	PBS	10X	1X	1.0
<input type="checkbox"/>	Total	-	-	10.0

- c. Prepare four total slide mailers (and one optional mailer) for Quencher Removal steps.

For H&E Staining		
Items	Preparation & Handling	
<input type="checkbox"/>	Quencher Removal Solution	Label one slide mailer as Quencher Removal Mailer. Dispense 10 ml Quencher Removal Solution.
<input type="checkbox"/>	Milli-Q Water	Label three slide mailers as Milli-Q Water Mailer 1, 2, and 3. Dispense 10 ml Milli-Q Water in each.
<input type="checkbox"/>	1X PBS (optional)	Label one slide mailer as 1X PBS Mailer. Dispense 10 ml 1X PBS.

1.2 Cassette Removal

- a. Retrieve the assembled Xenium Cassette containing FFPE or fresh frozen tissue sections from the Xenium Analyzer.
- b. Using a pipette, remove all PBS-T from well corners.
- c. Remove slide from Xenium Cassette.



See *Tips & Best Practices* for guidance on Xenium Cassette Removal.

Remove Slide from Xenium Cassette




1.3 Quencher Removal

Quencher Removal steps should be performed in a fume hood due to the hazardous nature of sodium hydrosulfite. Quencher Removal is required for any staining protocol.

- a. Gently immerse slide in the Quencher Removal Solution Mailer and incubate for **10 min** at **room temperature**.

When immersing slides, ensure that the tissue sections are completely submerged.

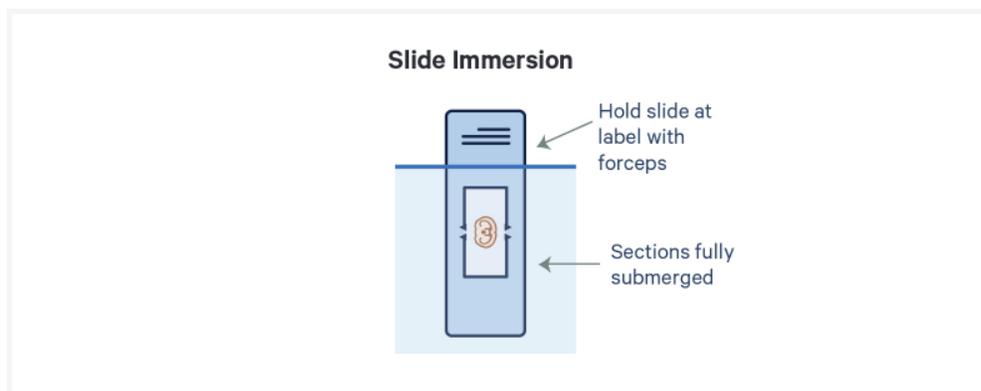

- b. Gently immerse slide in Milli-Q Water Mailer 1 and incubate for **1 min** at **room temperature**.
- c. Gently immerse slide in Milli-Q Water Mailer 2 and incubate for **1 min** at **room temperature**.
- d. Gently immerse slide in Milli-Q Water Mailer 3 and incubate for **1 min** at **room temperature**.
- e.  Proceed to Post-Xenium Analyzer H&E Staining protocol on the following page or alternative staining protocol, if desired. Slides may be stored temporarily in a slide mailer containing 1X PBS at **4°C** for up to two days until ready to stain. If proceeding to Post-Xenium In Situ Applications: Immunofluorescence, H&E, and Visium CytAssist Spatial Gene Expression Technical Note (CG000709), then proceed **immediately**.

Get Started - H&E Staining

H&E Staining Items		10x PN	Preparation & Handling	Storage
Obtain				
<input type="checkbox"/>	Hematoxylin	-	-	Ambient
<input type="checkbox"/>	Eosin	-	-	Ambient
<input type="checkbox"/>	Bluing Reagent	-	-	Ambient
<input type="checkbox"/>	Mounting Media	-	-	Ambient
<input type="checkbox"/>	Xylene	-	-	Ambient
<input type="checkbox"/>	Ethanol	-	Prepare Ethanol dilutions using Milli-Q water.	Ambient
<input type="checkbox"/>	Forceps	-	-	Ambient
<input type="checkbox"/>	Coplin Jars/Staining Dishes	-	-	Ambient
<input type="checkbox"/>	Milli-Q Water	-	-	Ambient

1.4 H&E Staining Preparation



Prepare buffers in appropriate sized conical tube or bottle and transfer carefully to corresponding coplin jar.

- a. Filter Hematoxylin & Eosin solutions using filter paper before starting H&E Staining protocol.
- b. Prepare sixteen total coplin jars for H&E Staining steps.

For H&E Staining		
Items	Preparation & Handling	
<input type="checkbox"/> Hematoxylin Solution	Label one coplin jar as Hematoxylin Jar. Fill to capacity with Mayer's Hematoxylin Solution.	
<input type="checkbox"/> Bluing Solution	Label one coplin jar as Bluing Solution Jar. Fill to capacity with Bluing Solution.	
<input type="checkbox"/> 70% Ethanol	Label one coplin jar as 70% Ethanol Jar. Fill to capacity with 70% ethanol.	
<input type="checkbox"/> 95% Ethanol	Label three coplin jars as 95% Ethanol Jar 1, 2, and 3. Fill to capacity with 95% ethanol.	
<input type="checkbox"/> Eosin Solution	Label one coplin jar as Eosin Solution Jar. Fill to capacity with Eosin Solution.	
<input type="checkbox"/> 100% Ethanol	Label two coplin jars as 100% Ethanol Jar 1 and 2. Fill to capacity with 100% Ethanol.	
<input type="checkbox"/> Xylene	Label two coplin jars as Xylene Jar 1 and 2. Fill to capacity with Xylene.	
<input type="checkbox"/> Milli-Q Water	Label five coplin jars as Milli-Q Water Jar 1, 2, 3, 4, and 5. Fill to capacity with Milli-Q Water.	



Alternatively, a slide staining dish can be used in place of a coplin jar. Adjust volumes accordingly. Use xylene-resistant dishes, gloves, and forceps during workflow.

1.5 H&E Staining

H&E Staining steps should be performed in a fume hood due to the hazardous nature of xylene. Xylene jars should be covered at all times to prevent evaporation.

- a. Gently immerse slide in the Milli-Q Water Jar 1 for **2 min** at **room temperature**.



Water immersions may be performed in glass beakers containing Milli-Q water, if preferred.

- b. Gently immerse slide in the Hematoxylin Solution Jar for **20 min** at **room temperature**.



Substitution of Hematoxylin with a brand other than the recommended option may lead to tissue staining that is darker than anticipated following 20 min incubation.


- c. Gently immerse slide in the Milli-Q Water Jar 2 for **1 min** at **room temperature**.
- d. Gently immerse slide in the Milli-Q Water Jar 3 for **1 min** at **room temperature**.
- e. Gently immerse slide in the Milli-Q Water Jar 4 for **1 min** at **room temperature**.
- f. Gently immerse slide in the Bluing Solution Jar for **1 min** at **room temperature**.
- g. Gently immerse slide in the Milli-Q Water Jar 5 for **1 min** at **room temperature**.
- h. Gently immerse slide in the 70% Ethanol Jar for **3 min** at **room temperature**.
- i. Gently immerse slide in the 95% Ethanol Jar 1 for **3 min** at **room temperature**.
- j. Gently immerse slide in the Eosin Solution Jar for **2 min** at **room temperature**.
- k. Gently immerse slide in the 95% Ethanol Jar 2 for **30 sec** at **room temperature**.
- l. Gently immerse slide in the 95% Ethanol Jar 3 for **30 sec** at **room temperature**.

- m.** Gently immerse slide in the 100% Ethanol Jar 1 for **30 sec** at **room temperature**.
- n.** Gently immerse slide in the 100% Ethanol Jar 2 for **30 sec** at **room temperature**.
- o.** Gently immerse slide in the Xylene Jar 1 for **3 min** at **room temperature**.
- p.** Gently immerse slide in the Xylene Jar 2 for **3 min** at **room temperature**.

1.6 Coverslipping

Before mounting the coverslip, ensure that the slide is dry. Moisture on the surface of the slide may result in faulty mounting. Wipe away any residual droplets with a lint-free laboratory wipe.

- a. Place slide on a flat, clean, non-absorbent work surface.
- b. Using a **wide-bore** pipette tip, add **150-200 μ l** mounting media to uniformly cover all tissue sections on the slide.
- c. Apply the coverslip at an angle on one end of the slide. Slowly lower the coverslip, without introducing bubbles. Allow mounting media to spread and settle.
- d. If needed, remove any large excess of mounting media by carefully wicking away from the edge of the coverslip with a lint-free laboratory wipe. Be careful not to move coverslip and disturb the tissue.
- e. Dry the coverslipped slide for **30 min** at **room temperature**.
- f. Once coverslipping is complete, proceed with imaging. Samples may be stored temporarily at **4°C** and imaged within 3 days after coverslipping. For long-term storage after imaging, maintain at **room temperature in the dark**.

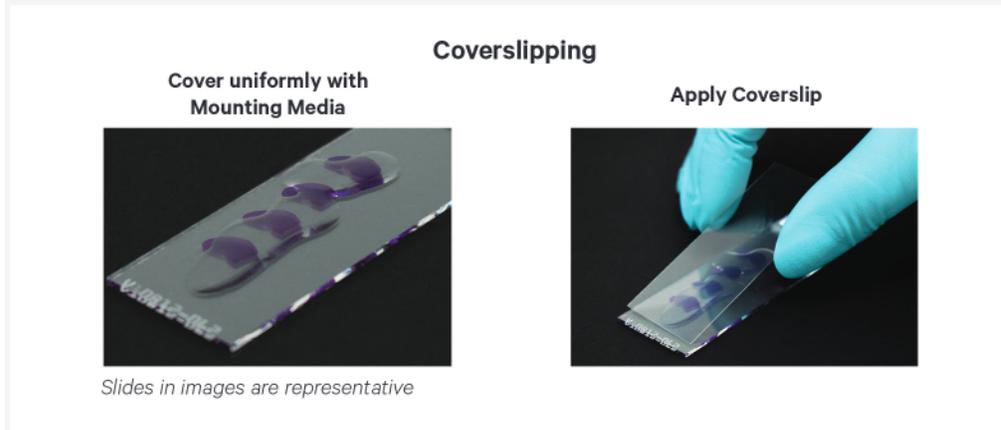

2. Troubleshooting

Dry Tissue

- Tissue shrinking is likely due to slides drying during the H&E staining workflow, and not related to Quencher Removal.
- Tissue drying can occur when:
 - Slides are not kept hydrated in between the Xenium workflow and H&E staining.
 - Slides are kept too long in Ethanol during Ethanol washes.

Using tap water vs. Milli-Q water will not impact slide drying during H&E staining.

Document Revision Summary

Document Number	CG000613
Title	Xenium In Situ Gene Expression - Post-Xenium Analyzer H&E Staining Demonstrated Protocol
Revision	Rev B
Revision Date	December 2023

Specific Changes	<p>Updated sodium hydrosulfite Part Number (157953-5G) in Specific Reagents & Consumables (page 4). Added alternative mounting media option (Cytoseal Mountant; 60, 280 and XYL, Fisher Scientific, PN 22-050-262) to Specific Reagents & Consumables (page 4). Updated Slide Storage guidance (page 7). Added new Quencher Removal section to Tips & Best Practices (page 9). Added new H&E Staining section to Tips & Best Practices (page 9). Added Troubleshooting section for Dry Tissue (page 24).</p>
General Changes	<p>Updated for general minor consistency of format, language, and terms throughout.</p>

© 2023 10x Genomics, Inc. (10x Genomics). All rights reserved. Duplication and/or reproduction of all or any portion of this document without the express written consent of 10x Genomics, is strictly forbidden. Nothing contained herein shall constitute any warranty, express or implied, as to the performance of any products described herein. Any and all warranties applicable to any products are set forth in the applicable terms and conditions of sale accompanying the purchase of such product. 10x Genomics provides no warranty and hereby disclaims any and all warranties as to the use of any third-party products or protocols described herein. The use of products described herein is subject to certain restrictions as set forth in the applicable terms and conditions of sale accompanying the purchase of such product. A non-exhaustive list of 10x Genomics' marks, many of which are registered in the United States and other countries can be viewed at: www.10xgenomics.com/trademarks. 10x Genomics may refer to the products or services offered by other companies by their brand name or company name solely for clarity, and does not claim any rights in those third-party marks or names. 10x Genomics products may be covered by one or more of the patents as indicated at: www.10xgenomics.com/patents. The use of products described herein is subject to 10x Genomics Terms and Conditions of Sale, available at www.10xgenomics.com/legal-notice, or such other terms that have been agreed to in writing between 10x Genomics and user. All products and services described herein are intended FOR RESEARCH USE ONLY and NOT FOR USE IN DIAGNOSTIC PROCEDURES.

The use of 10x Genomics products in practicing the methods set forth herein has not been validated by 10x Genomics, and such non-validated use is NOT COVERED BY 10X GENOMICS STANDARD WARRANTY, AND 10X GENOMICS HEREBY DISCLAIMS ANY AND ALL WARRANTIES FOR SUCH USE. Nothing in this document should be construed as altering, waiving or amending in any manner 10x Genomics terms and conditions of sale for the Chromium Controller or the Chromium Single Cell Controller, consumables or software, including without limitation such terms and conditions relating to certain use restrictions, limited license, warranty and limitation of liability, and nothing in this document shall be deemed to be Documentation, as that term is set forth in such terms and conditions of sale. Nothing in this document shall be construed as any representation by 10x Genomics that it currently or will at any time in the future offer or in any way support any application set forth herein.

Contact:

Email: support@10xgenomics.com
10x Genomics
6230 Stoneridge Mall Road
Pleasanton, CA